



Medical-Biological
Research & Technologies

WORLD OF BIOTECH-INNOVATICA



Catalogue 2023

WELCOME TO THE BIOSAN PRODUCT CATALOGUE 2023

For more than 30 years, the **mission of BIOSAN** has remained unchanged, reducing the risks of sample preparation by offering lines of laboratory equipment based on the most advanced scientific and industrial technologies.

The **catalogue of 2023** traditionally represents a wide range of laboratory equipment for carrying out both the simplest stages of sample preparation, such as mixing, shaking, centrifuging, incubation, and for implementing more complex and advanced stages — cell cultivation, nucleic acid isolation, as well as instruments for final result analysis.

Biosan look forward introducing to you a new range of sample preparation instruments: new multi-functional centrifuge **LMC-56** with automatic rotor recognition, reciprocal homogenizer **RCP-24**, new shaker-incubator **ES-20/80C** with integrated cooling capability and new CO₂ tube roller **CTR-6**.

For more convenience of selection, we have created a number of workflow description, which, we hope, will help you select the most appropriate product for your needs. ELISA, Nucleic Acid extraction and amplification, Microbial cell cultivation — these are only a few of them. From **Assist** pipettes to **3D-IW8** automated ELISA plate washers, from **MagSorb-16** magnetic stand for nucleic acid isolation to automated 12 channel nucleic acid processor **Biomagpure 12 Plus**. Our devices will allow you to properly and effectively prepare your sample and provide you with reliable results — like **Bioquant 96** Real-Time PCR instrument or **HiPo MPP-96** ELISA plate reader.

Instrument models are constantly being improved in accordance with the wishes and comments of our customers. We pay special attention to reliability, ergonomics and ease of use.

SAMPLE PREPARATION SOLUTIONS

We are focused on solving the problems of reproducibility of experimental data, the cause of which is usually the sample preparation of a biomaterial. The absence of strict rules for sample preparation leads to the greatest number of errors in implementing of the methodology.

Errors accumulate in connection with:

1. A significant decrease in the volumes of reactants (from millilitres to microliters);
2. The lack of intermediate temperature logistics of the technological process (temperature shelf);
3. The characteristics of mixing micro-quantities of reagents;
4. The lack of laboratory air decontamination systems in the process;
5. Inappropriate storage of cellular material not regulated by the method.

We offer solutions for sample preparation only after we make sure that they are ideal for eliminating the above-mentioned errors.

CUSTOMER SUPPORT

We are attentive to all customer requests. Specialists of the company promptly provide warranty and post-warranty service, and solve problems that may arise from users both at the stage of ordering equipment and during operation and maintenance. We are always happy to assist you in developing skills in the operation and maintenance of our products.

Additional information about the products can be found on the website www.biosan.lv, including video of the products demonstrating the functional characteristics. Electronic brochures, catalogue and user manuals are also available for download.

VISION

We plan to continue our continuous improvement in order to remain your partner and expert in both Life Science research and medical diagnostics. This will allow us to develop new promising products, while remaining in the same row with the world's leading bioengineering companies.

We will be sincerely happy if you are interested in BIOSAN products. Thank you for your cooperation!

Vasily Bankovsky

Ph.D., Biology

President of BIOSAN



World of BIOTECH-INNOVATICA

“Most of the great innovations arose from the interaction of creative personalities with teams that managed to realize their ideas.”

— Walter Isaacson



Vasily Bankovsky, Ph.D. (Biology), Head of R&D Department, Chairman of the Board at Biosan

The concept of development for Biosan called World of Biotech-innovatica. Four planetary systems with satellites — devices revolve around Terra Innovatica (biomaterial under research). We have marked out four planets — 4 contemporary levels of life science research, medical and veterinary diagnostics:

- 1. Terra Genomica** — level of genes (DNA-analysis, oligonucleotide and mononucleotide polymorphism — ONP, SNP);
- 2. Terra Immunologica** — level of immunology (detection of polymorphism of antibodies and immune response);
- 3. Terra Biochemica** (metabolomics) — level of metabolism products and ferment activity;
- 4. Terra Cellomica** — level of cellular morphogenesis (cellular polymorphism).

The distance from the planet orbitals to Terra Innovatica corresponds to the time of disease detection at each level (from one week, as in the case of DNA-analysis, to several years, when the changes can be tracyomorphism at the level of genes leads to the manifestation of polymorphism at all higher levels, it results in the ambiguity (if not more) of any decision made based on the obtained data. The definition comprising the polymorphism of norm and abnormality (disease) is not yet available, this experience is still being gained, hence, the multianalysis technology, though expensive, is the only solution as of today.

Although the classic determinism in life science research and diagnostics has finally yielded its position to the stochastic one, there are still no instrumental solutions, that allow channelling our new knowledge into informed and unambiguous decisions. This is the real situation; these are the temporary sacrifice of progress.

Biosan is the only company in the World of Biotech-innovatica, which develops, produces and distributes instrument lines for all four levels. These satellites of four planets are specialised devices providing the instrumental basis for multilevel analysis, whereas the reagent sets make these satellites move. Consequently, by the World of Biotech-innovatica we mean the direction of Biotechnology, responsible for the development of multilevel analysis sets (product lines). In the future perspective, multianalysis chips may appear with the development of chip technologies, allowing to unify all the technologies described above in one chip.

I am pleased to point out that many of our ideas and products have been developed as a result of the long-standing cooperation between scientists and developers of Biosan with universities, as well as with academic institutes and institutes of applied sciences and our company customers worldwide.

All our inventions resulted from joint efforts, and today we are still open for collaboration. We will be delighted if the result of our work — which has already received wide recognition of the scientific community — would also be of interest for you, particularly if it would serve as yet another starting point for the development of innovative biotechnologies and appearance of new planets and their satellites in the sky of the World of Biotech-innovatica.

Sincerely,
Vasily Bankovsky, Ph.D. (Biology)
Head of R&D Department
Biosan, Chairman of the Board

World of BIOTECH-INNOVATICA

TERRA CELLOMICA

Research and diagnostics at the level of cellular morphogenesis (cellular polymorphism);

TERRA BIOCHEMICA (metabolomics)



Research and diagnostics of metabolism products and enzyme activity;

TERRA IMMUNOLOGICA

Research and diagnostics at the level of immunology (detection of polymorphism of antibodies and immune response);

TERRA GENOMICA

Research and diagnostics at the level of genes (DNA-analysis, sequencing oligonucleotide and mononucleotide polymorphism — ONP, SNP).

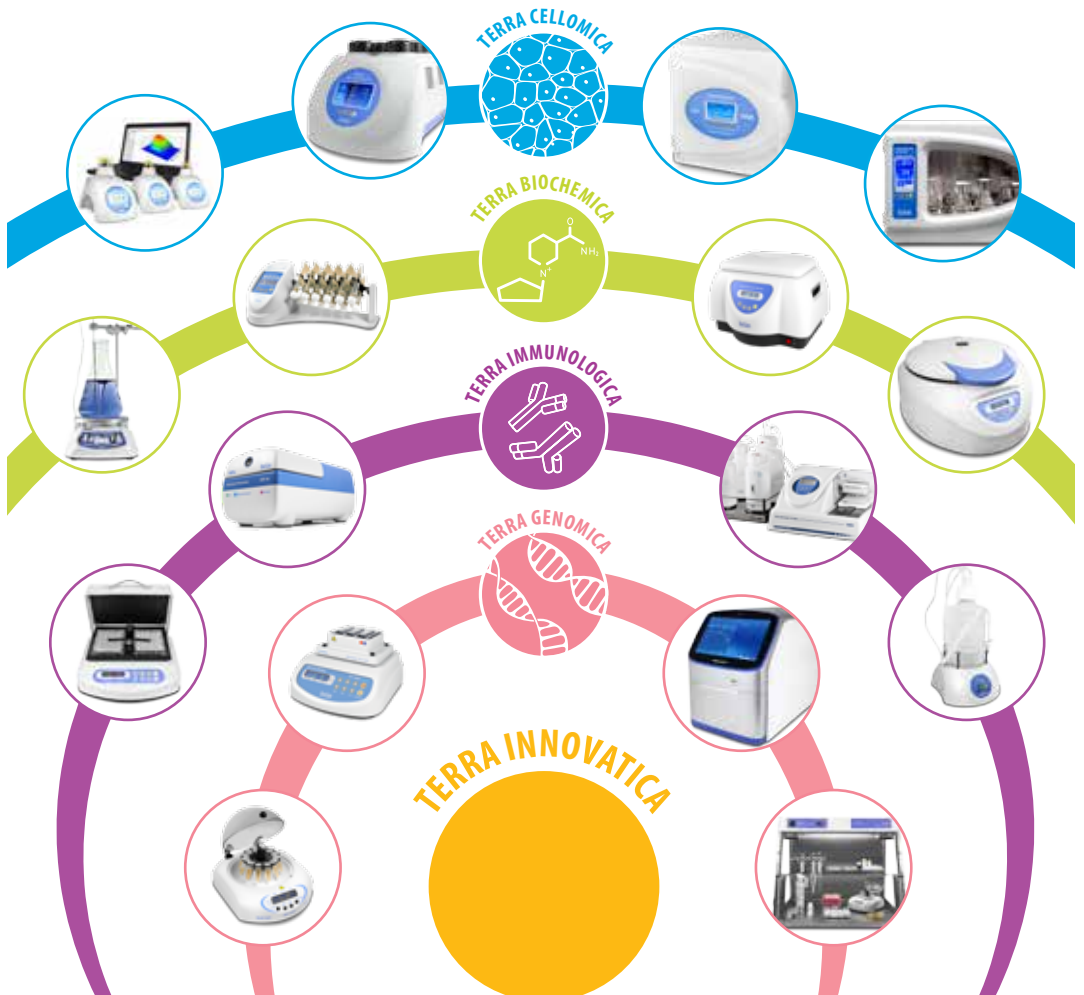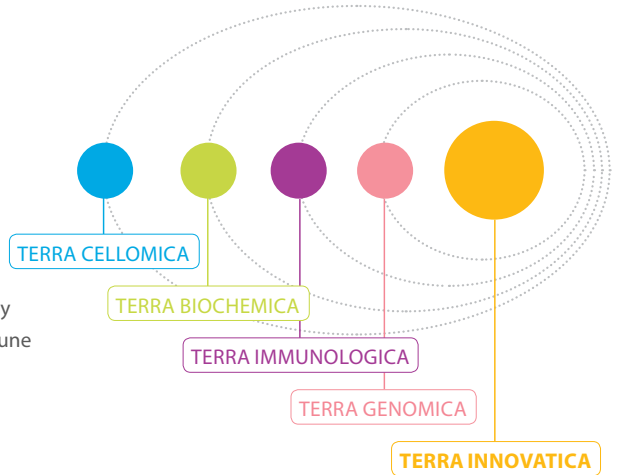

NEW PRODUCTS AND ANNOUNCEMENT

NEW PRODUCTS

- Photometer:** DEN-600
- Centrifuge:** LMC-56
- Shaker-incubator with cooling:** ES-20/80C
- Homogenizer:** RCP-24
- CO₂ shaker:** CPS-20
- Shelves for DNA/RNA UV-cleaner boxes:** P-5, F-1
- Sticky mats:** SPML, SPM
- Combined Colorimetric MicroArray and ELISPOT Imager:** ImagerBio C
- Fluorescent MicroArray Imager:** ImagerBio F
- CO₂ Tube Roller:** CTR-6



ANNOUNCEMENTS

Please visit the news section of our website for the latest product announcements.

IVD REGISTERED PRODUCTS

- Photometer:** DEN-600m
- Washer:** IW-8m

PRODUCT LINE

COVID-19

GENERAL LAB EQUIPMENT

ROCKERS, SHAKERS, ROTATORS, VORTEXES, HOMOGENIZER 15

- Rockers:** MR-1, MR-12
- Shakers:** 3D, Multi Bio 3D, PSU-10i, PSU-20i, MPS-1, PSU-2T
- Rotators:** Bio RS-24, Multi Bio RS-24, Multi RS-60
- Vortexes:** V-1 plus, V-32, MSV-3500
- Homogenizer:** RCP-24

THERMO-SHAKERS 35

- PST-60HL, PST-60HL-4, PST-100HL, TS-DW, TS-100, TS-100C, TS-100C Smart

MINICENTRIFUGES-VORTEXES 43

- FV-2400, FVL-2400N, MSC-3000, MSC-6000, CVP-2

CENTRIFUGES 50

- Microspin 12, LMC-3000, LMC-4200R
- Rotors for LMC-3000 and LMC-4200R

DRY BLOCK THERMOSTATS 59

- Bio TDB-100, TDB-120, CH-100, CH 3-150
- Dry block heating systems with interchangeable blocks

WATER BATHS 67

- Stirred water baths:** WB-4MS
- Unstirred water baths:** JB Academy, JB Nova, SUB Aqua Pro
- Circulating baths:** LSB Aqua Pro range, OLS26, Optima™ series, Optima™ R series, LT ecocool™, SBB series

MAGNETIC STIRRERS, OVERHEAD STIRRERS 85

- Magnetic Stirrers:** MS-3000, MMS-3000, MSH-300, Intelli-Stirrer MSH-300i
- Overhead Stirrer:** MM-1000

BIOSAFETY AIR, SURFACE – UV-cleaner boxes, UV-cleaner recirculators 96

- UV-cleaner box:** UVC/T-AR, UVC/T-M-AR, UVT-B-AR, UVT-S-AR
- UV Cleaner--Recirculators:** UVR-M, UVR-Mi
- DNA/RNA Decontamination Solution, Spray:** PDS-250, PDS-10

WATER PURIFICATION SYSTEMS 102

- Labacqua Trace, Labacqua HPLC, Labacqua Bio

DENSITOMETERS, PHOTOMETER	105
---------------------------------	-----

Densitometers: DEN-1, DEN-1B

Photometer: DEN-600

ASPIRATORS, PIPETTES	109
----------------------------	-----

Aspirator: FTA-1, FTA-2i

Pipettes: Assist series; Assistboy

WASHERS	144
---------------	-----

IW-8, 3D-IW8

BIOPROCESSING

CO ₂ INCUBATOR	130
---------------------------------	-----

S-Bt Smart Biotherm, Compact CO₂ Incubator

CPS-20, CO₂ shaker

PERSONAL BIOREACTORS	108
----------------------------	-----

RTS-1, RTS-1C, RTS-8, RTS-8 Plus

SHAKERS–INCUBATORS	124
--------------------------	-----

ES-20, ES-20/60, ES-20/80

LAB DIAGNOSTICS

DNA/RNA PURIFICATION	136
----------------------------	-----

DNA/RNA Extraction Lines

BioMagPure 12 Plus, Workstation for automated nucleic acid purification and reagents

MagSorb-16, Magnetic rack for manual nucleic acid extraction

REAL TIME PCR DETECTION	140
-------------------------------	-----

PCR Analysis Line

BioQuant-96 and reagents

IMMUNODIAGNOSTICS	143
-------------------------	-----

ELISA Line

IW-8, Intelispeed Washer

3D-IW8, Inteliwasher

ImagerBio C, Combined Colorimetric MicroArray and ELISPOT Imager

ImagerBio F, Fluorescent MicroArray Imager

MPP-96 HiPo, Microplate Photometer

APPLICATIONS

Article links	157
---------------------	-----

How to choose a Proper Shaker, Rocker, Vortex	158
---	-----

Product line examples	160
-----------------------------	-----

GENERAL INFORMATION

GENERAL INFORMATION ABOUT BIOSAN POLICY	156
---	-----

COVID-19

The global coronavirus (COVID-19) outbreak marks the necessity of fast and reliable sample preparation as well as safe working environment. To facilitate products selection, here are shown ready product lines to start COVID-19 sample analysis rapidly.

Preventive

Decontamination



UVR-M/UVR-Mi
UV Cleaner-Recirculators



PDS-250, PDS-10L
DNA/RNA decontamination solution

Air

Surface

ELISA based detection methods

Steps

Incubation



PST-60HL-4
Plate Shaker-Thermostat



PST-60HL
Plate Shaker-Thermostat



PSU-2T
Mini-shaker
for immunology

Washing



Assist
Pipette series



Assistboy
Pipette
controller



FTA-2i
Aspirator
with Trap Flask



3D-IW8
Inteliwasher

Detection and analysis



HiPo MPP-96
Microplate
Photometer



Quant Assay
Software for
HiPo MPP-96

NA based detection methods

Manual preparation



Working area

UVT-5-AR
DNA/RNA UV-cleaner box

Resuspension of probes and reagents



V-1 Plus
Personal Vortex



MSC-3000
Centrifuge/Vortex
Multispin

Sample lysis




TDB-120
Dry block thermostat



TS-100C Smart
Programmable Thermo-Shaker
with cooling for microtubes
and PCR plates

Sample wash

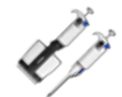

Assist
Pipette series



FTA-2i
Aspirator with Trap Flask

Elution of NA



MagSorb 16
Magnetic Rack for
Manual Nucleic Acid
Extraction



Microspin 12
High-speed
Mini-centrifuge



Biomagpure 12 Plus
Compact Bench-Top
Robotic Workstation for
Automated Nucleic Acid
Purification (+ reagents)

Sample detection and analysis

Preparation



CVP-2
Plate centrifugation
and vortexing



MSC-6000
Tube centrifugation
and vortexing

Real time PCR detection



BioQuant-96
Real time
PCR instrument



abTes COVID-19 qPCR kit
(Reagents are CE-IVD certified
– among WHO listed kits.)



Bucket rotor, only for LMC-56

LMC-56 NEW see page 54

Centrifuge

LMC-56 is a bench-top, low-speed multifunctional centrifuge for routine laboratory tasks. Perfect for biochemical, cytological, immunological analyzes. It can be used for sample preparation in research, biomedical and genetic laboratories. The precursor of this centrifuge is Well-known among users in Europe laboratory centrifuge LMC-3000.

New centrifuge has significantly improved features:

- Higher speed — maximum 6,000 rpm or 3,160 g;
- New bucket rotors, having maximum total capacity up to 400 ml;
- A wide selection of rotors — up to 13;
- Automatic recognition of rotors and setting the maximum speed for each rotor type;
- Improved ventilation of the working chamber to reduce tube heating during centrifugation;
- Several types of acceleration and braking;
- Auto-stop in case of imbalance.





CPS-20 **NEW** [see page 132](#)

CO₂ Shaker

CO₂ Shaker **CPS-20** provides regulated orbital motion of the platform and is designed for use in Biosan's S-Bt CO₂ Incubator. **CPS-20** is specifically designed in harsh environments such as CO₂ and humidity and provides reproducible results for cell culture growth. A choice of five interchangeable platforms provide the possibility of performing various procedures and techniques in various cultivation vessels.

CPS-20 incorporates a brushless motor with a guaranteed service life of up to 35,000 hours and a direct-drive system, ensuring the most reliable, stable long-time operation (up to 30 day nights). The unit is equipped with a triple eccentric mechanism for platform motion that provides supreme balancing characteristics, superior reliability and quiet operation. The specially designed remote controller allows for the protection of electronics from a CO₂ incubator environment, as well as, the remote control minimizes interference with the incubator environment and the ongoing experiment.



CTR-6 **NEW** [see page 129](#)

CO₂ Tube Roller

CO₂ Tube Roller **CTR-6** provides regulated rocking and rolling of maximum up to 6 rollers and is designed for use specifically in CO₂ incubators. **CTR-6** is specifically designed for use in harsh environments such as CO₂ and humidity and provides reproducible results for cell culture growth. Possibility to remove rollers makes the unit flexible and allows for performing various procedures and techniques in various cultivation vessels. The specially designed remote controller allows for protection of electronics from CO₂ incubator environment, as well as, does not interfere with the experiment.





ES-20/80C **NEW** [see page 126](#)

Shaker-incubator with cooling

Orbital Shaker-Incubator with cooling **ES-20/80C** for biotechnological and pharmaceutical laboratories is a next-generation, professional category equipment.

ES-20/80C duplicates the functionality of ES-20/80 and also uses a Peltier element to cool the camera. A built-in heat-resistant brushless fan provides precise temperature distribution inside the chamber (from 12.5 °C below ambient up to +80 °C). Additionally, excellent sample temperature uniformity of ± 0.2 °C at 37 °C is achieved.



RCP-24 **NEW** [see page 34](#)

Homogenizer

Reciprocal Homogenizer **RCP-24**, a bench-top mechanical device designed for mixing, grinding, homogenizing and emulsifying biological objects in microtubes by vigorously mixing by reciprocal motion with various beads for sample preparation for subsequent academic, pharmaceutical, biotechnological or biomedical studies.





DEN-600 **NEW** see page 108

Photometer

DEN-600 is a compact, portable, rechargeable battery-powered photometer. It comprises of 600 nm wavelength optical system, which enables to apply:

- 1) OD₆₀₀ method that estimates total number of cells.
- 2) McFarland (McF) turbidity measurement method.
- 3) Bradford protein assay method for protein concentration measurement.
- 4) Other methods that can be adjusted or optimized using 600 nm wavelength.

The device serves as an inexpensive alternative to a spectrophotometer, which is commonly used for these applications. Because **DEN-600** is battery powered and compact, it can be comfortably located in a biosafety cabinet, anaerobic chamber or quickly moved to another lab room. Additionally, the vessel holding mechanism allows accommodating round bottom, conical vials or falcon tubes, therefore enabling to measure the absorbance (Abs) and turbidity in Abs, OD and McFarland units.

USB connectivity and DEN software allow for data transfer, data processing and calculation, software calibration for Bradford protein assay method or a custom calibration for a specifically suitable vessel.



P-5



F-1


P-5, F-1 **NEW** see page 100

Shelves for DNA/RNA UV-cleaner boxes

Two types of shelves have been developed for DNA/RNA UV-cleaner boxes to increase the box's effective area:

P-5 — shelf-holder for 5 pipettes and **F-1** flat shelf.

On the **F-1** shelf, you can place laboratory glassware, reagents and other items that are convenient to keep in close proximity.

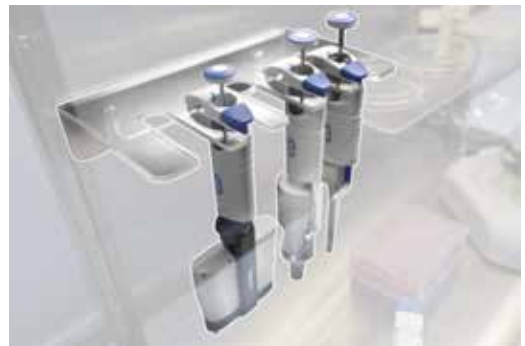



SPM



SPML



IVD REGISTERED PRODUCTS

Market continues to evolve and Biosan following the trends now can proudly announce that we have completed the implementation and successfully certified the Quality Management System according to ISO 13485 standard. As from now on we are proudly offering two medical devices:



DEN-600m

Densitometer

Medical photometer for bacterial culture normalisation. Helps many laboratories across the globe to standardise microbiology examinations.

SPM, SPML **NEW** see page 22

Double-sided adhesive strips and mat

Convenient alternative to traditional steel holders, an easy way to fix tubes, plates, flasks and other laboratory wares on platforms for cultivation, incubation and mixing.

SPML, set of 3 double-sided adhesive strips



SPM, double-sided adhesive mat



IW-8m

Microplate washer

Medical microplate washer. Whenever you need to standardise your ELISA protocol our microplate washer will assist you with unrivalled performance and lowest possible in-well residual volume.



ImagerBio C **NEW** [see page 148](#)

Combined Colorimetric MicroArray and ELISPOT Imager

ImagerBio C is a compact and portable instrument designed for imaging colorimetric Microarrays and Elispots. The instrument by standard operates in brightfield and darkfield illumination modes with bottom and top illuminations respectively. The instrument is able to work with 96-well plate, a 12×8 well strip, as well as microscopic glasses (placed in a special adapter).




ImagerBio F **NEW** [see page 150](#)

Fluorescent MicroArray Imager

ImagerBio F is a compact and portable fluorescent imager designed for imaging fluorescent microarrays. The instrument by standard operates on 2 channels: CY3 and CY5. The instrument is able to work with 96-well plate, a 12×8 well strip, as well as microscopic glasses (placed in a special adapter). The device is available in 2, 3 or 4 channel configuration: DAPI, FITC, CY3, CY5. Qdot options are also available.

MIXING DEVICES:

**ROCKERS, SHAKERS,
ROTATORS, VORTEXES,
HOMOGENIXER**

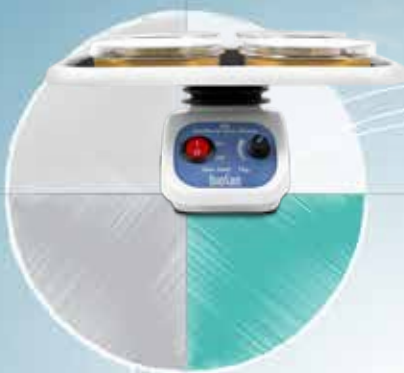


3D

Sunflower Mini-shaker



MPS-1

High-Speed Multi Plate Shaker

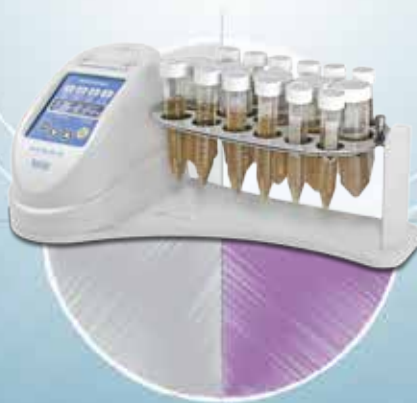

Multi Bio RS-24

Programmable rotator

MR-1, Mini Rocker–Shaker

DESCRIPTION

Mini Rocker-Shaker **MR-1** provides regulated gentle rocking motion of the platform and is ideal for mini gel destaining after electrophoresis, conducting Northern, Southern and Western blot analysis.

Shaker is a compact, noiseless device designed for personal use. Drive and brushless motor allow continuous mixing up to 7 days and ensures reliable, trouble-free operation for more than 2 years.

Non-slip, temperature resistant, silicone mat located on the rocker's platform provides a stable position for vessels during shaking. Optional dimpled PDM mat fixes tubes of different sizes.

The unit is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C.

ACCESSORIES FOR THE STANDARD PLATFORM:

Optional dimpled mat **PDM** prevents different size tubes from rolling around the platform



MR-1



Rocking uni-rotation



Product video

MR-12, Rocker–Shaker

DESCRIPTION

MR-12 Rocker–Shaker provides both soft and intensive mixing of solutions or nutrient media in vessels or plastic bags placed on the platform. Adjustable speed and platform tilt angle allows setting parameters for optimal solution transfer and mixing.

The device is ideal for gel destaining after electrophoresis and homogenisation of bioextraction media. It is optimal for biomolecule hybridisation on strips and staining/destaining procedures. When installed inside a bioincubator it is ideal for growing cells and cell cultures in disposable plastic reactor-bags (working volumes up to 10 litres, media volumes up to 5 litres).

The unit is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C. Low voltage external power supply (12V) provides electrical safety in humid environment.



MR-12



Rocking uni-rotation with adjustable tilt

MR-1 and MR-12, Rocker-Shakers

	MR-1	MR-12
Mixing frequency range	1–30 oscill./min	1–99 oscill./min (increment 1 oscill./min)
Fixed tilt angle	7° (fixed)	0°–10° (increment 1°) (for 1–50 oscill./min) 10° (for 51–99 oscill./min)
Max. continuous operation time	168 h	
Digital time setting	1 min–24 h/non-stop	1 min–99 h 59 min (increment 1 min)/non-stop
Timer sound signal	—	yes
Non-slip silicone mat is supplied as standard	215 × 215 mm	480 × 380 mm
Maximum load	1 kg	5 kg
Display	LED	LCD, 2 × 16 signs
Platform working area	215 × 215 mm	480 × 380 mm
Overall dimensions (W×D×H)	220 × 205 × 120 mm	430 × 480 × 210 mm
Weight	2.1 kg	11.9 kg
Input current/power consumption	12 V, 320 mA/3.8 W	12 V, 1.1A/13 W
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V	

PDM, dimpled mat



MR-1 with PDM dimpled mat



MR-12



ORDERING INFORMATION:

MR-1 with standard platform Bio PP-4S

MR-12 with standard platform PP-480

Optional accessories: for MR-1:

PDM, dimpled mat

Cat. number



BS-010152-AAG

BS-010130-AAI

PDM

3D, Sunflower Mini-Shaker

DESCRIPTION

“Sunflower” **3D** Mini-Shaker provides adjustable three-dimensional smooth rotation of the platform and is designed for mixing blood samples, minigel staining and destaining, sample washing, blot hybridisation reactions.

Mini-Shaker is a compact device with low energy consumption. The use of direct drive and brushless motor allows continuous mixing up to 7 days and ensures reliable, trouble-free operation for many years. Non-slip, temperature resistant, silicone mat located on the shaker's platform provides a stable position for vessels during shaking. The platform is suitable for placing a versatile dimpled PDM mat for different size tubes.

Mini-Shaker can be used in cold rooms or incubators, operating at ambient temperature range +4°C to +40°C.

Basic Plus
Product Class



3D — uni-rotation



Product video

Multi Bio 3D, Programmable mini-shaker («Sunflower» type)

DESCRIPTION

Programmable mini-shaker **Multi Bio 3D** is designed for various applications: hybridization reactions, cell growing, gel washing, soft extraction and homogenisation of biological components in solutions.

Multi Bio 3D provides realization of several types of motion in one module. This option of Biosan instruments essentially extends possibilities and enhances the efficiency of preparation of test samples as well as allows selecting the mixing type according to individual requirements.


Microprocessor control allows performing **1 Orbital 3D rotation** of the platform and but also **2 Reciprocal 3D motion** (of ping-pong type) and **3 Soft vibrating rocking**. These three motion types can be performed separately, pairwise and in cycles, periodically repeating the sequence of three motion types. The shaker is designed for laboratories with increased demands for the quality of mixing, extraction and cell growing processes.

Non-slip, temperature resistant, silicone mat located on the shaker platform provides a stable position for vessels during shaking. Optional dimpled PDM mat fixes tubes of different sizes.

Programmable shaker can be used in cold rooms or incubators, operating at the ambient temperature range +4°C to +40°C.

Premium
Product Class

Ø 22 mm orbit

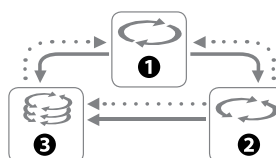

Multi-rotation



Product video

3D Mini-Shaker and Multi Bio 3D, Programmable 3D shaker («Sunflower» type)

	3D	Multi Bio 3D
❶ Speed control range (orbital and reciprocal motion)	5–60 rpm	1–100 rpm
❷ Turning angle (reciprocal motion)	—	0–360° (increment 30°)
❸ Rocking angle (vibro motion)	—	0–5° (increment 1°)
Fixed tilt angle	7°	
Orbit	—	22 mm
Platform working area	215 × 215 mm	
Non-slip silicone mat is supplied as standard		
Maximum continuous operation time	168 h	24 h
Time setting range for ❶ ❷	—	0–250 s
Time setting range for ❸	—	0–5 s
Number of cycles	—	0–125 times
Timer sound signal	—	yes
Maximum load	1 kg	
Overall dimensions (W × D × H)	235 × 235 × 140 mm	
Weight	1.2 kg	1.8 kg
Input current/power consumption	12 V, 260 mA/3.1 W	12 V, 380 mA/4.6 W
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V	

Accessories for the standard platform:

Optional dimpled mat PDM prevents different size tubes from rolling around the platform

PDM, dimpled mat



Multi Bio 3D with PDM mat



ORDERING INFORMATION:

3D with stand, platform **Bio PP-4S**

Multi Bio 3D with stand, platform **Bio PP-4S**

Optional accessories:

PDM dimpled mat

Cat. number 

BS-010151-AAG

BS-010125-AAG

PDM

PSU-10i, Orbital Shaker

DESCRIPTION

Shaker **PSU-10i** provides regulated orbital motion of the platform and is designed for use both in small specialized biotechnological laboratories and in large multidisciplinary laboratories: a choice of five (5) interchangeable platforms provides the possibility of performing various procedures and techniques.

Shaker **PSU-10i** incorporates a direct drive system, a brushless motor with a guaranteed service life of up to 35,000 hours and an automatic loading balancing system. These innovations allow for continuous mixing up to 7 days, ensure reliable, trouble-free operation for more than 2 years and significantly expand the device performance range in both high and low limits.

The unit is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C.

Premium
Product Class

Ø 10 mm
orbit



Product
video

PSU-20i, Orbital Shaker

DESCRIPTION

Shaker **PSU-20i** provides three motion types: **1 Orbital**, **2 Reciprocal** and **3 Vibrating**, which can be performed separately, pairwise and sequentially in repeated cycles.

Shaker is designed for applications both in small specialized laboratories and in large multidisciplinary laboratories. **PSU-20i** is an ideal instrument for laboratories researching biopharmaceutics and biomedicine.


Shaker **PSU-20i** is noiseless and reliable in operation, incorporates a direct drive system and brushless motor with a guaranteed service life up to 35,000 working hours. Direct drive and brushless motor allows for continuous mixing for up to 7 days and ensures reliable operation for more than two years.

A choice of nine different interchangeable platforms provides the possibility of performing various procedures and techniques. Special attention should be paid to a multilevel platform, which allows accommodation of a large number of various microplates, Petri dishes, cultural bags and other low containers.

The unit is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C.

Premium
Product Class

Ø 20 mm
orbit

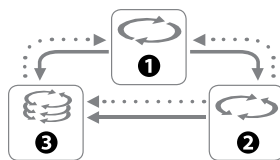

Multi-motion

Description and pictures of all platforms can be found on page 22-23

PSU-10i and PSU-20i, Orbital Shakers

SPECIFICATIONS

	PSU-10i	PSU-20i
Multi-motion	—	yes
Speed control range*	50–450* rpm (increment 10 rpm)	20–250* rpm (increment 5 rpm)
Digital speed control	yes	
Max. continuous operation time	168 h	
Orbit	10 mm	20 mm
Digital time setting	1 min–96 h/non-stop	
Timer sound signal	yes	
Maximum load	3 kg	8 kg
Overall dimensions (WxDxH)	255 × 255 × 100 mm	410 × 410 × 130 mm
Weight	3.4 kg	11.7 kg
Input current/power consumption	12 V, 800 mA/9.6 W	12 V, 3.2 A/40 W
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V	

* — max. speed depends on the load and vessels' shape

Platform P-6/250 for PSU-10i



Platform Bio PP-4 for PSU-10i



Platform for PSU-20i PP-20/4



Platform Bio PP-4 for PSU-10i



ORDERING INFORMATION:

Cat. number

PSU-10i, Shaker without platform BS-010144-AAN








PSU-20i, Shaker without platform BS-010145-ACI

PSU-20i motion types	Description	Speed range	Turning angle	Motion timer*	Digital time setting
1 Orbital	Orbital motion with an option of shifting direction	20–250 rpm	—	0–250 s	1 min–96 h (increment 1 min) or non-stop
2 Reciprocal	Orbital motion with shifting direction of rotation	20–250 rpm	0–360° (30° increment)	0–250 s	
3 Vibrating	High speed, low amplitude motion	—	0–5° (1° increment)	0–5 s	

* — for switching to the next motion in the cycle

Description and pictures of all platforms can be found on page 22-23

Platforms for PSU-10i and ES-20

Platform	Description	Dimensions (Working area)	Cat. number
UP-12 Used on PSU-10i, ES-20	 Universal platform with adjustable bars for different types of flasks, bottles and beakers with silicone mat	285 × 220 × 40 mm (270 × 185 × 40 mm)	BS-010108-AK
Bio PP-4 Used on PSU-10i	 Flat platform with silicone mat for Petri dishes, culture flasks, agglutination cards	255 × 255 mm (230 × 230 mm)	BS-010116-AK
PP-4 Used on ES-20, PSU-10i	 Metallic flat platform with silicone mat for Petri dishes, culture flasks, agglutination cards	220 × 220 mm (215 × 215 mm)	BS-010108-BK
NEW SPM	 Double-sided adhesive mat as an alternative for regular flask clamps (for PP-4)	220 × 220 mm	BS-010111-BK
P-12/100 Used on PSU-10i, ES-20	 Platform with clamps for flasks, 100–150 ml (12 places)	250 × 190 mm (250 × 190 mm)	BS-010108-EK
P-6/250 Used on PSU-10i, ES-20	 Platform with clamps for flasks, 250–300 ml (6 places)	250 × 190 mm (250 × 190 mm)	BS-010108-DK
P-16/88 Used on PSU-10i, ES-20	 Platform with spring holders for up to 88 tubes up to 30 mm diameter (e. g. 10 ml, 15 ml, 50 ml tubes)	275 × 205 × 75 mm (275 × 205 × 75 mm)	BS-010116-BK

SPML, SPM, Double-sided adhesive strips and mat NEW

DESCRIPTION

Convenient alternative to traditional steel holders, an easy way to fix tubes, plates, flasks and other laboratory wares on platforms for cultivation, incubation and mixing. Two size options are offered **SPML** can be used with UP-168 platform on Biosan orbital shaker PSU-20i and in ES 20/80, ES 20/60 orbital shakers.

SPM is compatible with PP-4 platform, which fits both on PSU-10i orbital shaker and in ES-20 Shaker-Incubator.

Made of polyurethane with adhesive, simple to clean and durable, able to withstand up to 1,000 times placement/removal or 12 months of use. Additional information about temperature, working volume and speed limitations available in the user manual.

SPECIFICATIONS

SPML Size (L×W×H)	390 × 80 × 3 mm (double sided 1.5 mm PU with PET adhesive)
SPM Size (L×W×H)	210 × 210 × 3 mm (double sided 1.5 mm PU with PET adhesive)
Colour	transparent
Duration of use	up to 1,000 times placement/removal or 12 months
Temperature range	+4C° to +80C°
Shaking speed	0–400 rpm

SPM on PP-4 platform



Video

SPML on UP-168 platform



ORDERING INFORMATION:

SPM, Double-sided adhesive mat


SPML, Set of 3 double-sided adhesive strips

Cat. number


BS-010111-BK

BS-010135-MK

Platforms for PSU-20i and ES-20/60

Platform	Description	Dimensions (Working area)	Cat. number
UP-330 Used on PSU-20i 	Universal platform with adjustable bars for different types of flasks, beakers	345 × 430 × 105 mm (300 × 400 × 80 mm)	BS-010145-AK
P-30/100 Used on PSU-20i, ES-20/60 	Platform with 30 clamps for 100–150 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010135-BK
P-16/250 Used on PSU-20i, ES-20/60 	Platform with 16 clamps for 250–300 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010135-CK
P-9/500 Used on PSU-20i, ES-20/60 	Platform with 9 clamps for 500 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010135-AK
P-6/1000 Used on PSU-20i, ES-20/60 	Platform with 6 clamps for 1,000 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010135-DK
PP-400 Used on PSU-20i, ES-20/60, ES-20/80 	Flat platform with non-slip silicone mat	360 × 400 mm (360 × 400 mm)	BS-010135-FK
UP-168 Used on ES-20/60, ES-20/80 	Universal platform for different flasks (Clamps ordered separately)	360 × 400 mm (360 × 400 mm)	BS-010135-JK
FC-50 FC-100 FC-250 FC-500 FC-1000 FC-2000 used on PSU-20i 	Clamp for 50, 100, 250, 500, 1000, 2,000 ml flask (for UP-168)	Ø 50 mm Ø 65 mm Ø 85 mm Ø 105 mm Ø 130 mm Ø 165 mm	BS-010126-MK BS-010126-HK BS-010126-JK BS-010126-LK BS-010126-IK BS-010126-NK
NEW SPML 	Set of 3 double-sided adhesive strips as an alternative for regular flask clamps (for UP-168)	390 × 80 × 3 mm	BS-010135-MK
TR-21/50 	Test tube rack for 50 ml with 21 drillings (for UP-168)	340 × 124 mm (2 per platform)	BS-010135-KK
TR-44/15 	Test tube rack for 15 ml with 44 drillings (for UP-168)	340 × 124 mm (2 per platform)	BS-010135-LK
PP-20/4 Used on PSU-20i 	Four-level flat platform with non-slip rubber mat	380 × 480 × 440 mm (355 × 455 mm × 4)	BS-010126-EK
PP-20/3 Used on PSU-20i 	Three-level flat platform with non-slip rubber mat	380 × 480 × 300 mm (355 × 455 mm × 3)	BS-010126-DK
PP-20/2 Used on PSU-20i 	Two-level flat platform with non-slip rubber mat	380 × 480 × 160 mm (355 × 455 mm × 2)	BS-010126-CK
PP-20 Used on PSU-20i 	One-level flat platform with non-slip rubber mat	380 × 480 mm (355 × 455 mm)	BS-010126-BK

MPS-1, High-Speed Multi Plate Shaker

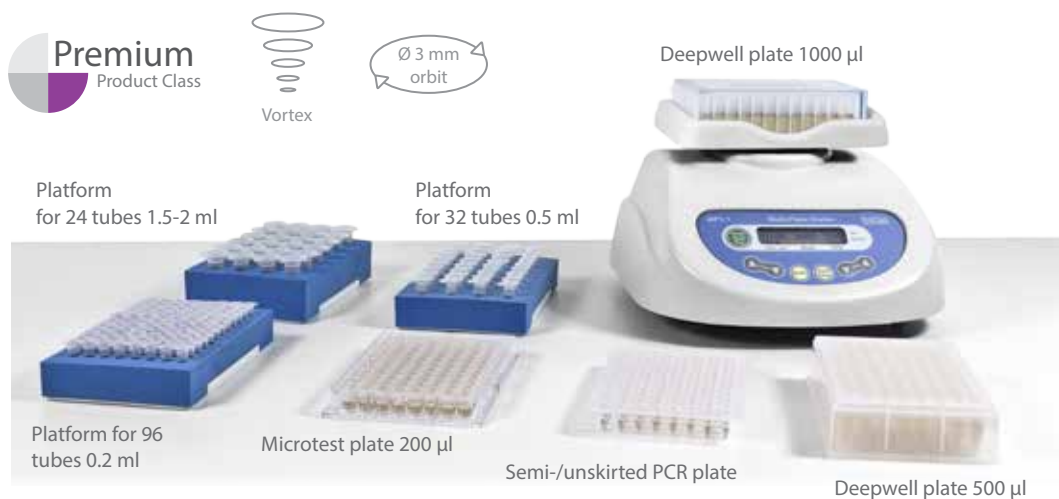**DESCRIPTION**

High-Speed Multi Plate Shaker **MPS-1** can be used in virtually any application by providing adjustable mixing of reagents in microtest plates, PCR plates, deepwell plates and test tubes (shaking tubes 0.2 to 2 ml and vortexing any volume up to 50 ml).

The shaker is compact and user-friendly. The shaker is ideal for personal use.

MPS-1 features a head for vortexing a single tube.

The unit is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C, decreasing linearly to 50% relative humidity at 40°C. Low voltage external power supply (12 V) provides electrical safety in a humid environment.

MPS-1 features **Pulse Mode** mixing function that works on the principle of giving a periodic impulse: the tube is accelerated to the set speed, holds it for 3 seconds and then drops the speed to zero. This motion is repeated until the timer runs out. This method provides a constant state particle resuspension inside a tube, as the acceleration is always changing. The advantage of this method is the high throughput of mixed samples compared to vortexing a single tube.



Product video

Features

- Speed control range 300–3,200 rpm
- Stable mixing with 3 mm orbit
- Five mixing presets
- **Pulse Mode** mixing function
- Quiet operation — low noise at maximum speed
- Universal platform holder for Deepwell plates and Microtest plates
- Additional four platforms for semiskirted and unskirted PCR plates 200 µl as well as for tubes from 0.2 to 2 ml

Platform for semi-/unskirted PCR plate 200 µl

Platform for 24 tubes 1.5–2 ml



MPS-1, High-Speed Multi Plate Shaker

Vortexing a 50 ml tube



Vortexing a 15 ml tube



Deepwell plate 96/1000 µl



Microtest plate 200 µl



Deepwell plate 96/500 µl



Mixing Speed control range	300–3,200 rpm
----------------------------	---------------

Platform options:

– For semi-/unskirted PCR plate or 96 microtest tubes 0.2 ml	P-02/96
– For 24 microtest tubes 1.5–2 ml	P-2/24
– For 32 microtest tubes 0.5 ml	P-05/32
– For 24 microtest tubes 0.5 ml and 48 microtest tubes 0.2 ml	P-02/05
– Universal platform for deepwell plates, 96-well microtest plates (U, V or flat bottomed), 384-well microtest plates	

Types of mixing presets:


VORTEX	3,200 rpm
HARD	2,600 rpm
MEDIUM	1,800 rpm
SOFT	1,000 rpm
CUSTOM	adjustable rpm

Features a **Pulse Mode** mixing function

Features a Vortex function

Maximum load	0.3 kg
Mixing Orbit	3 mm
Acceleration time to maximum speed	5 s
Digital time setting	0–60 min (15 s increment)/non-stop
Timer sound signal	yes
Maximum continuous operation time	8 h
Noise level, not more	65 dB
Weight	5.1 kg
Overall dimensions (WxDxH)	225 × 215 × 150 mm
Input current/power consumption	12 V, 800 mA / 10 W
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 12 V

ORDERING INFORMATION:

Cat. number 

MPS-1 , Multi Plate Shaker with built-in universal platform	BS-010216-A03
MPS-1 , Multi Plate Shaker with built-in universal platform and set of 4 platforms (P-02/96, P-2/24, P-05/32, P-02/05)	BS-010216-A11

Optional platforms:

Cat. number

1 P-02/96	For semi-/unskirted PCR plate or 96 microtest tubes 0.2 ml	BS-010216-CK
2 P-2/24	For 24 microtest tubes 1.5–2 ml	BS-010216-AK
3 P-05/32	For 32 microtest tubes 0.5 ml	BS-010216-BK
4 P-02/05	For 24 microtest tubes 0.5 ml and 48 microtest tubes 0.2 ml	BS-010216-DK

1 Platform P-02/96



2 Platform P-2/24



3 Platform P-05/32



4 Platform P-02/05



PSU-2T, Mini-Shaker

DESCRIPTION

Mini-Shaker **PSU-2T** is designed for immunoassays and provides adjustable mixing of reagents in microplates. The device ensures smooth movement of the platform even at low speeds.

Shaker is a compact and user-friendly device. It takes up little space on a desk and is ideal for personal use. Direct drive and brushless motor allow continuous mixing up to 7 days and ensures reliable, trouble-free operation for more than 2 years. Display of the device switches between time and speed readings.

The unit is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C.

SPECIFICATIONS

Speed control range	150–1,200 rpm
Digital time setting	1 min–24 h/non-stop
Digital setting and control of time and speed	
Max. continuous operation time	168 h
Direct drive mechanism	
Orbit	2 mm
Overall dimensions (WxDxH)	220 × 205 × 90 mm
Weight	2 kg
Input current/ power consumption	12 V, 280 mA/3.4 W
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V

ORDERING INFORMATION: Cat. number

PSU-2T with standard platform **IPP-2** BS-010155-AAG

Optional platforms

IPP-4 BS-010102-AK



Basic Plus
Product Class

Ø 2 mm orbit



Product video

A Platform **IPP-2**



B Platform **IPP-4**

**Platforms for microtest plates:**

A **IPP-2** (standard platform) 184 × 132 mm
for 2 microtest plates

B **IPP-4** (optional platform) 266 × 170 mm
for 4 microtest plates



Multi Bio RS-24 and Multi RS-60, rotators



Multi Bio RS-24



Product video



Multi RS-60




Product video

It is possible to choose the position of tubes for rocking motion – horizontal or vertical. The platform does not make an additional revolution before stopping in the horizontal plane.

Programmable Rotators performs several motion types in one module. Microprocessor control allows performing not only **1 Vertical overhead rotation** of the platform, but also **2 Reciprocal rotation (rocking motion)** as well as **3 Vibration**. These three motion types can be performed separately, pairwise and in cycles, periodically repeating the sequence of three motion types. Multi-Rotation option of Biosan instruments substantially expands possibilities and enhances the efficiency of sample preparation for the examined materials and allows adjusting the mixing procedure according to the individual tasks.

Programmable Rotators can be used for variety of applications in modern life science laboratories: for hybridisation reactions, cell growing, soft extraction and homogenisation of biological components in solutions, as well as for reactions of binding and washing of magnetic particles.

Multi Bio RS-24 and **Multi RS-60** are designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40 °C. Low voltage external power supply (12 V / 24V) provides electrical safety in a humid environment.

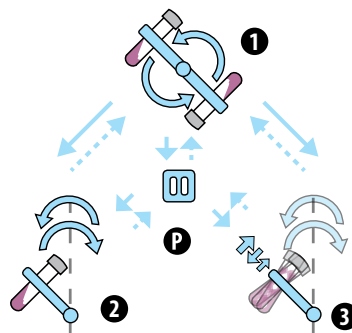

Programmable Rotator provides 3 rotation types and Pause:

- 1** Vertical overhead rotation
- 2** Reciprocal rotation (rocking motion)
- 3** Vibro
- P** Pause

Multi Bio RS-24 and Multi RS-60, rotator

SPECIFICATIONS

	Multi Bio RS-24	Multi RS-60
1 Vertical overhead rotation:		
Speed control range	1–100 rpm (increment 1 rpm)	
Vertical rotation movement	360°	
Time setting range	0–250 s	
2 Reciprocal rotation (rocking motion):		
Speed control range	1–100 rpm (increment 1 rpm)	
Tilt angle range	1–90° (increment 1°)	
Time setting range	0–250 s	
3 Vibro:		
Tilt angle range	0–5° (increment 1°)	
Pause/Vibro time setting range	0–5 s	
GENERAL SPECIFICATIONS:		
Digital time setting	1 min–24 h/non-stop (increment 1 min)	
Timer sound signal	yes	
Maximum load	0.5 kg	0.8 kg
Overall dimensions (W×D×H)	365 × 195 × 155 mm	430 × 230 × 230 mm
Weight	1.7 kg	3.8 kg
Input current/power consumption	12 V, 660 mA/8 W	24 V, 750 mA/18 W
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V	Input AC 100–240 V, 50/60 Hz; Output DC 24 V

Multi Bio RS-24 with optional platform PRSC-22



Multi RS-60 with standard platform PRS-48

**ORDERING INFORMATION:**

Cat. number

Multi Bio RS-24 with standard platform PRS-26

BS-010117-AAG

Multi RS-60 with standard platform PRS-48

BS-010118-AAI

Optional platforms for Multi Bio RS-24:

PRS-5/12

BS-010117-HK

PRS-10

BS-010117-IK

PRSC-22

BS-010117-LK

PRSC-10

BS-010117-JK

PRS-1DP

BS-010149-DK

M-8/50

BS-010117-PK

Optional platforms for Multi RS-60:

PRS-8/22

BS-010118-AK

PRS-14

BS-010118-BK

Description and pictures of all platforms can be found on page 29

Platforms for Multi Bio RS-24

Standard:	Capacity	Tube Volume	Tube Diameter	Cat. number
1 PRS-26	26	1.5–15 ml	10–16 mm	BS-010117-GK
Optional				
2 PRS-5/12	5 and 12	up to 50 and 1.5–15 ml	20–30 and 10–16 mm	BS-010117-HK
3 PRS-10	10	up to 50 ml	20–30 mm	BS-010117-IK
4 PRSC-22	22	15 ml	16 mm	BS-010117-LK
5 PRSC-10	10	50 ml	25–30 mm	BS-010117-JK
6 M-8/50	8	50 ml	25–30 mm	BS-010117-PK
7 PRS-1DP	Platform for microplates and racks for tall tubes 0.5 and 1 ml (e.g. Thermo 3741MTX, 3742MTX, 3744MTX)			BS-010149-DK

1 PRS-26

2 PRS-5/12

3 PRS-10

4 PRSC-22

5 PRSC-10

6 M-8/50

7 PRS-1DP


Clamps on PRSC-10



PRS series platforms are equipped with universal rubber clamps for different size tube fixation; PRSC series platforms have metal clamps able to hold heavier solutions (e.g. soil, sand).

Platforms for Multi RS-60

Standard:	Capacity	Tube Volume	Tube Diameter	Cat. number
1 PRS-48	48	1.5–15 ml	10–16 mm	BS-010118-CK
Optional:				
2 PRS-8/22	8 and 22	up to 50 and 1.5–15 ml	20–30 and 10–16 mm	BS-010118-AK
3 PRS-14	14	up to 50 ml	20–30 mm	BS-010118-BK

1 PRS-48

2 PRS-8/22

3 PRS-14


Bio RS-24, Mini-Rotator

DESCRIPTION

Mini-rotator **Bio RS-24** provides vertical rotation of the platform. The rotator is an ideal instrument for preventing blood coagulation in tubes and fulfilling of procedures of biological components extraction.

The device is simple to operate; it is designed as a low-cost solution.

The unit is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C. Low voltage external power supply (12 V) provides electrical safety in a humid environment.

SPECIFICATIONS

Speed control range	5–30 rpm
Vertical rotation movement	overhead, 360°
Digital time setting	1 min–24 h/non-stop (increment 1 min)
Timer sound signal	yes
Maximum continuous operation time	8 h
Overall dimensions (W×D×H)	325 × 190 × 155 mm
Weight	1.4 kg
Recommended load	75% of the rated volume
Input current/power consumption	12 V, 110 mA/1.3 W
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 12 V

PRS series platforms are equipped with universal rubber clamps for different size tube fixation;

PRSC series platforms have metal clamps able to hold heavier solutions (e.g. soil, sand).

Basic Plus
Product Class



Product video



Vertical rotation 360°

Bio RS-24 in operation



ORDERING INFORMATION:

Cat. number

Bio RS-24 with standard platform PRS-22	BS-010133-AAG
Optional platforms:	
PRS-4/12	BS-010117-AK
PRSC-18	BS-010117-EK

Platform	Capacity	Tube Volume	Tube Diameter, Ø
1 PRS-22 (standard)	22	1.5–15 ml	10–16 mm
2 PRS-4/12 (optional)	4 and 12	up to 50 and 1.5–15 ml	20–30 mm and 10–16 mm
3 PRSC-18 (optional)	18	15 ml	16 mm

1 PRS-22



2 PRS-4/12



3 PRSC-18



Basic Plus
Product Class



V-1 plus



Basic Plus
Product Class



V-32 with platform PV-32



V-1 plus and V-32, Vortexes

V-1 plus vortex and **V-32** multi vortex are intended for intensive mixing of samples in tubes with an eccentric mechanism.

Vortex can be used for different operations:

- Mixing tissue samples;
- Suspending cell samples;
- Mixing chemical samples;
- Mixing bacterial and yeast cells when washing from the culture medium;
- Extracting metabolites and enzymes from cells and cell cultures, etc.

Vortex can be used to perform various DNA/RNA operations, such as purification of low-molecular DNA/RNA fragments in PCR-diagnostics.

Vortex is applicable in all the fields of laboratory research in biotechnology, microbiology and medicine.

Vortexes have two operation modes:

- Continuous operation;
- Impulse operation. (**V1 plus** pressure activated)


Model **V-1 plus** is a personal vortex with a fluoroplastic head for single tube vortexing.

Model **V-32** is a universal vortex multipurpose device with different accessories. It is supplied with a 32-socket universal platform PV-32 for Eppendorf type tubes up to 1,5 ml (1.5/0.5/0.2 ml–16/8/8 sockets) and a PL-1 head for vortexing a single tube up to 50 ml. An optional 6-socket platform PV-6/10 for 10 ml tubes (maximum tube diameter 15 mm) or a platform PV-48 for six strips of eight 0.2 ml microtubes can be supplied on request.



Product video

Platform PL-1 for V-32

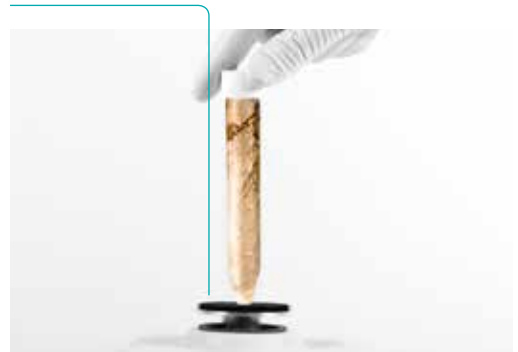

V-1 plus and V-32, Vortexes

SPECIFICATIONS

	V-1 plus	V-32
Mixing principle	Vibro Eccentric	
Speed control range	500–3,000 rpm	
Acceleration time	<1 s	3 s
Maximum continuous operation time	24 h	
Mixing module for tubes	from 0.2 to 50 ml	from 0.2 to 10 ml
Maximum mixing volume	30 ml	45 ml
Maximum load	30 g	70 g
Orbit	4 mm	2 mm
Dimensions (W×D×H)	90 × 150 × 80 mm	120 × 180 × 100 mm
Weight	0.8 kg	1.5 kg
Input current/power consumption	12 V, 320 mA/3.8 W	
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V	

V-1 Plus



V-1 Plus




V-1 Plus



Platform PV-6/10 for V-32



Platform PV-48 for V-32

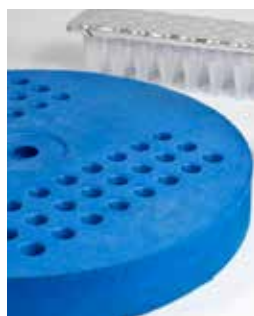

V-32 with platform PV-48



ORDERING INFORMATION:

V-1 plus

Cat. number

BS-010203-AAG

V-32 with standard platforms **PL-1** and **PV-32**

BS-010207-AAG

Optional platforms for V-32:

PV-6/10 platform for 6–10 ml tubes (max. Ø 15 mm)

BS-010207-BK

PV-48, platform for 6–8 × 0.2ml strips or 48 tubes of 0.2 ml

BS-010207-GK



MSV-3500 with platform SV-8/15



Product video

ORDERING INFORMATION:

Cat. number

MSV-3500 with all platforms	BS-010210-TAH
MSV-3500 without platform	BS-010210-AAH

MSV-3500, Multi Speed Vortex

Multi Speed Vortex **MSV-3500** is designed for soft or intensive mixing of reagents in different size and type plastic tubes (0.2 to 50 ml).

It is designed for operation in life science laboratories working in biochemistry, cell and molecular biology.

Unit has four types of interchangeable platforms: for Eppendorf type microtest tubes, 10/15/50 ml tubes (diameter 12/16/30 mm). Platforms can be ordered separately or as one set with **MSV-3500**.

Speed and time are under microprocessor control. LCD display indicates two lines of values: the set and actual values of speed and time.

Unit provides high maximum speed of platform rotation efficiently mixing microvolumes (less than 5 µl) of samples.

Speed control range	300–3,500* rpm
Digital time setting	0–60 min/non-stop (increment 1 min)
Timer sound signal	yes
Display	LCD, 2 × 16 signs
Orbit	4 mm
Maximum load	0.2 kg
Maximum continuous operation time	8 h
Dimensions (WxDxH)	180 × 170 × 145 mm
Weight	2.6 kg
Input current/power consumption	12 V, 1 A / 12 W
External power supply	Input AC 100–240 V, 50/60 Hz, Output DC 12 V

* — Maximum speed depends on load

DESCRIPTION

SPECIFICATIONS

Optional platforms:

		Cat. number
1 SV-16/8	Platform for 16 × 1.5 ml + 8 × 0.5 ml + 8 × 0.2 ml microtubes, Ø 11/8/6 mm	BS-010210-CK
2 SV-10/10	Platform for 10 × 10 ml tubes 12 mm diameter	BS-010210-BK
3 SV-8/15	Platform for 8 × 15 ml tubes 16 mm diameter	BS-010210-DK
4 SV-4/30	Platform for 4 × 50 ml tubes 30 mm diameter	BS-010210-AK

1 SV-16/8



2 SV-10/10



3 SV-8/15



4 SV-4/30



RCP-24, Homogenizer NEW

DESCRIPTION

Reciprocal Homogenizer **RCP-24**, a bench-top mechanical device designed for mixing, grinding, homogenizing and emulsifying biological objects in microtubes by vigorously mixing by reciprocal motion with various beads for sample preparation for subsequent academic, pharmaceutical, biotechnological or biomedical studies.

Homogenizer facilitates the formation of a supernatant containing nucleic acids and proteins suitable for subsequent purification, extraction or analysis. The device is optimized for extracting proteins, DNA, RNA or tRNA from various tissue sources, but it can also be used for other applications. **RCP-24** performs efficient homogenization of mammalian tissue, plant tissue or other biomaterials.

SPECIFICATIONS

Test tubes capacity	up to 24
Test tubes	2 ml
Speed control range	500–2,000 rpm (increment 100 rpm)
Digital time setting	1–15 min (increment 1 min)
Oscillation amplitude	44 mm, vertical
Dimension (WxDxH)	285 × 400 × 440 mm
Input current / power consumption	230 V, 50 Hz/ 220 W (1.3 A)
Weight	19.1 kg



Basic Plus
Product Class



Product video is available
on the website


ORDERING INFORMATION:
RCP-24, Homogenizer

Ceramic 1.4 mm Pre-Filled Bead Mill Tubes, sterile, set of 25 pcs.

Ceramic 2.8 mm Pre-Filled Bead Mill Tubes, sterile, set of 25 pcs.

Glass 0.1 mm Pre-Filled Bead Mill Tubes sterile, set of 25 pcs.

Glass 0.5 mm Pre-Filled Bead Mill Tubes sterile, set of 25 pcs.

Stainless steel 2.8 mm Pre-Filled Bead Mill Tubes, sterile, set of 25 pcs.

Cat. number

BS-010701-A02

BS-010701-AK

BS-010701-BK

BS-010701-EK

BS-010701-FK

BS-010701-KK

THERMO-SHAKERS



PST-60HL
Plate Shaker-Thermostat



TS-100C
Thermo-Shaker with cooling
for microtubes and PCR plates



TS-DW
Thermo-Shaker for Deep Well Plates

PST-60HL, PST-60HL-4 and PST-100HL, Thermo-Shakers

PST-60HL, PST-60HL-4 and **PST-100HL** Thermo-shakers are designed for shaking standard 96-well microtiter plates in the thermal regulation mode. Models **PST-60HL** and **PST-100HL** hold 2 plates, model **PST-60HL-4** has four plates.

A multisystem principle used in design of the Thermo-Shaker, allows operating it as three independent devices:

- Incubator;
- Microplate shaker;
- Thermo-Shaker.

A distinctive feature of Biosan Plate Thermo-Shakers is patented by the company **Two-Side Microplates Heating**, which allows achieving full correspondence of the set and actual temperature in the microplate wells.

Standard versions of Thermo-shakers provide heating up to 60°C, sufficient for carrying out ELISA tests.

Thermo-shaker **PST-100HL** with the ability to stabilize the temperature up to 100°C is specially designed for hybridisation reactions.

Plate Thermo-Shakers provide:

- Soft or intensive sample shaking;
- Rotation speed regulation, stabilization and indication
- Even rotation amplitude throughout the Thermo Automatic
- Setting and indication of the required temperature on the platform
- Automatic fault diagnostics (temperature sensor, platform heating, lid heating etc.)
- With the help of the temperature calibration function, the user can calibrate the unit to compensate for differences in the thermal behaviour of plates from different manufacturers; (**PST-60HL, PST-60HL-4**).

Application fields:

PST shakers can be used in various applications such as:

- **Immunochemistry** — Enzyme-Linked Immuno Sorbent Assay (ELISA). Unique bottom and top heating, while shaking, ensures the most efficient linkage of the target, thus providing the most reliable results;
- **Molecular biology** — Micro and Macro array applications — incubation with shaking provides more efficient hybridization of target nucleic acid with on the surface of Micro and Macro chip printed probes (**Specific holder is required**)



Ø 2 mm orbit

PST-60HL



Product video is available on the website



Ø 2 mm orbit

PST-60HL-4



Ø 2 mm orbit


PST-100HL



PST-60HL, PST-60HL-4 and PST-100HL, Thermo-Shakers

	PST-60HL	PST-60HL-4	PST-100HL
Temperature setting range	+25°C... +60°C		+25°C... +100°C
Temperature control range	+5°C above ambient... +60°C		+5°C above ambient... +100°C
Temperature setting resolution	0.1°C		
Temperature stability	±0.1°C		
Temperature uniformity @ +37°C	±0.25°C		±0.2°C
Temperature calibration coefficient range	0.936–1.063 (±0.063)		—
Heating	Two-side microplate heating (platform and lid)		Two-side microplate heating (platform and lid) + double heating contour of the platform
Orbit	2 mm		
Speed regulation range	250–1,200 rpm (increment 10 rpm)		
Digital time setting	1 min–96 h/non-stop (increment 1 min)		
Timer sound signal	yes		
Display	LCD, 2 × 16 signs		
Max. height of microtest plate	18 mm		
Number of microtest plates	2	4	2
Weight	6.1 kg	8.8 kg	5.9 kg
Platform dimensions (W×D)	250 × 150 mm	290 × 210 mm	250 × 150 mm
Overall dimensions (W×D×H)	270 × 260 × 125 mm	380 × 390 × 140 mm	270 × 260 × 125 mm
Input current/power consumption	12 V DC, 3.3 A/40 W	12 V DC, 4.15 A/50 W	12 V, 5 A/60 W
External power supply	Input AC 100–240 V 50/60 Hz, Output DC 12 V		



Heat up time PST-60HL and PST-60HL-4:

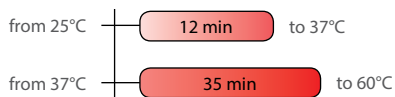

Heat up time PST-100HL:



PST-60HL-4 spring holders

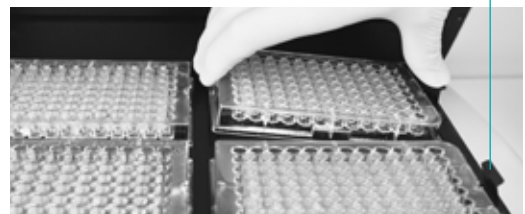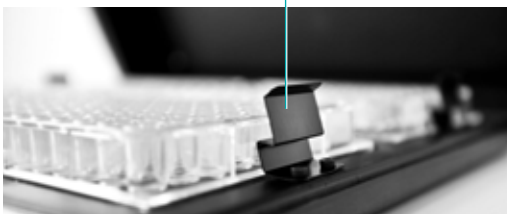

ORDERING INFORMATION:

PST-60HL

PST-60HL-4

PST-100HL

Cat. number

BS-010119-AAI

BS-010128-AAI

BS-010142-AAI

TS-100, TS-100C, and TS-100C Smart Thermo-Shakers

DESCRIPTION

TS-100 and **TS-100C** thermo-shakers are designed for intensive mixing of samples in microtest tubes or PCR plates in a temperature control environment. The **TS-100C** model of thermo-shaker differs from **TS-100** in the possibility of cooling samples down to +4°C.

Features of thermo-shakers meet the highest expectations of users according to many parameters:

- Fast reaching of specified mixing speed and maintenance of equal amplitude of rotation throughout the thermo-shaker block;
- Stability of maintaining the set temperature in a wide range throughout the block surface of thermo-shakers;
- With the help of the temperature calibration function, the user can calibrate the unit approximately $\pm 6\%$ of the selected temperature to compensate differences in the thermal behaviour of tubes from different manufacturers;
- LCD display indicates pre-set and current values of temperature, speed and time of operation;
- Quiet motor operation, compact size, prolonged service life.

Functions of heating and mixing can be performed either simultaneously or independently, which allows using the unit as three independent devices:

- **Thermostat;**
- **Shaker;**
- **Thermo-shaker.**

We offer five heating and cooling blocks for each model, including a block with a plastic lid for PCR-plates. Within one model of thermo-shaker, the blocks are mutually interchangeable and can be easily installed.



Mixing Efficiency
Video is available
on the website

The new model allows you to control the device in the following modes:

1. Manual using the front panel interface.
2. Through a computer program using Bluetooth® technology.

The software allows you to manage following parameters:

- Rotation speed
- Temperature
- Time
- Sound signal
- Creating Profiling programs using controlled parameters
- Visualization of temperature vs time and speed vs time graphs
- Data export to Excel and CSV formats
- Error messages/Fault diagnostics

Possibility of control up to seven units from PC. Independent parameter setting allows performing different tasks simultaneously on several units.

Premium
Product Class

Ø 2 mm
orbit

TS-100



TS-100
product video

Premium
Product Class

Ø 2 mm
orbit

TS-100C



TS-100C
product video

Smart Plus
Product Class

Ø 2 mm
orbit

TS-100C Smart




Bluetooth®
connection


TS-100, TS-100C, and TS-100C Smart Thermo-Shakers

	TS-100	TS-100C, TS-100C Smart
Temperature setting range	+25°C ... +100°C	+4°C ... +100°C
Temperature control range	5°C above ambient ... +100°C	15°C below ambient ... +100°C
Temperature setting resolution	0.1°C	
Temperature stability	±0.1°C	
Temperature accuracy @ +37°C	±0.5°C	
Average heating speed:	4°C/min from +25°C to +100°C	5°C/min from +25°C to +100°C
Average cooling speed:	—	from +100°C to +25°C 5°C/min from +25°C to +4°C 1.8°C/min
Temperature uniformity over the block:	@ +37°C ±0.1°C @ +60°C ±0.2°C @ +100°C ±0.2°C	@ +4°C ±0.6°C @ +37°C ±0.1°C @ +100°C ±0.3°C
Temperature calibration coefficient range	0.936–1.063 (±0.063)	
Speed control range	250–1,400 rpm	
Acceleration time	3 s	
Orbit	2 mm	
Display	LCD, 2 × 16 signs	
Microprocessor controlled temperature, mixing speed and operation time		
Digital time setting	1 min–96 h (1 min increment)	
Timer sound signal	yes	
Maximum continuous operation time	96 h	
Overall dimensions (W×D×H)	220 × 240 × 130 mm	
Weight	3.7 kg	
Input current/power consumption	12 V, 3.5 A/42 W	12 V, 4.9 A/60 W
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V	
PC software	—	only for TS-100C Smart

Heat up times for TS-100:

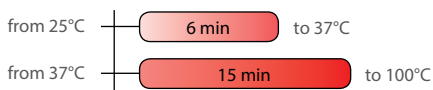

Heat up and cool down times for TS-100C and TS-100C Smart:

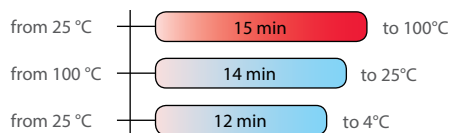

Thermo-Shakers are capable to support various application such as:

- **Molecular diagnostics** — Sample lysis for further Nucleic acid automated or manual extraction;
- **Genetic** — Amplicon denaturation for NGS Library preparation;
- **Biochemistry** — Enzymatic reaction;
- **Genomics** — Protein degradation studies;
- **Cellular biology** — Extraction of metabolites from cellular material.

ORDERING INFORMATION:

	Cat. number 
TS-100 without block	BS-010120-AAI
TS-100C without block	BS-010143-AAI
TS-100C Smart with software, without block	BS-010171-A01

Photos and descriptions of all blocks can be found on next page

Interchangeable Blocks for TS-100

Optional Blocks:		Tube's volume	Cat. number	
1	SC-18	20 and 12 microtubes	0.5 ml and 1.5 ml	BS-010120-AK
2	SC-18/02	20 and 12 microtubes	0.2 ml and 1.5 ml	BS-010120-CK
3	SC-24	24 microtubes	2 ml	BS-010120-EK
4	SC-24N	24 microtubes	1.5 ml	BS-010120-GK
5	SC-96A	96-well unskirted or semi-skirted microplate (0.2 ml) for PCR or 12 × 8–0.2ml strips or 96 tubes of 0.2 ml.		BS-010120-FK

1 SC-18



2 SC-18/02



3 SC-24



4 SC-24N



5 SC-96A



Interchangeable Blocks for TS-100C and TS-100C Smart

Optional Blocks:		Tube's volume	Cat. number	
1	SC-18C	20 and 12 microtubes	0.5 ml and 1.5 ml	BS-010143-AK
2	SC-18/02C	20 and 12 microtubes	0.2 ml and 1.5 ml	BS-010143-CK
3	SC-24C	24 microtubes	2 ml	BS-010143-EK
4	SC-24NC	24 microtubes	1.5 ml	BS-010143-GK
5	SC-96AC	96-well unskirted or semi-skirted microplate (0.2 ml) for PCR or 12 × 8–0.2ml strips or 96 tubes of 0.2 ml.		BS-010143-FK

1 SC-18C



2 SC-18/02C



3 SC-24C



4 SC-24NC



5 SC-96AC



TS-DW, Thermo-Shaker for deep well plates



Product
video

Deep Well Plate Thermo-Shaker provides:

- Soft or intensive sample shaking;
- Rotation speed regulation, stabilisation and indication;
- Even rotation amplitude throughout the Thermo-Shaker platform;
- Exceptional temperature uniformity across the plate;
- Required operation time setting and indication;
- Automatic stopping of the platform movement after expiration of the set time;
- Setting and indication of the required temperature on the platform;
- A variety of changeable blocks that can accommodate most popular deepwell plates;
- Automatic fault diagnostics (temperature sensor, platform heating, lid heating etc.).

Separate blocks to accommodate different deepwell plates will be released. For example:

- Deep Well Plates NUNC® 96/2,000 µl
- Deep Well Plates Eppendorf® 96/0.5 ml

Application fields:


- **Cytochemistry** — for *in situ* reactions;
- **Immunochemistry** — for immunofluorescent reactions;
- **Biochemistry** — for enzyme and protein analysis;
- **Molecular biology** — for nucleic acid extraction.

TS-DW Thermo-Shaker is designed for shaking and incubating deep well plates.

A multisystem principle, used in the Thermo-Shaker design, allows operating it as three independent devices: Incubator, Plate shaker and Thermo-Shaker.

TS-DW provides excellent temperature uniformity across the plate due to patented two-sided heating of the block and the lid, contour heating of the block and close proximity of heating elements to plate walls.

There is a number of interchangeable blocks to suit different plates such as Eppendorf® 96/1,000 µl, Sarstedt® Megablock 96/2,200 µl, Porvair® 96/2,000 µl, Axygen® 96/2,200 µl. Also, we can manufacture a customized block on request.

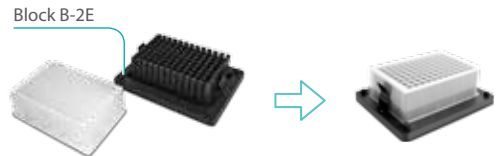

The block for deepwell plate is mountable, thus a custom plate module can be manufactured on request

Temperature setting range	+25 °C ... +100 °C
Temperature control range	5 °C above ambient ... +100 °C
Temperature setting resolution	0.1 °C
Temperature uniformity @ +37 °C	±0.1 °C*
Temperature accuracy @ +37 °C	±0.5 °C*
Temperature calibration coefficient range	0.936–1.063 (±0.063)
Time of platform heating from +25 °C to +37 °C	6 min*
Speed control range	250–1,400 rpm
Orbit	2 mm
Display	LCD, 2 × 16 signs
Digital time setting	1 min–96 h (1 min increment)
Timer sound signal	yes
Overall dimensions (W×D×H)	240 × 260 × 160 mm
Weight	5.1 kg
Input current/power consumption	12 V, 4.8 A/58 W
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 12 V

* — For B-2E block

ORDERING INFORMATION:

TS-DW without block

Cat. number

BS-010159-A02

Interchangeable blocks description and catalogue numbers can be found on next page

TS-DW, Thermo-Shaker for deep well plates

ORDERING INFORMATION

 Cat. number 

Interchangeable Blocks:		Volume	D-w plate cat. no.	
AP-96	Block for one deep-well plate Eppendorf®	96/500 µl	0030501101	BS-010159-QK
B-05PO	Block for one deep-well plate PlateOne®	96/0.5 ml	S1896-5000	BS-010159-QK
B-06A	Block for one deep-well plate Corning Axygen®	96/600 µl	P-DW-500C	BS-010159-KK
B-08AB	Block for one deep-well Abgene™ Storage Plate	96/800 µl	AB0765, AB0859	BS-010159-MK
B-12AB	Block for one deep-well plate Abgene™	96/1.2 ml	AB0564, AB0787	BS-010159-SK
B-1R	Block for one deep-well plate Riplate®	96/1,000 µl	43001-0101	BS-010159-UK
B-2A	Block for one deep-well plate Axygen®	96/2,200 µl	P-2ML-SQ-C	BS-010159-FK
B-2BBI	Block for one deep-well strips with rack SSIbio®	96/1.2 ml	703B00, 713B00	BS-010159-JK
B-2KF	Block for one PowerMag® Glass Bead 96 Plate	96/-	27600-4-KF-BP	BS-010159-LK
B-2M	Block for one KingFisher™ deep-well 96 Plate	96/-	95040450	BS-010159-RK
B-2N	Block for one deep-well plate Nunc®	96/2,000 µl	278743, 278752	BS-010159-DK
B-2PO	Block for one deep-well plate PlateOne®	96/2 ml	S1896-2000	BS-010159-NK
B-2R	Block for one deep-well plate Riplate®	96/2,000 µl	43001-0103	BS-010159-TK
B-2S	Block for one deep-well plate Sarstedt® Megablock	96/2,200 µl	82.1972.002	BS-010159-CK
B-2SD	Block for one deep-well plate Slicprep™ 96 Device	96/2.2 ml	AB0932	BS-010159-PK
B-2SL	Block for one deep-well plate Starlab®	96/1,200 µl	E2896-0120	BS-010159-IK
B-2E	Block for one deep-well plate Eppendorf®	96/1,000 µl	0030505204, 0030506200, 0030502205	BS-010159-AK
B-2P	Block for one deep-well plate Porvair®, Nunc®, Biotage®	96/2,000 µl	219009, 95040452, 121-5203	BS-010159-EK


Block B-2S



Block B-2N



Block B-2E

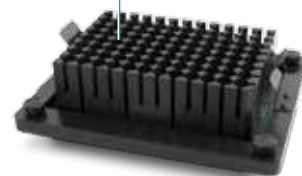

MINICENTRIFUGES-VORTEXES, MINI-CENTRIFUGE, CENTRIFUGES



FVL-2400N
Mini-Centrifuge/Vortex



MSC-6000
Centrifuge/Vortex Multispin



CVP-2
Centrifuge/Vortex for PCR plates

FV-2400, Microspin and FVL-2400N, Combi-Spin


DESCRIPTION

Minicentrifuges-Vortexes Microspin **FV-2400** and Combi-Spin **FVL-2400N** is specially designed for genetic engineering research (for PCR-diagnostics experiments). Units can be used in biomedical and biotechnological laboratories.

Minicentrifuges-Vortexes provide simultaneous mixing and separation of 12 samples, using centrifuge and mixing modules located on the common spin-module. Sequential combination of these operations allows you to collect all material at the bottom of the tube.

FV-2400 is an "open type" centrifuge (without lid), that increases the speed of centrifugation and resuspension operations.

FVL-2400N has a bioform design and equipped with a transparent protective lid accompanied by a protection mechanism that stops the rotor motion when the lid is opened.

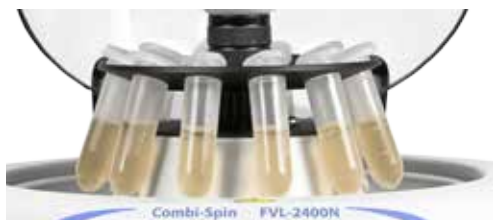

Rotor R-1.5



SPECIFICATIONS

	FV-2400	FVL-2400N	FV-2400	FVL-2400N
Rotation speed (fixed)	2,800 rpm		3,500 rpm	
Max. RCF	500×g		700×g	
Continuous and impulse operation modes				
Safety		Stop at open lid		Stop at open lid
Overall dimensions (W×D×H)	120×170×120 mm	190×235×125 mm	120×170×120 mm	190×235×125 mm
Weight	1.4 kg	1.7 kg	1.4 kg	1.7 kg
Nominal operating voltage	230 V, 50 Hz	230 V, 50 Hz	120 V, 60 Hz	120 V, 60 Hz
Power consumption (120 / 230 V)	30 W (0.13 A)		30 W (0.27 A)	

Rotors for FV-2400 and FVL-2400N

ORDERING INFORMATION:

Cat. number 

FV-2400 white with standard rotor R-1.5M and R-0.5/0.2M

BS-010201-AAA

FVL-2400N with standard rotors R-1.5 and R-0.5/0.2

BS-010202-AAA


Optional rotors: **see table below**

Rotors for FV-2400:		Capacity	Type	Cat. number
1 R-0.5/0.2M	12 × 0.5 ml and 12 × 0.2 ml microtubes	24	Standard	BS-010201-BK
2 R-1.5M	12 × 1.5/2 ml microtubes	12	Standard	BS-010201-AK
3 R-2/0.5	8 × 1.5/2 ml and 8 × 0.5 ml microtubes	16	Optional	BS-010205-CK
4 R-2/0.5/0.2	6 × 1.5/2 ml, 6 × 0.5 ml and 6 × 0.2 ml microtubes	18	Optional	BS-010205-DK
5 SR-16	Two 8-section strips for 0.2 ml microtubes	16	Optional	BS-010202-AK
6 SR-64*	Eight 8-section strips for 0.2 ml microtubes	64	Optional	BS-010201-EK

* — For any type of strips including paired



Tube vortexing on FV-2400

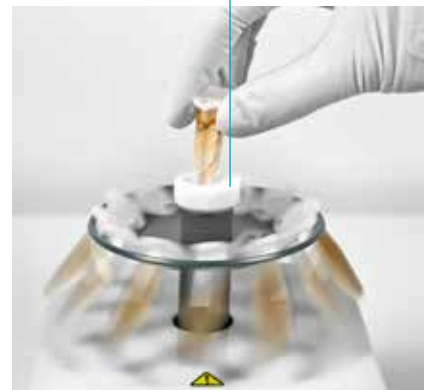

Rotors for FVL-2400N:		Capacity	Type	Cat. number
1 R-0.5/0.2	12 × 0.5 ml and 12 × 0.2 ml microtubes	24	Standard	BS-010205-BK
2 R-1.5	12 × 1.5/2 ml microtubes	12	Standard	BS-010205-AK
3 R-2/0.5	8 × 1.5/2 ml and 8 × 0.5 ml microtubes	16	Optional	BS-010205-CK
4 R-2/0.5/0.2	6 × 1.5/2 ml, 6 × 0.5 ml and 6 × 0.2 ml microtubes	18	Optional	BS-010205-DK
5 SR-16	Two 8-section strips for 0.2 ml microtubes	16	Optional	BS-010202-AK
6 SR-32*	Four 8-section strips for 0.2 ml microtubes	32	Optional	BS-010205-FK

* — Not compatible with Combi-Spins produced before 2015



MSC-3000 and MSC-6000, Multi-Spins

Centrifuge/vortex Multi-Spins **MSC-3000** and **MSC-6000** are products of extensively evolving Spin-mix-Spin technology that is intended for collecting micro volumes of reagents on the microtube's bottom (first centrifugation spin), following mixing (mix) and collecting the reagents again from the walls and cap of the microtube (second spin). Aim of this repetitive algorithm of operation is to reduce the mistakes during sample preparation for PCR analysis. We named it "sms-algorithm".

Multi-Spin is a fully automatic device for reproducing sms-algorithm for 12 tubes at one time, thus saving time considerably. A must-have instrument for PCR and DNA analyses laboratory.

Multi Spin is four devices combined in one:

1. Centrifuge — Maximum RCF:
 - MSC-3000:** up to $800 \times g$
 - MSC-6000:** up to $2,350 \times g$
2. Vortex (3 mixing modes — **soft, medium, hard**; regulated time; Vortexing regulation timer 1–20 s)
3. Centrifuge/Vortex;
4. SMS-cycler for realisation of the "sms-algorithm".



MSC-3000 product video



MSC-6000 product video

Saving time with multi-spin

Multi-Spin allows considerable time saving compared to Combi-Spin by automatically performing a cycling program of sample mixing and spinning according to the set spin-mix-spin cycle for 12 microtubes simultaneously.



	FVL-2400N	MSC-3000	MSC-6000
Speed control max.	2,800 rpm	3,500 rpm	6,000 rpm
RCF max.	$500 \times g$	$800 \times g$	$2,350 \times g$
Number of tubes vortexing	1 individually	12 simultaneously	
Time for completing "spin-mix-spin" cycle:			
for 2 microtubes	60 s	25 s	15 s
for 12 microtubes	5–6 min	90 s	60 s
for 100 microtubes	60 min	15 min	10 min
Unit price ratio	1 ×	1.5 ×	1.6 ×

MSC-3000 and MSC-6000, Multi-Spins

	MSC-3000	MSC-6000
Speed regulation range (increment 100 rpm)	1,000–3,500 rpm	1,000–6,000 rpm
RCF max.	800 × g	2,350 × g
Spin timer	1 s–99 min	1 s–30 min
Timer sound signal	yes	
Vortexing intensity	Soft, medium, hard	
Vortexing time	0–20 s (increment 1 s)	
SMS-cycle regulation	1–999 cycles	
Display	LCD, 2 × 16 signs	
Safety	Autostop at open lid	Lid lock
Overall dimensions (W×D×H)	190 × 235 × 125 mm	
Weight	2.1 kg	2.5 kg
Input current/power consumption	12 V, 11 W (0.9 A)	24 V, 24 W (1 A)
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 12 V	Input AC 100–240 V 50/60 Hz; Output DC 24 V

ORDERING INFORMATION:

Cat. number

MSC-3000 with standard rotors R-1.5, R-0.5/0.2


BS-010205-AAN

MSC-6000 with standard rotors R-1.5, R-0.5/0.2

BS-010211-AAL

Optional rotors: see table below

Rotor R-1.5

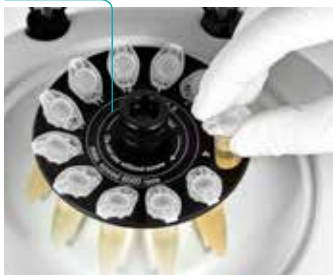

MSC-3000



MSC-6000

Optional rotors:	Capacity	Type	Cat. number	
1 R-0.5/0.2	12 × 0.5 ml and 12 × 0.2 ml microtubes	24	Standard	BS-010205-BK
2 R-1.5	12 × 1.5/2 ml microtubes	12	Standard	BS-010205-AK
3 R-2/0.5	8 × 1.5/2 ml and 8 × 0.5 ml microtubes	16	Optional	BS-010205-CK
4 R-2/0.5/0.2	6 × 1.5/2 ml, 6 × 0.5 ml and 6 × 0.2 ml microtubes	18	Optional	BS-010205-DK
5 SR-16	Two 8-section strips for 0.2 ml microtubes	16	Optional	BS-010202-AK
6 SR-32*	Four 8-section strips for 0.2 ml microtubes	32	Optional	BS-010205-FK

* — Not compatible with Multi-Spins produced before 2015

1 R-0.5/0.2



2 R-1.5



3 R-2/0.5



4 R-2/0.5/0.2



5 SR-16



6 SR-32



CVP-2, Centrifuge vortex for PCR plates

DESCRIPTION

After many years of Combined Centrifuge/Vortex concept success, we are proud to introduce the long-awaited Centrifuge vortex for PCR plates, **CVP-2**, to the sample preparation market. The Spin–Mix–Spin technology is intended to spin-down micro volumes of reagents on the well's bottom (first centrifugation spin), following mixing (mix) and spin-down the reagents again from the walls and cap of the well (second spin). Aim of this repetitive algorithm of operation is to reduce the mistakes during sample preparation for PCR analysis.

CVP-2 is a fully automatic device for reproducing sms–algorithm for 2 PCR plates at the same time, thus saving time considerably. A must-have instrument for PCR and DNA analyses laboratory.

CVP-2 is 4 devices combined in 1:

1. Centrifuge — Maximum RCF: $245 \times g$ (1,500 rpm)
2. Vortex (300–1,200 rpm; Vortexing regulation timer 0–60 sec)
3. Centrifuge vortex
4. SMS–cyclers for realization of the “sms–algorithm”

Tested plate types for use with CVP-2 centrifuge:

- Full-skirted 96-well standard micro-plates (without adapter)
- Half-skirted 96-well standard micro-plates (with adapter AP-96)
- Unskirted 96-wel standard I micro-plates (with adapter AP-96)
- Applied Biosystems™ MicroAmp™ Optical 96-well reaction plate (with adapter AP-96)
- Applied Biosystems™ MicroAmp™ Optical 384-well reaction plate (with adapter AP-384)
- For specific plate usage, please contact us for evaluation.

SPECIFICATIONS

Speed regulation range	300–1,500 rpm
Min. RCF at 1,500 rpm	175 × g
Vortex regulation range	300–1,200 rpm
Setting resolution	100 rpm
Plate type:	
• Without adapter:	96-well skirted PCR plates, PCR strips in a frame;
• With adapter AP-96 :	96-well semi-skirted and non-skirted PCR plates;
• With adapter AP-384 :	384-well PCR plates;
Display	LCD, 2 × 16 signs
Centrifugation mode time range	0–30 min
Centrifugation mode time increment	1 s; after 1 min–1 min
Vortex mode time range	0–60 s
Timer sound signal	yes
Number of programmable cycles	1–999
Chamber diameter	210 mm
Overall dimensions (W×D×H)	285 × 350 × 190 mm
Weight	6.15 kg
Input current/power consumption	12V, 1.5 A/18 W
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 12 V

Premium
Product Class



Product
video



Adapter AP-96 for 96-well semi-skirted and unskirted PCR plates



Adapter AP-384 for 384-well PCR plates



ORDERING INFORMATION

Cat. number



CVP-2

BS-010219-A02

With rotor for two PCR plates, protection lid and adapters AP-96* (a set of 2 adapters for 96-well semi-skirted and unskirted PCR plates)

Optional accessories:

AP-384*

BS-010219-EK

A set of 2 adapters for 384-well PCR plates

* — Adapters are made of Ertacetal® C and are autoclavable

High-speed Mini-centrifuge **Microspin 12**

DESCRIPTION

SPECIFICATIONS

Basic Plus
Product Class



Product video



Protection lid

1 A-02 Adapters



2 A-05 Adapters



High-speed Mini-centrifuge **Microspin 12** is a compact desktop centrifuge designed for biomedical laboratories.

Microspin 12 is used to extract RNA/DNA samples, sedimentation of biological components, biochemical and chemical analysis of microsamples.

A display simultaneously shows actual and set values for:

1. Centrifugation time;
2. Set and actual speed values;
3. Relative centrifugal force.

A brushless motor provides noiseless performance at maximal speed and long service life. An angular rotor is designed to accommodate 12 Eppendorf microtubes and spin columns (autoclavable adapters for 0.2, 0.5 ml tubes included). The rotor is made of aluminium, it is equipped with a fixing lid and included in the standard specification of the centrifuge. Constant airflow around the rotor reduces the risk of samples overheating during operation.

Metal protective inserts inside the casing and lid, automatic imbalance switch-off and locking of a lid provide safe operation. A sound signal indicates the completion of centrifugation.

The external power supply unit allows operating **Microspin 12** in cold rooms (at ambient temperatures from +4°C to +40°C).

Speed control range	1,000–14,500 rpm (100 rpm increment)
Relative centrifugal force control range	50–12,400 × g
Digital time setting	15 s–30 min
Timer sound signal	yes
Time setting resolution	1 min–15 s; after 1 min–1 min
Acceleration time up to 14,500 rpm	20 s
Slowdown time, not more	10 s
Display	LCD, 2 line
Safety: Rotor imbalance diagnostics: automatic stop, "IMBALANCE" warning	
Overall dimensions (W×D×H)	200 × 240 × 125 mm
Weight	3.5 kg
Input current/power consumption	24 V, 2.5 A / 60 W
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 24 V

ORDERING INFORMATION: Cat. number

Microspin 12 BS-010213-AA1
Built-in rotor MSR-12 (12 places for microtubes 1.5/2 ml)
with protection lid MSL-SC and adapters A-02, A-05 (autoclavable)

Additional/replacement parts:

- | | |
|--|--------------|
| MSL-SC, protection lid for rotors | BS-010213-EK |
| 1 A-02, 12 pieces for microtubes 0.2 ml | BS-010213-BK |
| 2 A-05, 12 pieces for microtubes 0.5 ml | BS-010213-AK |

LMC-3000, Laboratory Centrifuge

DESCRIPTION

LMC-3000 is a modern low-speed bench-top centrifuge designed for operation with microtest plates and centrifuge tubes up to 50 ml, Gel Cards. This device is widely used in biomedical profile laboratories.

Features:

- Soft start and run-down of the rotor;
- User-friendly setting of centrifugation parameters and simultaneous display of both set and actual values;
- Safe operation at any speed is provided by metal protection chamber and case cover, automatic stop at imbalance and a lock keeping the lid closed while the centrifuge is running;
- Low noise level;
- Rotor selection;
- Setting rotor speed in RPM or RCF (Relative Centrifugal Force);
- Multiple accelerations (Slow, Normal, Fast) and deceleration (0, Slow, Normal, Fast) modes and possibility to switch off forced braking;
- Wide choice of accessory rotors.

SPECIFICATIONS

Speed regulation range for centrifuge tubes	100–3,000 rpm (1,610 × g)
Speed regulation range for microtitre plates	100–2,000 rpm (560 × g)
Setting resolution	100 rpm

Rotor imbalance diagnostics
(automatic stop, "IMBALANCE" warning)

Display	LCD, 2 × 16 signs
Digital time setting	1–90 min (increment 1 min)
Timer sound signal	yes
Chamber diameter	340 mm
Overall dimensions (W×D×H)	420×495 × 235 mm
Weight	11.8 kg


Nominal operating voltage	230 V, 50/60 Hz or 120 V, 50/60 Hz
Power consumption (230/120 V)	110 W (0.5 A)/ 120 W (1 A)



ORDERING INFORMATION: Cat. number

LMC-3000 without rotors BS-010208-AAA

Basic Plus
Product Class

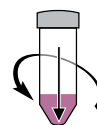

Product
video



Rotor R-12/15



Rotors description, pictures and catalogue numbers can be found on page 52–53

LMC-4200R, Laboratory Refrigerated Centrifuge



Product video



DESCRIPTION

Laboratory bench-top centrifuge with refrigeration **LMC-4200R** provides temperature control of biomaterial during centrifugation. Temperature control of the so-called "cold-shelf" is a gold standard for enzymologists and cell biologists because it ensures conditions necessary for reproducibility of the sample preparation stage. Temperature control absence at this stage can cause unpredictable results.

LMC-4200R is a modern centrifuge designed for operation with microtest plates, Gel Cards and tubes from 2 to 50 ml.

SPECIFICATIONS

Temperature control range	-10°C ... +25°C
Stable temperature maintenance range	25°C below ambient ... to +25°C
Temperature setting resolution	1°C
Speed regulation range for centrifuge tubes	100–4,200 rpm (3,160 × g)
Speed regulation range for microtitre plates	100–2,000 rpm (560 × g)
Speed setting resolution	100 rpm
Rotor imbalance diagnostics (automatic stop, "IMBALANCE" warning)	
Slowdown time, not more	30 s
Display	LCD, 2 lines
Digital time setting	1–90 min (increment 1 min)
Timer sound signal	yes
Chamber diameter	360 mm
Dimensions (WxDxH)	635 × 580 × 335 mm
Weight	56 kg
Nominal operating voltage	230 V, 50 Hz
Power consumption (230 V)	990 W (4.3 A)

Features:

- Effective way of acceleration and deceleration: Run-up time 20 sec; Run-down time, not more 30 sec;
- Efficient rate of chamber refrigeration: under 10 min;
- Maintenance of stable temperature during operation;
- User-friendly setting of centrifugation parameters (speed, temperature, time) and simultaneous display of both set and actual values;
- Safe operation is provided by a metal protection chamber and a case cover, automatic stop at imbalance (emergency shutdown, "IMBALANCE" displayed) and a lock keeping the lid closed while the centrifuge is running;
- Low noise level;
- Possibility to switch off forced braking;
- Wide choice of accessory rotors;
- Rotor selection;
- Setting rotor speed in RPM or RCF (Relative Centrifugal Force);
- Multiple accelerations (Slow, Normal, Fast) and deceleration (0, Slow, Normal, Fast) modes and possibility to switch off forced braking;

Rotor R-24/10



ORDERING INFORMATION:

Cat. number

LMC-4200R without rotors

BS-010212-AAA

Rotors description and pictures can be found on next pages

Interchangeable Rotors for LMC-3000 and LMC-4200R

Rack RR-U



	Rotor R-12/10	Rotor R-24/10	Rotor R-6	Rotor R-6P
Rotor type	Swing-out			
Dimensions	Ø16 × 105 mm		Ø29 × 115 mm	
Capacity	12	24	6	
Tube's volume	10–15 ml		50 ml	
Max. speed	4,200 rpm	4,000 rpm	4,200 rpm	
Max. RCF:	LMC-3000 1,610 × g	<i>Not applicable</i>	1,610 × g	
	LMC-4200R 3,160 × g	2,860 × g	3,160 × g	
Cat. number:	BS-010208-BK	BS-010212-JK	BS-010208-DK	BS-010208-XK



Plastic conical bottom centrifuge tube

Manufacturers: **Falcon, Greiner Bio-one, Sarstead, Corning, Nunc, TPP, etc.**

HOW TO CHOOSE ROTOR



15 ml



R-12/15



50 ml




R-6



R-6P


Plastic round bottom centrifuge tube, Vacutainers

Manufacturers: **Nunc, Greiner, Greiner Bio-one, TPP, etc.**

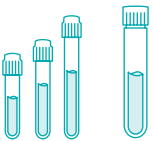

2–9 ml

10–15 ml

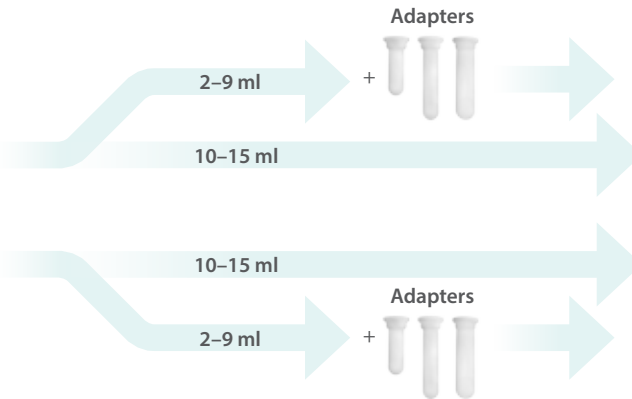

Adapters



R-12/10



R-24/10

Only for LMC-4200R


Interchangeable Rotors for LMC-3000 and LMC-4200R



Rotor R-2	Rotor R-24GC	Rotor R-12/15	Rotor RMT-24
Swing-out		Angled Swing-out	Fixed-angle
128 × 85.6 mm	53 × 74 mm	Ø17 × 120 mm	Ø11.1
2	24	12	24
up to 45 mm	—	15 ml	1.5–2 ml
2,000 rpm	1,500 rpm	4,200 rpm	6,000 rpm
560 × g	280 × g	1,610 × g	<i>Not applicable</i>
560 × g	280 × g	3,160 × g	1,837 × g
BS-010208-AK	BS-010208-VK	BS-010208-EK	BS-010221-BK

Standard 96-well microtitre plates, skirted PCR plates and deepwell plates up to 45 mm

Manufacturers: **Nunc, Greiner, Greiner Bio-one, etc.**

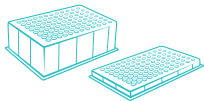

R-2


NEW

1.5–2 ml

RMT-24

96-well semi-/ unskirted PCR plate

Manufacturers: **Nunc, Greiner, Greiner Bio-one, etc.**

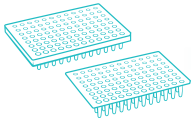

Adapter

AP-96


Material: **Ertacetal® POM-C** and is autoclavable



R-2

Gel Cards

Manufacturers: **Grifols®, DiaMed®, Bio-Rad® etc.**

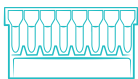

R-24GC, Rotor for Gel Cards for blood group serology testing (Forward Group, Reverse Group, RhD Type and 3 cell antibody screen). Recommended centrifugation time – 9 minutes



R-24GC

ORDERING INFORMATION: optional accessories for rotors

Cat. number

Adapter* for R-2:

AP-96	2 adapters for 96-well semi-skirted and non-skirted PCR plates	BS-010219-DK
Adapters** for R-12/10, R-24/10:	Vacutainers dimensions (Ø × length)	
BN-13/75	for vacutainers® 2–5 ml 13 × 80 mm	BS-010208-PK
BN-13/100	for vacutainers® 4–8 ml 13 × 105 mm	BS-010208-QK
BN-16/100	for vacutainers® 8–10 ml 16 × 105 mm	BS-010208-RK

Rack for rotors

RR-U	BS-010208-UK
-------------	--------------

* — Set of 2 adapters, made of **Ertacetal® POM-C** and is autoclavable

** — Set of 12 adapters, made of **Kocetal® POM**. Max. temperature +100°C



HOW TO CHOOSE ROTOR

LMC-56, Laboratory Centrifuge NEW

LMC-56 is a modern benchtop, low-speed laboratory centrifuge designed for convenient sedimentation, centrifugation and collection of necessary samples. It provides operation with tubes, blood collection systems, gel cards, microtest and ELISA plates. Our centrifuges are designed for safe work (metal protecting housing), easy maintenance and wide application range in medical, biochemical, chemical, industrial and other type of laboratories.

Speed regulation range 100–6,000 rpm (3,750 × g)

Speed resolution 100 rpm

Digital time setting 1–90 min (inc.1 min)

Timer sound signal +

Rotor imbalance diagnostics +
(automatic stop, "IMBALANCE" warning)

Rotor detection +

Display LCD, 2 × 16 signs

Chamber diameter 355 mm

Overall dimensions 560 × 480 × 315 mm
(W×D×H)

Weight 16.5 kg


Nominal operating voltage 230 V, 50/60 Hz
or 120 V, 50/60 Hz

Power consumption (230V / 120 V) 320 W / 330 W

Bucket rotor **BR-4U** with **BA-14/2U**
adapter for 2 ml tubes



Premium
Product Class

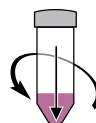


Features:

- User-friendly centrifugation parameter input and simultaneous display of the set and actual parameter values.
- Safe assay performance: metal protective housing and metal lid, automatic imbalance switch-off, lid lock during the centrifuge operation provide safe operation at all speeds.
- Rotor imbalance automatic diagnostics (emergency stop, imbalance indication).
- Automatic rotor detection with active rotational speed limit.
- Relatively high speed – maximum 6,000 rpm or 3,750 g.
- Wide choice of accessory rotors and adapters.
- A wide selection of rotors – up to 13.
- Improved chamber to reduce sample heating during centrifugation.
- Different modes of acceleration and deceleration, including deceleration mode with switched off forced braking.
- Possibility to set the speed both in revolutions per minute and by relative centrifugal force.



Interchangeable bucket adapters for LMC-56

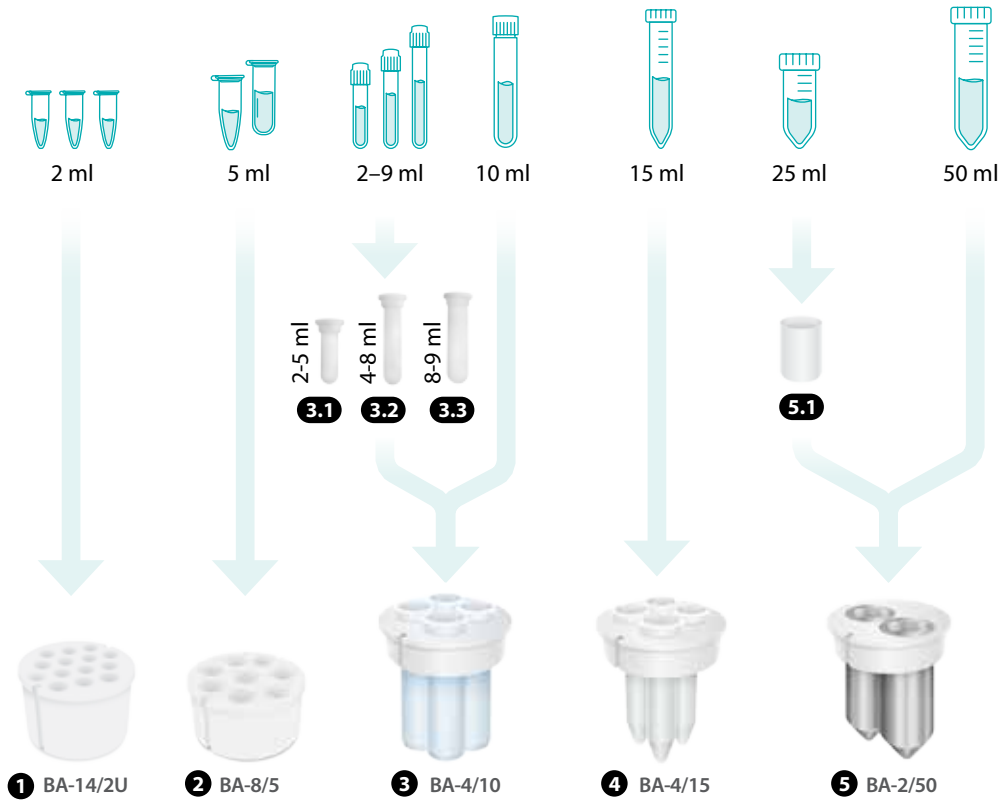

HOW TO CHOOSE BUCKET ADAPTER

ORDERING INFORMATION:



Rotor	Description	Capacity	Cat. number
BR-4U	Bucket rotor without bucket adapters	4	BS-010221-AK

Bucket adapters	Description	Capacity	Cat. number
1 BA-14/2U	Bucket adapters for 14 microtest tubes of 2 ml (set of 4 pcs.)	14 × 4	BS-010221-JK
2 BA-8/5	Bucket adapters for 8 tubes of 5 ml (set of 4 pcs.)	8 × 4	BS-010221-FK
3 BA-4/10	Bucket adapters for 4 tubes of 10 ml (set of 4 pcs.)	4 × 4	BS-010221-DK
4 BA-4/15	Bucket adapters for 4 tubes of 15 ml (set of 4 pcs.)	4 × 4	BS-010221-EK
5 BA-2/50	Bucket adapters for 2 tubes of 50 ml (set of 4 pcs.)	2 × 4	BS-010221-CK

Adapters for BA-4/10	Description	Cat. number
3.1 BN-13/75-16	For bucket adapters BA-4/10 optional adapter set (16 pcs) for vacutainers 2-5 ml	BS-010221-OK
3.2 BN-13/100-16	For bucket adapters BA-4/10 optional adapter set (16 pcs) for vacutainers 4-8 ml	BS-010221-PK
3.3 BN-16/100-16	For bucket adapters BA-4/10 optional adapter set (16 pcs) for vacutainers 8-9 ml	BS-010221-QK

Adapters for BA-2/50	Description	Cat. number
5.1 BI-25-8	For bucket adapters BA-2/50 optional adapter set (8 pcs) for 25 ml conical tubes	BS-010221-UK

Interchangeable Rotors for LMC-56

Rack RR-U




	Rotor Ri-12/10	Rotor Ri-24/10	Rotor Ri-6	Rotor Ri-6P
Rotor type	Swing-out			
Dimensions	Ø16 × 105 mm		Ø29 × 115 mm	
Capacity	12	24	6	
Tube's volume	10–15 ml		50 ml	
Max. speed	4,200 rpm	4,000 rpm	4,200 rpm	
Max. RCF:	3,160 × g	2,860 × g	3,160 × g	
Cat. number:	BS-010221-MK	BS-010221-LK	BS-010221-HK	BS-010221-IK



Plastic conical bottom centrifuge tube


Manufacturers: **Falcon, Greiner Bio-one, Sarstead, Corning, Nunc, TPP, etc.**

HOW TO CHOOSE ROTOR

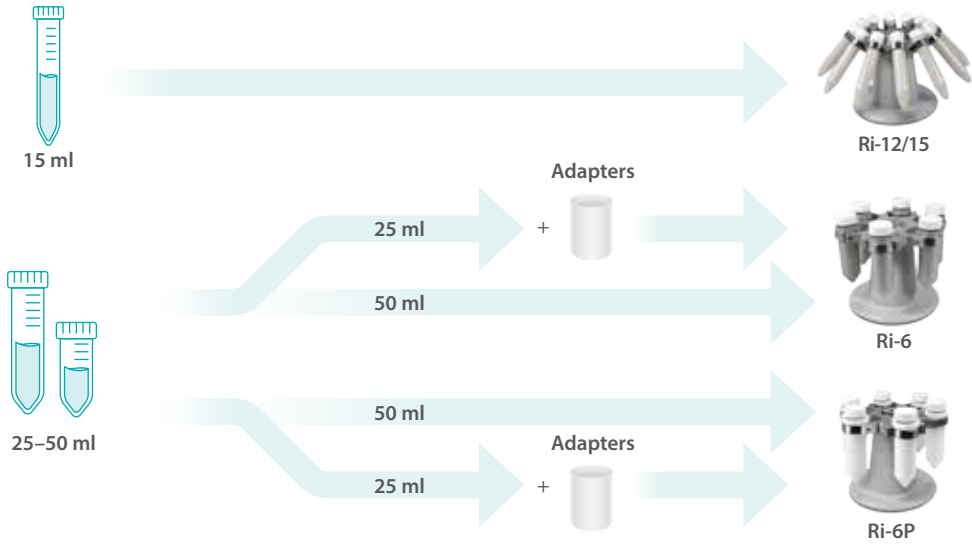

Plastic round bottom centrifuge tube, Vacutainers

Manufacturers: **Nunc, Greiner, Greiner Bio-one, TPP, etc.**

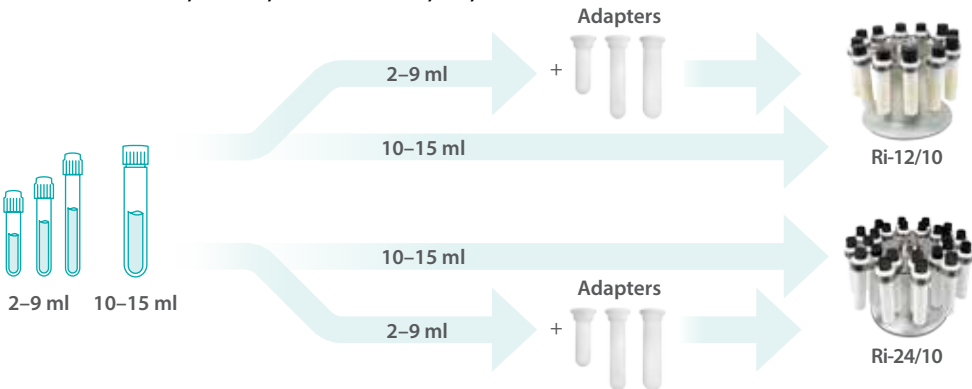

Interchangeable Rotors for LMC-56




Rotor Ri-2	Rotor Ri-24GC	Rotor Ri-12/15	Rotor RMT-24	Bucket rotor BR-4U
Swing-out		Angled Swing-out		Swing-out
128 × 85.6 mm	53 × 74 mm	Ø17 × 120 mm	Ø11.1	—
2	24	12	24	4 buckets
up to 45 mm	—	15 ml	1.5–2 ml	2–50 ml
2,000 rpm	1,500 rpm	4,200 rpm	6,000 rpm	4,200 rpm
560 × g	280 × g	3,160 × g	3,750 × g	2070–3260 × g
BS-010221-GK	BS-010221-NK	BS-010221-KK	BS-010221-BK	BS-010221-AK



Standard 96-well microtitre plates, skirted PCR plates and deepwell plates up to 45 mm


Manufacturers: **Nunc, Greiner, Greiner Bio-one, etc.**

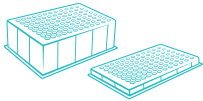

Ri-2

96-well semi-/ unskirted PCR plate

Manufacturers: **Nunc, Greiner, Greiner Bio-one, etc.**

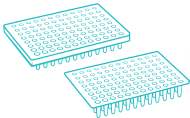


Material: **Ertacetal® POM-C** and is autoclavable



Ri-2

Gel Cards

Manufacturers: **Grifols®, DiaMed®, Bio-Rad® etc.**

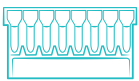

Ri-24GC, Rotor for Gel Cards for blood group serology testing (Forward Group, Reverse Group, RhD Type and 3 cell antibody screen). Recommended centrifugation time – 9 minutes



Ri-24GC



1.5–2 ml



RMT-24



2–50 ml













BR-4U

← Detailed information on **page 55**

Interchangeable Rotors for LMC-56

 ORDERING INFORMATION:

Adapters	Description	Cat. number
BI-25-6 	For rotors Ri-6 , Ri-6P optional insert set (6 pcs) for 25 ml conical tubes	BS-010221-VK
BN-13/75 	For rotors R-12/10, Ri-12/10 optional adapter set (12 pcs) for vacutainers 2-5 ml	BS-010208-PK
BN-13/100 	For rotors R-12/10, Ri-12/10 optional adapter set (12 pcs) for vacutainers 4-8 ml	BS-010208-QK
BN-16/100 	For rotors R-12/10, Ri-12/10 optional adapter set (12 pcs) for vacutainers 8-9 ml	BS-010208-RK
BN-13/75-24 	For rotors R-24/10, Ri-24/10 optional adapter set (24 pcs) for vacutainers 2-5 ml	BS-010221-RK
BN-13/100-24 	For rotors R-24/10, Ri-24/10 optional adapter set (24 pcs) for vacutainers 4-8 ml	BS-010221-SK
BN-16/100-24 	For rotors R-24/10, Ri-24/10 optional adapter set (24 pcs) for vacutainers 8-9 ml	BS-010221-TK
AP-96 	2 adapters for 96-well semiskirted and non-skirted PCR plates - made of Ertacetal® C. Autoclavable.	BS-010219-DK
AP-384 	2 adapters for 384-well plates made of Ertacetal® C. Autoclavable.	BS-010219-EK
RR-U 	Rack for rotors.	BS-010208-UK

Rotor RMT-24



Bucket rotor BR-4U with BA-2/50 adapter



THERMOSTATED EQUIPMENT:


THERMOSTATS – DRY BLOCK, HEATING/COOLING SYSTEMS



CH-100
Heating/Cooling Dry Block



CH 3-150
Heating and cooling thermostat

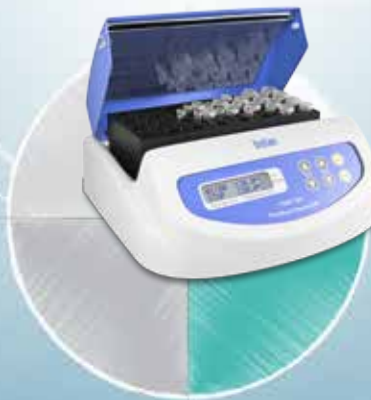

TDB-120
Dry block thermostat

Bio TDB-100 and TDB-120, Dry Block Thermostats

DESCRIPTION

Bio TDB-100 / TDB-120 — compact, easy-to-use thermostat for Eppendorf type micro tubes. It is specially designed for long incubation at different temperatures. The thermostat has an undeniable advantage in working with microquantities of reagents in microtubes. The thermostat possesses unprecedentedly high precision and uniformity of temperature distribution over the block.

With the help of the software-enabled temperature calibration function, the user can calibrate the unit in the range of several percent of the selected temperature to compensate for differences in the thermal behaviour of tubes from different manufacturers.



Blocks (built in) specifications:

Bio TDB-100

- 1 Block 24 × 2/1.5 ml + 15 × 0.5 ml + 10 × 0.2 microtubes

TDB-120

- 2 Block A-53 21 × 0.5 ml + 32 × 1.5 ml microtubes
- 3 Block A-103 21 × 0.5 ml + 32 × 1.5 ml + 52 × 0.2 ml microtubes

- 1 Block for Bio TDB-100



Basic Plus Product Class



Bio TDB-100



Bio TDB-100 product video

Heat up times for Bio TDB-100:



Basic Plus Product Class



TDB-120 with block A-103



TDB-120 product video

Heat up times for TDB-120:



Bio TDB-100 and TDB-120, Dry Block Thermostats

	Bio TDB-100	TDB-120
Temperature setting range	+25°C ... +100°C	+25°C ... +120°C
Temperature control range	5°C above ambient ... +100°C	5°C above ambient ... +120°C
Temperature setting resolution	0.1°C	
Temperature stability	±0.1°C	
Temperature uniformity @ +37°C	±0.1°C	
Temperature calibration coefficient range	0.936–1.063 (± 0.063)	0.968–1.031 (± 0.031)
Digital time setting	1 min – 96 h/no-stop (increment 1 min)	
Timer sound signal	yes	
Display	LCD, 2 × 16 signs	
Block capacity	24 × 2/1.5 ml + 15 × 0.5 ml + 10 × 0.2 ml microtubes	A-53 21 × 0.5 ml + 32 × 1.5 ml microtubes A-103 21 × 0.5 ml + 32 × 1.5 ml + 52 × 0.2 ml microtubes
Overall dimensions (W×D×H)	210 × 230 × 115 mm	230 × 210 × 110 mm
Weight	2.8 kg	
Nominal operating voltage	230 V, 50/60 Hz or 120 V, 50/60 Hz	
Power consumption	200 W (870 mA)	

ORDERING INFORMATION:

Cat. number 

Bio TDB-100 with built-in block

BS-010412-AAA

TDB-120 with built-in block A-103

BS-010401-QAA

TDB-120 with built-in block A-53

BS-010401-PAA



2 Block A-53



3 Block A-103



CH-100, Heating/Cooling Dry Block

DESCRIPTION

CH-100 is the result of combining two popular Biosan instruments:

1. Heating Dry block and
2. Cooling Dry block thermostat

The combined construction of the aluminium block and Peltier element module cooled with the forced ventilation radiator provides fast switching of the cooling and heating modes.

CH-100 is a very effective instrument for sample preparation during enzyme reactions, hybridization reactions, DNA analysis.

Microprocessor controlled time and temperature. Simultaneous indication of the and actual temperature and time.

SPECIFICATIONS

Temperature setting range	-10 °C ... +100 °C
Temperature control range	30°C below ambient ...+100°C
Temperature setting resolution	0.1°C
Temperature stability	±0.1°C
Temperature uniformity @ +37 °C	±0.25°C
Temperature calibration coefficient range	0.936–1.063 (±0.063)
Digital time setting	1 min–96 h/non-stop (increment 1 min)
Timer sound signal	yes
Display	LCD, 2 × 16 signs
Overall dimensions (W×D×H)	240 × 260 × 165 mm
Weight	3.2 kg
Input current/power consumption	12 V, 4.4 A / 55 W
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 12 V

Blocks (built in) capacity:

Block CH-1	20 × 0.5 ml +12 × 1.5 ml microtubes
Block CH-2	20 × 1.5 ml microtubes
Block CH-3	20 × 2 ml microtubes

ORDERING INFORMATION: Cat. number

CH-100 with block CH-1	BS-010410-BAI
CH-100 with block CH-2	BS-010410-CAI
CH-100 with block CH-3	BS-010410-UAI


Ice on block CH-2



Basic Plus
Product Class



Heat up and cool down times for CH-100:

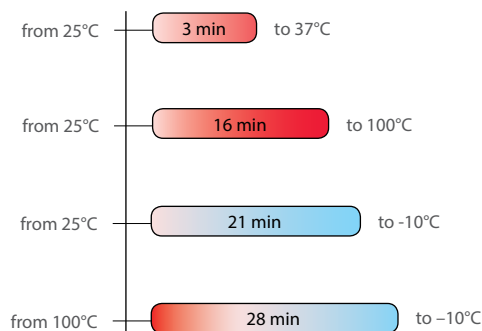

Product video






Product
video



Heat up and cool down times for CH3-150:

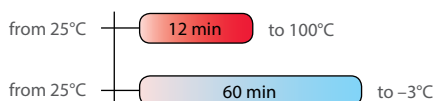

INTERCHANGEABLE THERMOBLOCKS:

- ❶ **B2-50** Ø48 mm × 2 sockets, depth 58 mm
- ❷ **B10-16** Ø16 mm × 10 sockets, depth 56 mm
- ❸ **B6-25** Ø25 mm × 6 sockets, depth 40 mm
- ❹ **B23-1.5** 23 sockets for 1.5 ml microtubes, depth 35 mm
- ❺ **B10-13** Ø13 mm × 10 sockets, flat bottom, depth 30 mm
- ❻ **B5-29** Ø29 mm × 5 sockets, flat bottom, depth 40 mm
- ❼ **B18-12** 18 sockets for Ø12 mm round bottom tubes, depth 58 mm

Different block types can be provided on request



ORDERING INFORMATION:

Cat. number

CH 3-150 without blocks BS-010418-AAA

Optional blocks:

- B2-50** BS-010418-AK
- B10-16** BS-010418-BK
- B6-25** BS-010418-CK
- B23-1.5** BS-010418-DK
- B10-13** BS-010418-LK
- B5-29** BS-010418-KK
- B18-12** BS-010418-EK

CH 3-150, Combitherm-2

Combitherm-2 **CH3-150** is specially designed to thermostabilise materials at temperatures from -3°C to $+150^{\circ}\text{C}$ according to analysis methods. To obtain useful functionality and decrease foot-print of instruments Combitherm-2 thermoblocks combined in a common electronic circuit board as well as inside a common external body. The front keyboard's left part is responsible for setting parameters for cooling plug-in blocks and the right part — for heating plug-in blocks. Both of them are regulated independently and can realise up to 16 programs, including temperature and time in each program. Peltier technology is used for cooling below room temperature; PCB is used for heating till $+150^{\circ}\text{C}$.

Separation of cooling and heating parts from each other increases durability of the instrument and speed of temperature changing after setting a new program.

Heating Block Specifications:

Temperature setting range	$+25^{\circ}\text{C} \dots +150^{\circ}\text{C}$
Temperature control range	5°C above ambient $\dots +150^{\circ}\text{C}$
Setting resolution	1°C
Stability	$\pm 0.1^{\circ}\text{C}$
Temperature calibration coefficient range	$0.936 \dots 1.063 (\pm 0.063)$

Cooling Block Specifications:

Temperature setting range	$-3^{\circ}\text{C} \dots +20^{\circ}\text{C}$
Temperature control range	23°C below ambient $\dots 5^{\circ}\text{C}$ below ambient
Setting resolution	0.1°C
Stability	$\pm 0.1^{\circ}\text{C}$

General Specifications

Digital time setting	1 min–99 h 59 min (increment 1 min)
Timer sound signal	yes
User adjustable programs (temperature and time)	16 (heating) +16 (cooling)
Display	LCD
Overall dimensions (W × D × H)	295 × 285 × 220 mm
Weight (without block)	5.6 kg
Nominal operating voltage	230 V, 50/60 Hz
Power consumption	430 W (1.8 A)

DESCRIPTION

SPECIFICATIONS

❶ B2-50



❷ B10-16



❸ B6-25



❹ B23-1.5



❺ B10-13



❻ B5-29



❼ B18-12



QB Series, Dry Block Heating Systems with Interchangeable Blocks

Equipment presented on pages 64–65 is produced by Grant Instruments (Cambridge) Ltd. Biosan is the sole distributor of Grant Instruments products in Russia, CIS and the Baltic States (Latvia, Lithuania, Estonia) and the official distributor for a number of other regions.

DESCRIPTION

A market-leading range of versatile, high-quality dry block heating systems with excellent temperature control, providing a source of precision heating for many sensitive analytical procedures.

A premium product range at an affordable price:

- Accurate, reproducible and safe heating of your samples — advanced temperature control combined with high quality, precision-engineered blocks providing excellent thermal contact;
- Versatile range of interchangeable heating blocks to fit any tube or plate you are using for your samples;
- Full range of models and options to cater for basic through to more sophisticated applications;
- Wide range of accessories.



Product video



QB4 with a lid



QBH2



Accessories






SPECIFICATIONS

Model (Cat. Num.)	QBD1 / QBD2 / QBD4	QBH2
Type	Digital	Digital
Number of blocks	1/2/4	2
Temperature range	amb. +5°C to 130°C	amb. +5°C to 200°C
Temperature setting range	+15°C to 130°C	+15°C to 200°C
Temperature stability @ 37°C	±0.1	±0.1
Temperature uniformity within the block @ 37°C	±0.1	±0.1
Display / Resolution	LED / 0.1°C	LED / 0.1°C
Safety: Overtemperature	Thermal fuse	
Timer with a sound alarm	1 min up to 72 h	
Heat up time from 25°C to 100°C	15 min	
Power consumption	150/300/600 W	300 W
Power supply	120 V or 230 V	

ORDERING INFORMATION:

Catalogue number matches the name of the product

QB Series, Dry Block Heating Systems with Interchangeable Blocks: Accessories

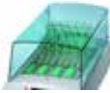Interchangeable blocks (Cat. Num.)		QBD1	QBD2	QBD4	QBH2	QBA1	QBA2
No. of blocks		1	2	4	2	1	2
QB-0 Plain block without holes		+	+	+	+	+	+
QB-10 24 × 10 mm Ø holes, 50 mm hole depth		+	+	+	+	+	+
QB-12 24 × 12 mm Ø holes, 50 mm hole depth		+	+	+	+	+	+
QB-13 12 × 13 mm Ø holes, 50 mm hole depth		+	+	+	+	+	+
QB-16 12 × 16 mm Ø holes, 50 mm hole depth		+	+	+	+	+	+
QB-17H for 10 × Falcon tubes tall 17 mm diam, 75 mm deep		+	+	+	+	+	+
QB-18 12 × 18 mm Ø holes, 50 mm hole depth		+	+	+	+	+	+
QB-24 5 × 24 mm Ø holes and universal bottles, 50 mm hole depth		+	+	+	+	+	+
QB-50 4 × 50 ml centrifuge tubes, glass universals, 50 mm hole depth		+	+	+	+	+	+
QB-H 56 × 0.2 ml microtube, 14 mm hole depth		+	+	+	+	+	+
QB-E0 24 × 0.5 ml microtube, 30 mm hole depth		+	+	+	+	+	+
QB-E1 24 × 1.5 ml microtube, 35 mm hole depth		+	+	+	+	+	+
QB-E2 24 × 2.0 ml microtube, 35 mm hole depth		+	+	+	+	+	+
QB-E5 12 × 5.0 ml microtube, 53.5 mm hole depth, 16.7 mm diameter		+	+	+	+	+	+
QB-DN Dolphin nose tube 24 × Ø 11.13 mm to Ø 6.1 mm		+	+	+	+	+	+
External Pt1000 temperature probe							
	Standard probe. For in-sample or in-block temperature control; encased in stainless steel sheath, Ø 3 mm × 30 mm long, with 350 mm of cable	+	+	+	+	-	-
	Short-form probe. For in-sample or in-block temperature control; encased in stainless steel sheath, Ø 3 mm × 14 mm long, with 350 mm of cable	+	+	+	+	-	-
Microplate blocks of molecular biology and biotechnology applications							
Double-size blocks 140 × 100 × 75 mm supplied with additional extraction tool							
	96 holes in microplate configuration for 0.2 ml microplates, strips or individual tubes. Uniformity ± 0.3°C within tubes across the block; 6.2 mm Ø holes, 14 mm hole depth	-	+	-	+	-	+
	Universal block for standard 96-well plates (u-well, v-well, flat bottom, high temperature) Uniformity ± 0.5°C between wells; supplied with hinged, double layer lid to create an insulated incubation chamber	-	+	-	+	-	+
Safety covers (not required with QDP-FL Microtiter blocks)							
	Made from tough clear acrylic for maximum visibility whilst preventing accidental touching of a hot block or contamination of samples from splashes. Clearance height 85 mm	QBL1	QBL2	QBL4	QBL2	QBL1	QBL2



THERMOSTATIC EQUIPMENT:

**WATER BATHS,
ORBITAL/LINEAR SHAKING BATHS,
UNSTIRRED WATER BATHS,
HEATING/COOLING CIRCULATORS**



Optima™ Series
Stirred Thermostatic Baths
and Heating Circulators



Ecocool
Energy Efficient Refrigerated /
Heating Circulating Baths




WB-4MS
Stirred water bath



How to choose thermostat?

How to choose thermostat?

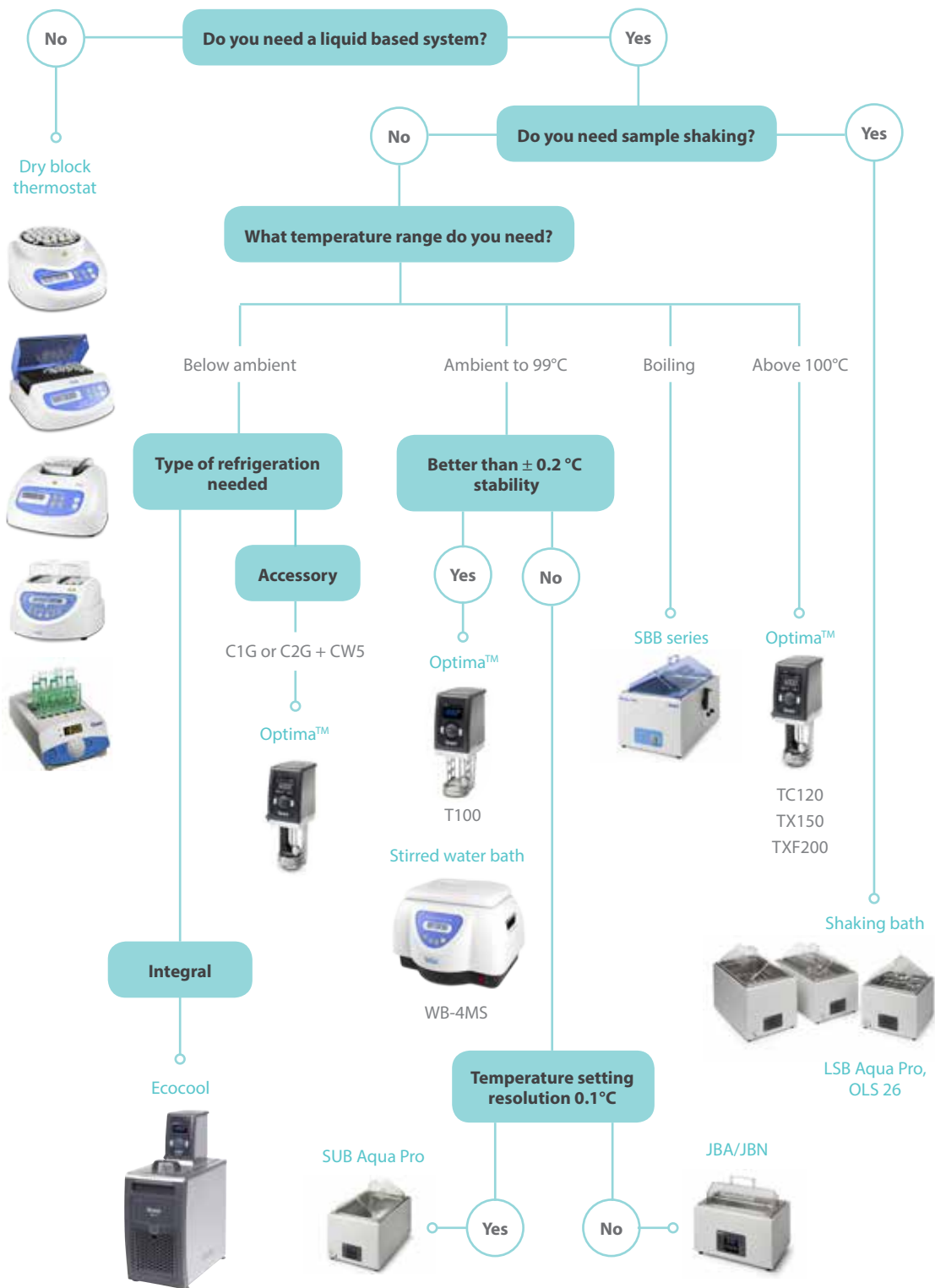

WB-4MS, Stirred water bath

DESCRIPTION

Stirred water bath **WB-4MS** is designed for chemical, pharmaceutical, medical and biological laboratory research, for processes requiring constant temperature ranging from ambient temperature to 100°C.

WB-4MS provides increased temperature stabilization (up to 0.1°C) due to a built-in magnetic stirrer (speed control range 250–1000 rpm).

Easy setup, high-temperature maintenance accuracy, compact size and attractive modern design make this water bath widely used.

SPECIFICATIONS

Tank capacity	4 l
Temperature setting range	+25°C ... +100°C
Temperature control range	5 °C above ambient ... +100°C
Temperature setting resolution	0.1°C
Temperature stability	±0.1°C
Temperature uniformity @ +37°C	±0.1°C
Stirring speed control range	250–1000 rpm
Digital time setting	1 min–96 h/non-stop (increment 1 min)
Timer sound signal	yes
Display	LCD, 2 × 16 signs
Digital setting of temperature, time and mixing speed	
Plastic lid with stainless steel interior included	
Quiet operation	
Working volume	235 × 135 × 110 mm
Overall dimensions (W×D×H)	340 × 270 × 250 mm
Weight	3.4 kg
Nominal operating voltage	230 V, 50/60 Hz or 120 V, 50/60 Hz
Power consumption	230 V, 50 Hz/600 W (2.6 A) 120 V, 60 Hz/670 W (5.6 A) 100 V, 50/60 Hz/600 W (6.0 A)
Maximum continuous operation time	24 h

Basic Plus

Product Class



WB-4MS with two adapters



Rack TR-16/19

ORDERING INFORMATION:

WB-4MS with base **BP-1** and lid

Optional accessories:

Cat. number

BS-010406-AAA

Racks:	Diameter/tube volume	Capacity	Cat. number
① TR-5/30	Ø 30 mm	5 tubes	BS-010406-KK
② TR-16/19	Ø 16–19 mm	16 tubes	BS-010406-FK
③ TR-30/13	Ø 10–13 mm	30 tubes	BS-010406-IK
④ TR-44/11	2/1,5 ml	44 tubes	BS-010406-JK

Combined Orbital/Linear Shaking Bath OLS26

DESCRIPTION



Patented, combined orbital and linear shaking mechanism of the **OLS26** allows optimisation of aeration and shear forces mixing for reproducible results.

- Precision digital temperature control;
- 0°C to 99°C operating range;
- Stability $\pm 0.1^\circ\text{C}$;
- Easy changeover from linear to orbital shaking;
- Adjustable shaking speed and stroke length;
- Polycarbonate lid included as standard;
- Drain tap for convenient emptying;
- 3 year warranty;
- TU26 included, other trays sold separately.

SPECIFICATIONS

Tank size	26 l
Minimum working depth	70 mm
Temperature control range	ambient +5 to 99°C. 0 to 99°C with accessory cooling
Temp. uniformity (DIN 12876-3) @ 70°C	$\pm 0.1^\circ\text{C}$
Temp. stability (DIN 12876-3) @ 70°C	$\pm 0.1^\circ\text{C}$
Display	2 x LED (individual displays and controls for temperature and shaking speed)
Orbital and Linear shaking speed	20 to 200 rpm (depending on load)
Orbital shaking radius	9 mm
Shaking speed display resolution	1 rpm
Linear shaking stroke length	18, 28, and 36 mm
Shaking tray area	380 x 235 mm
Timer	1 to 999 min
Dimensions (WxDxH)	335 x 590 x 475 mm
Heater power 120 V/230 V	1.05/1.4 kW
Drain tap	yes
Safety	over temperature protection/ low liquid level cut-out
Supply voltage	110–120 V or 220–230 V

ORDERING INFORMATION:

Cat. number



OLS26 with TU26 tray

OLS26

Equipment presented on pages 64–65, 71–83 is produced by Grant Instruments (Cambridge) Ltd. Biosan is the sole distributor of Grant Instruments products in Russia, CIS and the Baltic States (Latvia, Lithuania, Estonia) and the official distributor for a number of other regions.

Description and pictures of all available accessories can be found on page 73

Linear shaking bath — LSB Aqua Pro range

DESCRIPTION

World-renowned shaking water baths. High quality, robust design with unique magnetically coupled shaking mechanism for maximum reliability, consistency and quiet operation. Extensive range of accessories to provide the right solution for your application. Varied vessels types can be securely held using high-quality springs, clamps or racks.

FEATURES

- Ambient +5°C to 99°C operation;
- Stability $\pm 0.1^\circ\text{C}$;
- Choice of two models — 12 and 18 litres;
- Drain tap for convenient emptying;
- 3 year warranty;
- Polycarbonate lid included;
- Extensive choice of accessory shaking trays. Tray sold separately.



SPECIFICATIONS

	LSB12	LSB18
	9.2 kg W: 360 mm D: 385 mm H: 425 mm	11.2 kg W: 335 mm D: 565 mm H: 425 mm
Tank size	12 l	18 l
Minimum working depth	60 mm	
Temperature range	ambient +5 to 99°C	
Uniformity (DIN 12876-3) @ 70°C	$\pm 0.1^\circ\text{C}$	
Stability (DIN 12876-3) @ 70°C	$\pm 0.1^\circ\text{C}$	
Display	LED	
Linear shaking speed	20 to 200 strokes/min (depending on load)	
Shaking speed display resolution	1 strokes/min	
Linear shaking stroke length	20 mm	
Shaking tray area	240 × 235 mm	420 × 235 mm
Timer	1 to 999 min	
Heater power 120/230V	0.8/0.8 kW	1.05/1.4 kW
Drain tap	yes	
Safety	over-temperature protection/low liquid cut-out	
Supply voltage	110–120 V or 220–230 V	



ORDERING INFORMATION:

Cat. number

LSB12, Linear shaking bath 12 l with TU12 tray

LSB12

LSB18, Linear shaking bath 18 l with TU18 tray

LSB18

Description and pictures of all available accessories can be found on next page

Accessories for Shaking Baths: LSB 12, LSB 18 & OLS 26

Accessories LSB and OLS Aqua Pro Product/description		OLS26	LSB12	LSB18
		Catalogue number		
	Universal tray with adjustable springs. Highly versatile for a variety of vessel types.	TU26	TU12	TU18
	Flask/plate tray — with threaded holes to accept flask clamps or holder for deep well plates (≥ 2 ml). See option below.	TF26	TF12	TF18
	Test tube tray — compatible with SR racks or can be used alone to accommodate bags and miscellaneous vessels. See rack option below.	TS26 (holds up to 5 SR racks)	TS12 (holds up to 3 SR racks)	TS18 (holds up to 5 SR racks)
	Base tray — perforated stainless steel, allows bath to be used as an unstirred bath.	SBT26	SBT12	SBT26
	Cooling coil — the source of constant cooling to enable bath to be operated at or below ambient, down to 0 °C. LS200 lid (with an access hole for cooling coil) recommended.	CC26R	—	
	Heat exchange coil — attach to a cold water supply or refrigerated circulator. Can be used down to 2°C above the temperature of the coolant. LS200 lid (with an access hole for cooling coil) recommended.	CW26	—	
	Stainless steel sloping lid, gabled.	LS200	LU14	LU28
	Replacement polycarbonate lid, clear, gabled.	AQL26	AQL12	AQL26

Flask clamps and plate holder for TF tray

Cat. number	Description	OLS26 Capacity	LSB12 Capacity	LSB18 Capacity
SC-25	for 25 ml flask	28	18	33
SC-50	for 50 ml flask	24	14	26
SC-100	for 100 ml flask	15	9	17
SC-250	for 250 ml flask	8	5	14
SC-500	for 500 ml flask	6	4	6
SC-1000	for 1,000 ml flask	3	2	4
SH-DWP	1 × deep well plate	4	1	4

Test tube racks / microtube racks for TS tray

Cat. number	Tube diameter (mm)	Rack capacity
SR-10	10	48
SR-13	13	44
SR-16	16	24
SR-19	19	21
SR-25	25	12
SR-30	30	10
Cat. number	Microtube size (ml)	Rack capacity
SR-SE	0.5	119
SR-LE	1.5	48

ORDERING INFORMATION:

Catalogue number matches the name of the product 

Unstirred Water Bath



SUB Aqua Pro — advanced water bath range with a choice of eight models. Supplied with a base tray, lid and drain on larger bath.



JB Nova — general purpose water bath range with a choice of four models. Supplied with a base tray, lid and drain on larger bath.



JB Academy — basic range with a choice of three models. Supplied with a base tray.



SBB Aqua Plus boiling bath range — basic range with a choice of three models. Supplied with a base tray.

- The reliability, quality and consistent performance of Grant products have made Grant a leading manufacturer of water baths for decades;
- A new era for Grant water baths — now all models from basic to advanced with digital controls;
- Proven performance — technology to deliver temperature control you can rely on;
- Set and Forget™ technology — minimal bath setup, maximum time for your work.

More information about offered Grant Instruments products can be found on Biosan web-site www.biosan.lv

SUB Aqua Pro Digital Unstirred Water Bath











Built to the highest standard and specifications and incorporating the latest technology, the SUB Aqua Pro advanced water bath range supports even the most demanding applications, that require accurate temperature control. Choose from eight models with base tray and lid included as standard.

- Ambient +5°C to 99°C operation;
- Set and Forget™ technology — fast heat-up, accurate temperature control;
- Stability ±0.2°;
- Adjustable over temperature alarm — protect samples from overheating;
- Advanced dry start and run-dry protection;
- Three programmable temperature presets;
- 3 year warranty.

DESCRIPTION

SPECIFICATIONS

	SUB Aqua Pro digital unstirred water bath range – summary of specifications							
	SAP2	SAP2S	SAP5	SAP12	SAP18	SAP26	SAP34	SAPD
								
	2.5 kg W: 185 mm D: 200 mm H: 200 mm	3 kg W: 335 mm D: 215 mm H: 150 mm	3 kg W: 335 mm D: 215 mm H: 200 mm	6 kg W: 360 mm D: 380 mm H: 225 mm	9.5 kg W: 335 mm D: 590 mm H: 275 mm	9 kg W: 335 mm D: 590 mm H: 275 mm	14.5 kg W: 335 mm D: 590 mm H: 275 mm	9 kg W: 545 mm D: 380 mm H: 225 mm
Tank capacity	2 l	2 l (shallow)	5 l	12 l	18 l	26 l	34 l	5 l & 12 l
Temperature range	ambient t°C + 5 to 99							
Temp. display and setting resolution	0.1°C							
Temp stability (DIN 12876) @ 70 °C	±0.2°C							
Temperature setting/energy regulation	digital							
User adjustable over temp. alarm	+							
Fixed thermal cut-out	+							
Dry start/boil dry protection	+							
Programmable temp. presets	3							
Countdown timer with audible alarm	1 to 999 min							
Working area D×W (mm)	117 × 131	139 × 289	131 × 281	281 × 306	485 × 281	481 × 278	635 × 281	131 × 281 & 281 × 306
Minimum fill level	50 mm	32 mm	50 mm	50 mm	50 mm	70 mm	70 mm	50 mm
Maximum fill level	25 mm below the top of the tank							
Drain tap included	-	-	-	+	+	+	+	+
Heater power 120 V/230 V kW	0.25/0.25	0.35/0.35	0.35/0.35	0.8/0.8	1.4/1.05	1.4/1.05	1.8/1.3	1.15/1.15
Supply voltage V	120 or 230							

ORDERING INFORMATION:

Catalogue number matches the name of the product



SUB Aqua Pro Digital Unstirred Water Bath

OPTIONS AND ACCESSORIES

SAP2	SAP25	SAP5	SAP12	SAP18	SAP26	SAP34	SAPD
2 l	2 l (shallow)	5 l	12 l	18 l	26 l	34 l	5 l & 12 l

Replacement polycarbonate transparent lids*

AQL2	AQL5	AQL5	AQL12	AQL26	AQL26	—	AQL5, AQL12

Directs condensation away from immersed vessels, avoids contamination, reduces evaporation and saves energy

Stainless steel sloping lids*

—	LU6	LU6	LU14	LU28	LU28	LU36	LU6 & LU14

Flat lids*

—	—	LF6 (2 ring sets)	LF14 (4 ring sets)	LF28 (6 ring sets)	LF28 (6 ring sets)	LF36 (8 ring sets)	LF6 / LF14

With ring sets of variable hole diameter to accommodate tall vessels whilst reducing evaporation

Polypropylene spheres* (packs per bath)

1 × PS20	1 × PS20	1 × PS20	1 × PS20	2 × PS20	2 × PS20	3 × PS20	2 × PS20

Useful alternative to a lid, minimises evaporation and heat loss whilst allowing easy access to vessels in the bath; particularly useful for tall vessels

Raised shelves – reversible, allows two shelf depths. h = shelf height above tank base (mm)

—	—	—	RS14H (h 40 or 78) shelf covers half area of SAP12	RS18H (h 40 or 135) shelf covers half area of SAP18	RS28H (h 45 or 135) shelf covers half area of SAP26	RS36H (h 45 or 135) shelf covers half area of SAP34	RS14H (h 40 or 78) shelf covers half area of SAPD

Racks (no. per bath)

—	—	1 × J2	2 × J2	4 × J2	4 × J2	6 × J2	1 + 2 × J2

Choice of eight variants to accommodate different tube diameters and microtubes (see below)

Replacement base trays

AQBT2	AQBT5	AQBT5	AQBT12	AQBT26	AQBT26	SBT36	AQBT5 & AQBT12

Required if flat-bottomed flasks are to be placed directly on the base of the bath and to promote thermal convection in the bath

* — Lid or spheres recommended for use above 60°C

Unstirred Bath Racks

J2 Racks	Tube size Ø	Capacity	J2 Racks	Tube size Ø	Capacity
J2-10	10 mm	84	J2-25	25 mm	18
J2-13	13 mm	55	J2-30	30 mm	12
J2-16	16 mm	36	J2-SE	0.5 ml	105
J2-19	19 mm	32	J2-LE	1.5 ml	65

Optima™ Series, Stirred Thermostatic Baths and Heating Circulators



A cost-effective range of multi-purpose systems combining Grant's legendary quality and reliability. Precise temperature control for a wide range of laboratory applications.

- **Accurate and safe temperature control** — for samples and users;
- **Intuitive programming and thoughtful design features** — makes working with Grant heated baths and circulators easy;
- **Robust, durable construction** — for longevity, reliability and long-term low cost of ownership;
- **A complete range** — 32 models to cover basic through to sophisticated needs, each model represents excellent value for money.

APPLICATIONS

Grant stirred baths and circulators provide a source of precision heating and cooling for many routines and sensitive analytical procedures including sample incubation, calibration and quality control testing. All models from the TC120 upwards are suitable for unnecessary both open and closed-loop circulators (i.e. remote vessel open or closed).

For more powerful heating requirements, i.e. above 200 °C, contact marketing@biosan.lv for advice.

DESCRIPTION

Heating Circulators Specifications on page 78 and all available accessories on page 80

Model selection (see next page):

Any of the four **Grant Optima™** digital thermostats can be combined with any of eight Grant tanks (five stainless steel and three plastic) to provide a choice of 32 models.

Optima™ Series, Heating Circulators Specifications


SPECIFICATIONS

Grant Optima™ Heating Circulators Specifications		General purpose Digital		Digital High Performance	
		T100	TC120	TX150	TXF200
Stability (DIN 12876) @ 70°C	°C	±0.05	±0.05	±0.01	±0.01
Uniformity (DIN 12876) @ 70°C	°C	±0.1	±0.1	±0.05	±0.05
Setting resolution	°C	0.1	0.1	0.1 (0.01 with Labwise™)	
Display		4 digit LED		full colour QVGA TFT	
Timer function		—	1 to 6,000 min	1 min to 99 h 59 min	
No. preset temperatures		3	3	3	3
Re-calibration points		2	2	5	5
Offset adjustment		—	—	+	+
Socket for external probe (TXPEP, TXSEP)		—	—	+	+
Communication interface		—	—	USB & RS232	USB & RS232
Programmable		—	—	remote via PC/laptop 1 program/ 30 segments	direct via user interface or remote via PC/laptop 10 programs / 100 segments
Relays		—	—	1	1
Safety	overtemperature	fixed	adjustable cut-out		
Safety	fluid level — float switch	+	+	+	+
Alarms (can be configured to switch a relay)		—	high, without relay	high and low	high and low
Heater power 230 V	kW	1.3	1.3	1.9	1.9
Electrical power 230 V	kW	1.4 (50–60 Hz)	1.4 (50 Hz)	2.0 (50 Hz)	2.0 (50–60 Hz)
Height above tank rim	mm	200	200	200	200
Depth below tank rim	mm	135	135	135	135
Grant Optima™ thermostat pumps (integral)					
Maximum pressure	water, mbar	—	210	310	530
Maximum flow	water, l/min	—	16	18	23 (adjusted flow rate)
Pipe bore	inlet/outlet, mm	—	6/11	6/11	6/11
Dimensions (HxDxW)	mm	315 × 145 × 115			

ORDERING INFORMATION:

Cat. number:

T100 EURO

TC120 EURO

TX150 EURO













TXF200 EURO

Optima™ Series, Water Bath Combinations and Accessories

Capacity (l)	Outer tank dimensions 1. Dimensions (HxDxW) Weight (kg) 2. Working area (DxW) 3. Min/max fluid depths 4. Inner tank dimensions (HxDxW)	T100	TC120	TX150	TXF200
		Temperature setting range	Temperature setting range	Temperature setting range	Temperature setting range
ST5 – 5 l Stainless steel	1. 215 × 335 × 187 mm, 2.9 kg 2. 150 × 260 mm 3. 85/140 mm 4. 150 × 300 × 150 mm	T100–ST5 amb.+15 to 100 °C	TC120–ST5 0 to 120 °C	TX150–ST5 0 to 150 °C	TXF200–ST5 0 to 200 °C
ST12 – 12 l Stainless steel	1. 215 × 332 × 360 mm, 4.5 kg 2. 205 × 300 mm 3. 85/140 mm 4. 150 × 325 × 300 mm	T100–ST12 0 to 100 °C	TC120–ST12 0 to 120 °C	TX150–ST12 0 to 150 °C	TXF200–ST12 0 to 200 °C
ST18 – 18 l Stainless steel	1. 215 × 545 × 340 mm, 7.3 kg 2. 385 × 300 mm 3. 75/130 mm 4. 150 × 505 × 300 mm	T100–ST18 0 to 100 °C	TC120–ST18 0 to 120 °C	TX150–ST18 0 to 150 °C	TXF200–ST18 0 to 200 °C
ST26 – 26 l Stainless steel	1. 270 × 535 × 340 mm, 7.7 kg 2. 385 × 300 mm 3. 125/180 mm 4. 200 × 505 × 300 mm	T100–ST26 0 to 100 °C	TC120–ST26 –15 to 120 °C	TX150–ST26 –15 to 150 °C	TXF200–ST26 –15 to 200 °C
ST38 – 38 l Stainless steel	1. 260 × 733 × 338 mm, 11.9 kg 2. 575 × 300 mm 3. 125/180 mm 4. 200 × 690 × 300 mm	T100–S38 0 to 100 °C	TC120–S38 –15 to 120 °C	TX150–S38 –15 to 150 °C	TXF200–S38 –15 to 200 °C
P5 – 5 l Plastic	1. 180 × 323 × 220 mm, 2.2 kg 2. 120 × 150 mm 3. 85/140 mm 4. 155 × 240 × 160 mm	T100–P5 amb.+15 to 99 °C	TC120–P5 amb.+15 to 99 °C	TX150–P5 amb.+15 to 99 °C	TXF200–P5 amb.+15 to 99 °C
P12 – 12 l Plastic	1. 180 × 412 × 340 mm, 3.4 kg 2. 210 × 280 mm 3. 85/140 mm 4. 155 × 325 × 280 mm	T100–P12 amb.+5 to 99 °C	TC120–P12 amb.+5 to 99 °C	TX150–P12 amb.+5 to 99 °C	TXF200–P12 amb.+5 to 99 °C
P18 – 18 l Plastic	1. 180 × 589 × 340 mm, 5.1 kg 2. 375 × 280 mm 3. 85/140 mm 4. 155 × 510 × 290 mm	T100–P18 amb.+5 to 99 °C	TC120–P18 amb.+5 to 99 °C	TX150–P18 amb.+5 to 99 °C	TXF200–P18 amb.+5 to 99 °C
OPTIONS AND ACCESSORIES					
Labwise™ PC software (optional)					
Allows two-way communication for status display, programming and data capture		—	—	+	+
External probes (optional)					
TXPEP flexible plastic probe, 3 m cable		—	—	+	+
TXSEP stainless steel probe, 3 m cable		—	—	+	+
Remote switching device (optional)					
For switching appliances on and off (up to max. 8 Amps)		—	—	1	2
Vertical turbine pumps (optional)					
Low noise, compact design. Supplied with pipe connections and special lid for fitting to tank, pipe bore 12.7 mm					
VTP 1	max. pressure 1,000 mbar max. flow 9 l/min	+	Required only where application demands a higher pressure than that delivered by the internal pump to maintain flow		
VTP 2	max. pressure 1,650 mbar max. flow 12 l/min	+			



Optima™ Series, Water Bath Accessories

ACCESSORIES							
	Lids to help reduce evaporation/heat loss and avoid sample contamination	Polypropylene spheres (no. of packs required, 300 spheres in one pack)	Rack systems to optimise use of available bath capacity (no. of racks accommodated)	Raised shelves to allow shallow vessels to be accommodated	Accessory cooling systems to allow systems to operate at or below room temperature by means of cooling coil dipped into the bath; designed for minimal impact on working area		
					Refrigerated immersion coolers Consist of a cooling coil connected to a refrigeration unit by a flexible pipe. Extract heat continuously, with the bath control unit controlling temperature	Heat exchange coil Designed to be attached to a supply of cooling tap water or a refrigerated circulator	
					C1GR (0 to 40°C)	C2GR (-15 to 40°C)	CW5 (2°C above coolant temperature)
ST5 – 5 L stainless steel	STL5 flat stainless steel	1 × PS20	1 × QR	—		—	
ST12 – 12 L stainless steel	STL12 gabled, hinged (removable) stainless steel	1 × PS20	2 × VR	RS14		—	
ST18 – 18 L stainless steel	STL26 gabled, hinged (removable) stainless steel	2 × PS20	4 × VR	RS22		—	
ST26 – 26 L stainless steel	STL26 gabled, hinged (removable) stainless steel	2 × PS20	4 × VR	RS28			
ST38 – 38 L stainless steel	STL38 gabled, hinged (removable) stainless steel	3 × PS20	6 × VR	RS28 or RS38			
P5 – 5 L plastic	PL5 flat, stainless steel	1 × PS20	1 × QR	—	—	—	—
P12 – 12 L plastic	PL12 curved plastic	1 × PS20	2 × VR	RS14	—	—	—
P18 – 18 L plastic	PL18 curved plastic	2 × PS20	4 × VR	RS22	—	—	—

Ecocool

Energy Efficient Refrigerated / Heating Circulating Baths



- Choice of two models, temperature range -30 °C to +150 °C (model dependent);
- 3 year warranty with renowned service and support, no registration required;
- Active cooling through the whole temperature range;
- True energy saving of up to 80% against standard compressor units.

A new range of innovative, eco-friendly, refrigerated heating circulating baths offering significant running cost savings whilst delivering powerful cooling.

All products in the Ecocool range are supplied assembled as ready to use kits, complete with accessory hoses, clips and connectors as standard.

DESCRIPTION

SPECIFICATIONS

H: 640 mm D: 430 mm W: 245 mm		ecocool 100R	ecocool 150R
Temperature range	°C	-20 to 100	-30 to 150
Temperature stability	±°C	0.05	0.02
Refrigerant		R290	
Flow rate (max)	L/min	17	14 – 22 (adjustable)
Pump pressure (max)	mbar	250	530
Tank volume	L	4.5	5.5
Working Area (D x W)	mm	120 x 150	
Min/Max liquid level	mm	105/140	160/190
Cooling power	@ 20°C W	250	450
	@ 0°C W	200	350
	@ -10°C W	100	300
	@ -20°C W	50	200
Programs		—	1 x 30 segments via Labwise™
Communication interface		—	USB
Temperature probe socket		—	6 pin mini DIN
Display		4 digit LED	Full colour QVGA TFT
Languages		—	5 (EN, FR, DE, IT, ES)
Timer		1 min to 99 hrs 59 mins	
Temperature presets			3
Alarms		High	High and low
Electrical power (max) kW	120V/230V	2.16/2.07 (50–60 Hz)	2.28/2.76 (50–60 Hz)
Safety		Adjustable over temperature cut-out	
Ready to use kits		Assembled and supplied with standard tubing, insulation, clips and connectors	





Ecocool

Energy Efficient Refrigerated / Heating Circulating Baths

APPLICATIONS

- **PHARMACEUTICAL** — Mini pilot plant reactors
- **EDUCATION** — Rotary evaporator cooling, replacement of running tap water cooling, immersing small samples, photometry, chromatography systems
- **INDUSTRIAL** — QC testing, sample preparation, general cooling, reaction chemistry, temperature control, semi-conductor manufacturing, rheometry
- **FOOD** — Refractometry
- **LIFE-SCIENCE** — Electrophoresis cooling
- **HIGH TEMPERATURE COOLING** — Active up to 150°C



Options and accessories	ecocool 100R	ecocool 150R
Labwise™ PC software (optional)		
Allows two-way communication for status display, programming and data capture + USB cable provided	—	
External probes (optional)		
PEP plastic probe	—	+
SEP stainless steel probe	—	+
Vertical turbine pumps (optional) when pump is fitted, available working area is reduced.		
Low noise, compact design. Supplied with pipe connections and special lid for fitting to tank, pipe bore 12.7 mm		Required only where application demands higher pressure than that delivered by the internal to maintain flow.
VTP1-LT max. pressure 1,000 mbar; max. flow 9 l/min		Note: The optional VTP pumps will transfer additional heat to the baths and reduce the net cooling power of the refrigeration unit. The above figures must be taken into consideration when choosing the refrigeration unit, when ordering a VTP pump, please specify which refrigeration base unit it is to be used with.
VTP2-LT max. pressure 1,650 mbar; max. flow 12 l/min		Note: Other sizes of heat exchange coil can be made to your specification. Contact us for further information.
Hose Kits		
HOSE100 General purpose hose kit: -40 to 100 °C HOSE200 High temperature hose kit: -50 to 200 °C		Hose kit 2 x 2 m, assembled with Optima™ pump outlet plate and simple hose clips, no tools required.

Optima™ R series, Refrigerated Thermostatic Bath and Circulators



TXF200-R4R

Cost-effective and efficient multi-purpose systems for low-temperature applications.

- Powerful precision cooling, whether used in open-loop or closed-loop format
- Combining legendary quality, reliability and design for everyday usage — useful features, straightforward maintenance, compact design
- Robust, durable construction for longevity, reliability and long-term low cost of ownership
- Up to 4 years warranty

We recommend using the following liquids with refrigerated thermostatic baths and circulators:

- **-50 to 50°C:** Silicone oil — low viscosity (Bayer silicone M3);
- **-30 to 30°C:** 50% water 50% antifreeze (inhibited ethylene glycol);
- **0 to 30°C:** 80% water 20% antifreeze (inhibited ethylene glycol);
- **5 to 99.9°C:** Water.

Grant low-temperature circulators provide a source of precision cooling for many sensitive analytical procedures, including spectrophotometry, viscometry, refractometry and electrophoresis. They are suitable for use in both open and closed-loop circulation (i.e. remote vessel open or closed).

Alternatively, Grant RC series of recirculating chillers (closed circulators) can be used. These are generally needed for more powerful cooling requirements, e.g. the removal of mechanical or electrical heat produced in apparatus or machinery. Please contact marketing@biosan.lv for advice.

Model selection:

The R4R refrigeration range consist of one refrigeration unit which can be combined with four heating circulators to offer a temperature range of -30°C to 100°C.

Capacity (l) Outer tank dimensions	• Working area (LxW) • Min/max liquid depths • Weight	T100 H: 333 mm D: 172 mm W: 120 mm	TC120 H: 333 mm D: 172 mm W: 141 mm	TX150 H: 342 mm D: 172 mm W: 141 mm	TXF200 H: 342 mm D: 172 mm W: 141 mm
R4R – 20 l stainless steel H: 550 mm D: 515 mm W: 393 mm; <i>Cat.num.: R4</i>	• 230 × 305 mm • 80/140 mm • 40.6 kg	T100-R4R 0°C to 100°C	TC120-R4R -25°C to 100°C	TX150-R4R -30°C to 100°C	TXF200-R4R -30°C to 100°C
Options and accessories					
Labwise™ PC software (optional)					
Allows two-way communication for status display, programming and data capture + USB cable provided		—	—	+	+
External probes (optional)					
TXPEP flexible plastic probe, 3 m cable		—	—	+	+
TXSEP stainless steel probe, 3 m cable		—	—	+	+
Remote switching device (optional)					
For switching mains power appliances on and off (up to max. 8 Amps)		—	—	1	1
Vertical turbine pumps (optional)					
Low noise, compact design. Supplied with pipe connections and special lid for fitting to tank, pipe bore 12.7 mm		Required only where application demands a higher pressure than that delivered by the internal pump to maintain flow			
VTP 1 max. pressure 1,000 mbar; max. flow 9 l/min	+				
VTP 2 max. pressure 1,650 mbar; max. flow 12 l/min	+				


MAGNETIC STIRRERS, OVERHEAD STIRRER



Intelli-Stirrer MSH-300i
Magnetic Stirrer with hot plate



MMS-3000
High speed magnetic stirrer

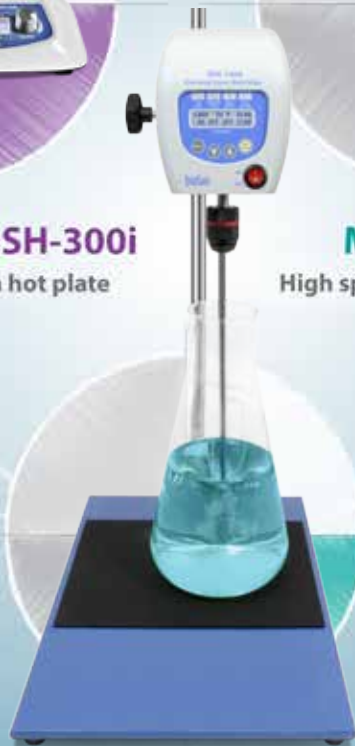

MM-1000
Overhead Stirrer Multi Mixer

MS-3000 and MMS-3000, Magnetic Stirrers

DESCRIPTION

MS-3000 and **MMS-3000** are compact magnetic stirrers with stainless steel working surface. Units provide stirring of liquids with the rotation speed of magnetic element up to 3,000 rpm. Up to date, it is the highest value of the maximal speed for magnetic stirrers of global producers.

Strong magnets hold the driven magnetic element firmly in the magnetic clutch. Stirring is performed without undesirable heating and noise.

Enclosures of stirrer **MS-3000** are made of strong steel and painted with powder enamel, which is chemically resistant to acids and alkali.

The stirrers are supplied with a cylinder-shaped magnetic stirring bar (6×25 mm) encapsulated in PTFE for universal use.

MMS-3000 is equipped with a detachable stand for supporting various sensor elements (temperature, pH and others) inside the stirred liquid.


Magnetic stirrer is an ideal laboratory instrument for PH-metering, extraction and dialysing with small quantities of substances.

Operation temperature range $+4^{\circ}\text{C}$ to $+40^{\circ}\text{C}$ (from cold rooms to incubators) at maximal relative humidity 80%.

Basic Plus
Product Class



Basic Plus
Product Class

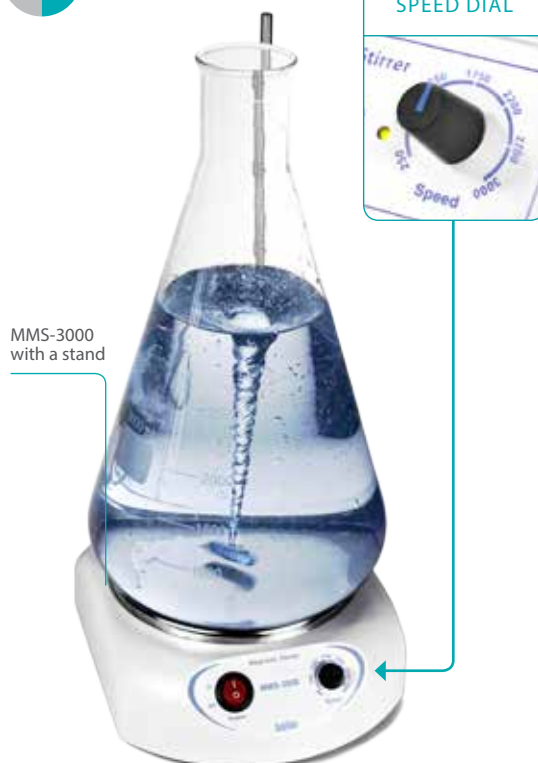

MS-3000 and MMS-3000, Magnetic Stirrers

	MS-3000	MMS-3000
Speed control range	0–3,000 rpm	
Stirring volume up to (H ₂ O)	5 L	20 L
Working surface material	Stainless steel	
SR-1, attachable stand size	—	Ø 8 × 320 mm
Max. length of magnetic stirring element (bar)	50 mm	70 mm
Stirring liquid viscosity	up to 1,170 mPa·s	
Maximum continuous operation time	24 h	
Operation in closed laboratory rooms	at ambient temperature from +4°C to +40°C	
Working plate size	110 × 110 mm	Ø 160 mm
Overall dimensions (W×D×H)	120 × 150 × 65 mm	185 × 230 × 75 mm
Weight	0.8 kg	1.5 kg
Input current/power consumption	12 V, 220 mA / 2.6 W	12 V, 250 mA / 3 W
External power supply	Input AC 100–240 V, 50/60 Hz; Output DC 12 V	

MMS-3000



MMS-3000



ORDERING INFORMATION:

MS-3000 white**MS-3000** blue (on request)**MMS-3000**

Optional accessories for MMS-3000:

HTP-1, Holder for temperature probe (see page 89)Cat. number 

BS-010301-AAF

BS-010301-ABF

BS-010305-AAF

BS-010309-FK

MSH-300 and Intelli-Stirrer MSH-300i, Magnetic Stirrers with hot plate

DESCRIPTION

MSH-300 and **Intelli-Stirrer MSH-300i** are magnetic stirrers of the new generation. Enclosures of stirrers are made of metal painted with powder enamel chemically resistant to acids and alkali. The stirrers are equipped with a detachable stand for supporting various sensor elements (temperature, pH and others) inside the stirred liquid.

The stirrers are supplied with a cylinder-shaped magnetic stirring bar (6×25 mm) for universal use covered with Teflon.

Units are equipped with overheat protection providing an automatic switch-off of the device when overheating for the set temperature difference occurs.

Magnetic stirrers with heating can be used for laboratory operations such as organic synthesis, extraction, analysis of oil products, pH-measurements, dialysis, soil suspending, preparing buffer solutions, etc.

Additional protection disables the heating if the temperature of the plate exceeds the set temperature for 30°C.

Operation temperature range +4°C to +40°C (from cold rooms to incubators) at maximal relative humidity 80%.



MSH-300i
product video

DESCRIPTION

Intelli-Stirrer MSH-300i is a digital version of magnetic stirrer with heating; it is designed for laboratories with higher requirements. It offers digital setting and control of temperature and rotation speed.

A powerful magnet allows mixing solutions with glycerine viscosity level. Maximum volume of stirred liquid (water) is 20 litres.

An external probe provides direct control of the stirred liquids temperature.

External temperature probe:

Probe type	thermocouple
Connection	type K

The cable is covered with Teflon, mechanically strong, elastic and chemically stable against oils, acids, aggressive reagents and liquids

Cable length	1 m
Operation temperature range	-50°C to +250°C

Basic Plus Product Class



Premium Product Class



MSH-300 and Intelli-Stirrer MSH-300i, Magnetic Stirrers with hot plate

	MSH-300	Intelli-Stirrer MSH-300i
Speed control range	250–1,250 rpm	100–1,250 rpm (10 rpm increment)
Max. stirring volume (H ₂ O)	15 L	20 L
Plate temperature regulation range	+30°C ... +330°C	+30°C ... +330°C (1°C increment)
Temperature control range with external probe	—	20 °C ... +150°C
Display	—	LCD
Temperature uniformity on the plate	±3°C	
Working plate heating time till 330°C	15 min	11 min
Diameter of working plate	160 mm	
Plate material	Aluminium alloy	
SR-1, attachable stand size	Ø 8 × 320 mm	
Length of magnetic stirring element	10–50 mm	20–70 mm
Max. stirring liquid viscosity	up to 1,170 mPa·s	
Maximum continuous operation time	24 h	168 h
Fault indication	Outputs sound signal and turns off the heating	Outputs an error code on display, turns off the heating
Overall dimensions (W×D×H)	190 × 270 × 100 mm	
Weight	2.9 kg	3.2 kg
Nominal operating voltage	230 V; 50/60 Hz or 120 V; 50/60 Hz	
Power consumption (Stirring)	8.5 W	
Power consumption (Heating)	550 W	

Connecting external probe to the Intelli-Stirrer MSH-300i




Plate heat up time for **MSH-300**:

from 25°C — **15 min** — to 330°C


Plate heat up time for Intelli-Stirrer **MSH-300i**:

from 25°C — **11 min** — to 330°C

External sensor installation:

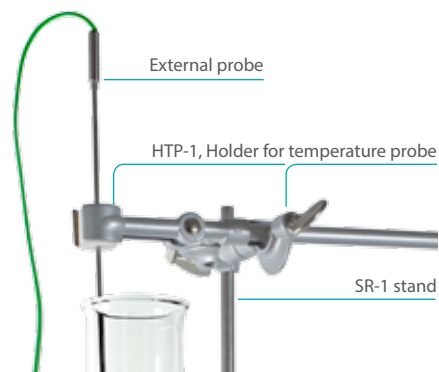

ORDERING INFORMATION:

Cat. number 

MSH-300 with stand BS-010302-OAA

Intelli-Stirrer MSH-300i with stand BS-010309-AAA

Optional accessories:

External temperature probe BS-010309-BK

HTP-1, holder for temperature probe BS-010309-FK

MM-1000, Overhead Stirrer Multi Mixer

DESCRIPTION

Overhead Stirrer Multi Mixer **MM-1000** is designed for stirring liquids up to 20 litres. Quiet and reliable mixer can provide stable continuous mixing up to seven days. It can realise three types of motion:

- 1 Rotational
- 2 Reciprocal
- 3 Vibration.

MM-1000 performs separate (mono-) (1; 2; 3), consecutive binary cycles (c) (1-2) × c; (1-3) × c and (2-3) × c and complex tri-cycles (1-2-3) × c.

Speed, angle and time of stirrer rotation are under micro-processor control.

Multi Mixer can be used for stirring solutions up to the "medium viscosity" range (from 1,000 to 10,000 mPa·s). It is an ideal instrument for biotechnology, organic synthesis, analytical laboratories.

The innovative combination of three motion types provides a high level of homogeneity due to consecutive combination of laminar and turbulent flows that cause substances to dissolve faster.

Electrically safe and energy efficient — powered by 12 V external power supply.

Specifications of movement types:

1 Rotation:

Speed regulation range 40–1,000 rpm

Time 0–250 s

2 Reciprocal motion:

Turning angle 0°–360° (increment 30°)

Time 0–250 s

3 Vibro motion:

Turning angle 0°–5° (increment 1°)

Timer 0–5 s

Timer sound signal yes

Stirring volume up to (H₂O) 20 L

Digital time setting 1 min–96 h/non-stop
(increment 1 min)


Overall dimensions (W×D×H) 140 × 135 × 250 mm

Weight 2.4 kg

Input current/power consumption 12 V, 700 mA / 8.4 W

External power supply Input AC 100–240 V 50/60 Hz,
Output DC 12 V

SPECIFICATIONS

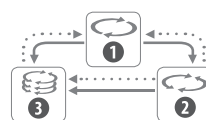

Multi mixing



Accessories for MM-1000


ORDERING INFORMATION:

Cat. number 

MM-1000 without stirrers

BS-010306-AAH

Optional accessories:	Type	Dimensions	Cat. number
A MP-1	Paddle stirrer	378 × (70 × 70) × 8 mm	BS-010306-AK
B MP-2	Propeller stirrer	2 folding blades (326 × 55 × 8 mm)	BS-010306-BK
C MP-3	Propeller stirrer	3 folding blades (325 × 50 × 8 mm)	BS-010306-CK
D MA-1	Anchor stirrer	332 × 90 × 8 mm	BS-010306-DK
E MC-1	Centrifugal stirrer	358 × 60 (110) × 8 mm	BS-010306-EK
F Double clamp	—	For device mounting	BS-010306-LK
G SRB, Support rod and base	—	For device mounting, 285 × 375 × 840 mm	BS-010306-KK

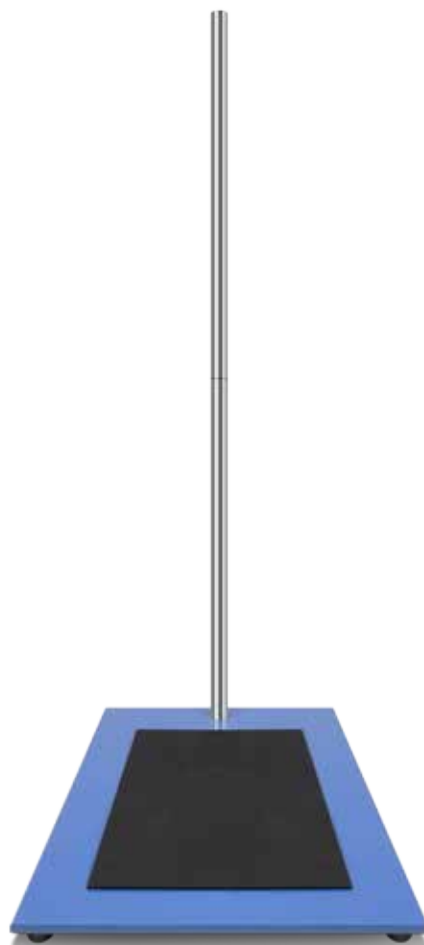A MP-1**B** MP-2**C** MP-3**G** SRB, Support rod and base**D** MA-1**E** MC-1**F** Double clamp

BIOSAFETY EQUIPMENT:

BIOSAFETY AIR, BIOSAFETY SURFACE, WATER PURIFICATION SYSTEMS



UVC/T-M-AR
DNA/RNA UV-cleaner box



UVT-S-AR
DNA/RNA UV-cleaner box



UVR-Mi
UV Cleaner-Recirculator

UVR-M and UVR-Mi, UV Cleaner-Recirculators

DESCRIPTION

How does UV-Air Flow Cleaner-Recirculator work?

Operation principle is based on a constant, forced air circulation through recirculator's chamber in close vicinity to UV lamps, thus ensuring maximal efficiency of disinfection. The inner mirror surface of the recirculator chamber reflects ultraviolet rays thereby increasing the UV radiation density and enhancing the disinfection effect.

What does UV Air Flow Cleaner-Recirculator consist of?

UV Air Flow Cleaner-Recirculator consists of a germicidal UV lamp, a fan unit equipped with dust filters and a control unit confined in a flow-through chamber.

What are the Benefits of UVR-M and UVR-Mi recirculators?

- UV Air Recirculators are ideal for air disinfection in hospitals (especially in outpatient departments, operating rooms, emergency rooms, delivery rooms etc.), kindergartens, research laboratories, veterinary clinics
- Recirculators are effective against common airborne diseases by disinfecting the air and efficiently destroying disease-causing agents (viruses, microorganisms) by UV radiation
- Provide complete protection from UV radiation
- Easy to install, operate and maintain. Very low noise level
- Built-in timer allows controlling the UV lamp operating time (UVR-Mi model)
- Digital control unit allows tracking overall UV lamp operating time (UVR-Mi model)


Recirculator fixation:

- Convenient fixation on walls (standard)
- Mounting on a movable stand (optional) **A**



UVR-Mi product video

Sensitivity of microorganisms to UV radiation intensity in UV air recirculators UVR-M and UVR-Mi

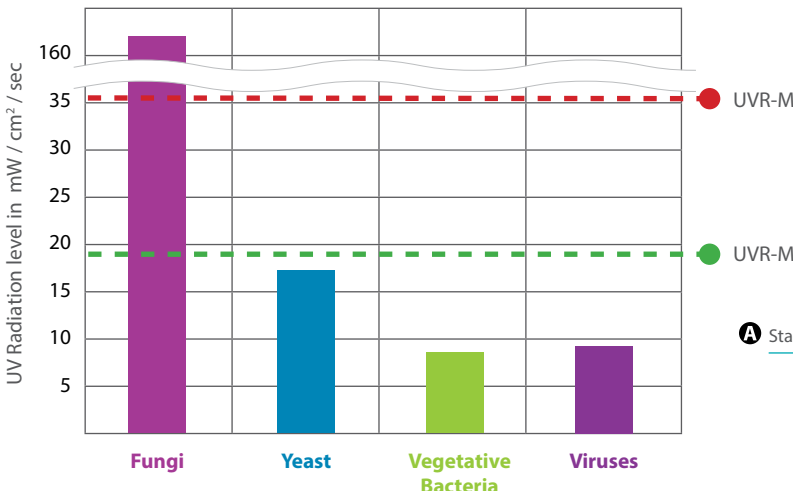


UVR-M and UVR-Mi, UV Cleaner–Recirculators

SPECIFICATIONS

	UVR-M	UVR-Mi
UV radiation source bactericidal UV-C, TUV 25W 1SL/25	1 lamp	2 lamps
UV radiation level	18 mW/cm ² /s	36 mW/cm ² /s
Air-flow productivity	21 m ³ /h	29 m ³ /h
Full user protection from direct UV light	Yes	
Display	—	LCD
UV lamp operation indicator	Yes	Yes
UV lamp lifetime counter	—	Yes
Timer	—	1 min–24 h/non-stop
Clock and ON/OFF scheduler	—	Yes
Automatic switch ON/OFF	—	Yes
Lamp fault detection	—	Yes
Overall dimensions (WxDxH)	110 × 135 × 660 mm	110 × 135 × 660 mm
Weight	2.7 kg	2.4 kg
Nominal operating voltage	230 V, 50 Hz or 120 V, 60 Hz	230 V, 50 Hz
Power consumption (230/120 V)	125 VA (540 mA)/160 VA (1.3 A)	110 W (0.5 A)

Operation principle

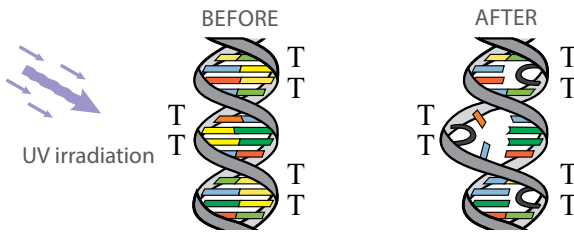

T — thymine formations

ORDERING INFORMATION: Cat. number

UVR-M BS-040105-AAA

UVR-Mi BS-040110-AAA

Optional adapters:

UVR-S (stand) BS-040105-AK

UVR-M and UVR-Mi, UV Air Recirculators Test Report.
biosan.lv/uvr-test



UVC/T-AR, UVC/T-M-AR, UVT-B-AR and UVT-S-AR DNA/RNA UV-Cleaner Boxes

DESCRIPTION

DNA/RNA UV-cleaner boxes (**UVC/T-AR**, **UVC/T-M-AR**, **UVT-B-AR** and **UVT-S-AR**) are designed for clean operations with DNA samples. They provide protection against contamination.

All models are bench-top type, made of metal framework, glass (or plexiglass) walls and working surface painted with powder enamel or made of stainless steel (See the specifications table on page 98).

UV-cleaner boxes are equipped with an open UV lamp installed in the upper hood. UV-radiation from the open lamps disinfects the working area inactivating DNA/RNA fragments during 15–30 min of exposure. A digital timer controls the duration of the direct UV irradiation. A daylight lamp provides proper illumination of the working surface.

UV-cleaner box is equipped with a flow-type bactericidal **UV cleaner–recirculator AR**, which provides constant decontamination inside the box during operation. They are recommended for operations with DNA/RNA amplicons.

UV cleaner-recirculator AR consists of a UV lamp, a fan and dust filters organized in a special body to protect a user working with a UV-cleaner box against UV light. Recirculator increases the maximum density of UV light, making it sufficiently effective for DNA/RNA inactivation. The UV–recirculator processes 100 UV-cleaner box volumes per hour, creating permanent aseptic operation conditions inside the UV-cleaner box.



Development and evaluation of DNA amplicon quantification video is available on the website

Advantages of Biosan UV-cleaner boxes:

- Ozone free high-density UV decontamination
- Long living UV lamps (9,000 hours average)
- Automatic switching off of UV-lamps when the protective screen is opened
- Bactericidal flow-type recirculator providing permanent decontamination inside UV-cleaner box during operation
- Shockproof glass walls
- Low noise, low energy consumption
- Tables for installation of UV-cleaner boxes
- UV-cleaner boxes with the bactericidal **UV cleaner–recirculator AR** is the patented Biosan solution

Premium
Product Class

UVC/T-M-AR



Basic Plus
Product Class

UVC/T-AR



Basic Plus
Product Class

UVT-B-AR



UVC/T-AR, UVC/T-M-AR, UVT-B-AR and UVT-S-AR DNA/RNA UV-Cleaner Boxes

Premium
Product Class

DNA/RNA UV-cleaner
box UVT-S-AR with equipment for nucleic acid extraction




Development and evaluation
of DNA amplicon
quantification.




biosan.lv/uv-box

Germicidal, shortwave (254 nm) ultraviolet energy is used for complete destruction of various biological agents

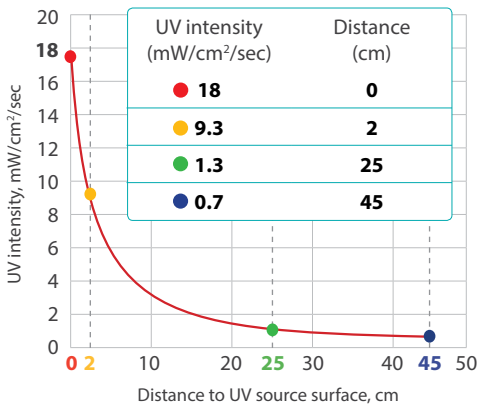

PER 1 SECOND

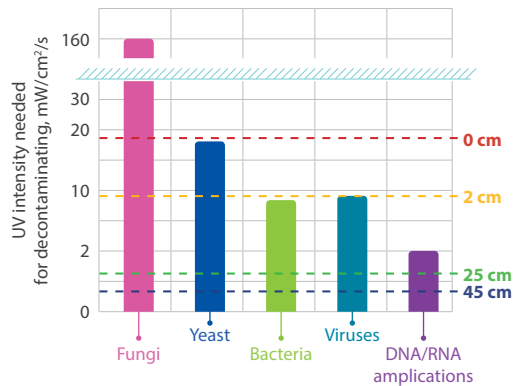

Average dosage for different surfaces

Surface	Dosage after 15 min	Dosage after 30 min
Working surface (40–50 cm)	570–680 mW/cm ²	1140–1360 mW/cm ²
Side walls (10–50 cm)	570–2500 mW/cm ²	1140–5000 mW/cm ²
Front window (10–50 cm)	570–2500 mW/cm ²	1140–5000 mW/cm ²

ORDERING INFORMATION:

Cat. number

UVC/T-AR with inlet	BS-040102-AAA
UVT-B-AR with internal socket and inlet	BS-040109-A06
UVC/T-M-AR with internal socket and inlet	BS-040104-A06
UVT-S-AR with internal sockets and inlet	BS-040107-AAA


UVC/T-AR, UVC/T-M-AR, UVT-B-AR and UVT-S-AR DNA/RNA UV-Cleaner Boxes



SPECIFICATIONS

Model	UVT-S-AR (double size)
Wall materials	Rear: stainless steel. Sides and front: glass (EUROGLASS, Germany)
Working surface material	Stainless steel
Open UV-lamp	2 × 30 W built-in bactericidal UV-C, TUV 30W 1SL/25
Recirculator UV radiation level	18 mW/cm ² /s
Radiation type	UV ($\lambda = 253.7$ nm), ozone-free
Digital time setting of direct UV exposure	1 min–24 h/non-stop (increment 1 min)
UV-recirculator	1 × 30W (efficiency >99% per 1 h)
Daylight lamp (for working area illumination)	1 × TLD-30W
Thickness of side panels	4 mm
Thickness of upper front panel	8 mm
Thickness of the front protective screen	5 mm
Optical transmission	95%
UV protection	>96% UV-protection film, type 4 mil, clear
Working area dimensions	1,210 × 500 mm
Opening size (W×H, fully raised protective screen)	1,185 × 190 mm
Safety features	Automatic open UV-lamp switching off when screen is open
Power outlets inside the unit (230/120 V)	3 built-in sockets max. 1000 W/600 W, Inlet for power cords
Nominal operating voltage	100–240 V, 50/60 Hz
Power consumption	135 W
Overall dimensions(W×D×H)	1250 × 600 × 590 mm
Weight (net/gross)	58/68,5 kg
Optional table	T-4L (W×D×H : 1290 × 600 × 770 mm)

UVT-S-AR

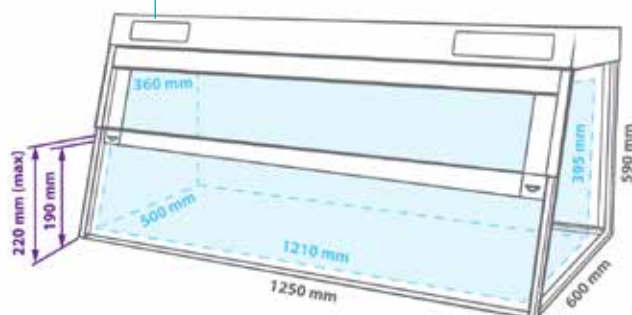


UVC/T-AR, UVC/T-M-AR, UVT-B-AR and UVT-S-AR DNA/RNA UV-Cleaner Boxes




SPECIFICATIONS

UVC/T-AR (compact)	UVC/T-M-AR (compact)	UVT-B-AR (compact)
Plexiglass: Polymethyl methacrylate (ALTUGLAS EX)	Rear: stainless steel. Sides and front: glass (EUROGLASS, Germany)	Rear: stainless steel. Sides: steel with chemical resistant powder coating. Front: glass (EUROGLASS, Germany)
Steel with chemical resistant powder coating	Stainless steel	
1 × 25 W built-in bactericidal UV-C, TUV 25W 1SL/25		
18 mW/cm ² /s		
UV (λ = 253.7 nm), ozone-free		
1 min–24 h/non-stop (increment 1 min)		
1 × 25W (efficiency >99% per 1 h)		
1 × TLD-15W		
4 mm	4 mm	2 mm
	8 mm	
8 mm	4 mm	4 mm
92%	95%	
>99.90% Polymethyl methacrylate ALTUGLAS EX	>96% UV-protection film, type 4 mil, clear	
645 × 490 mm		
645 × 165 mm	630 × 190 mm	
Automatic open UV-lamp switching off when screen is open		
Inlet for power cords	Inlet for power cords and 1 built-in socket, max. 1,000 W/600 W	
100–240 V, 50/60 Hz		
67 W		
720 × 535 × 555 mm	700 × 580 × 555 mm	
23/33 kg	28,8/39 kg	31,2/42 kg
T-4 (W×D×H : 800×600×745 mm)		

UVC/T-AR

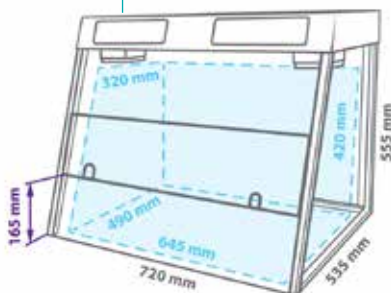

UVC/T-M-AR, UVT-B-AR

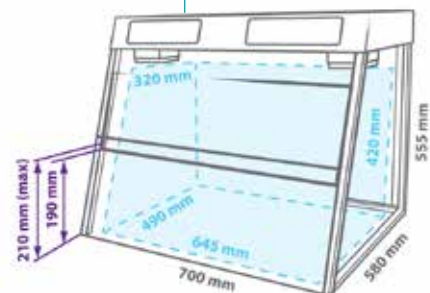

P-5, F-1, Shelves for DNA/RNA UV-cleaner boxes

DESCRIPTION

Two types of shelves have been developed for DNA/RNA UV-cleaner boxes to increase the effective area of the box: **P-5** — shelf-holder for five pipettes and **F-1** flat shelf.

On the **F-1** shelf, you can place laboratory glassware, reagents and other items that are convenient to keep in close proximity.


A P-5, shelf for pipettes:

Dimension (WxD) 230 × 140 mm

Capacity 5 pipettes

B F-1, shelf:

Dimension (WxD) 400 × 140 mm

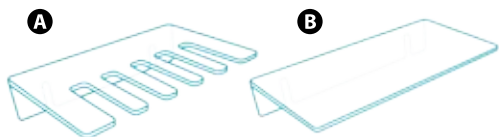

A P-5, shelf for pipettes



B F-1, shelf



ORDERING INFORMATION:

Cat. number

P-5, shelf for pipettes BS-040104-DK

F-1, shelf BS-040104-CK

PDS-250 and PDS-10L, DNA/RNA decontamination solution

DESCRIPTION

Contamination is especially problematic in the highly sensitive PCR technique. Originating from aerosolized fragments, contaminant DNA can lead to cross-contamination, thus resulting in inaccurate data and, as a result, misinterpreted analysis.

PDS is a ready-to-use solution for eliminating DNA and RNA from the surface prior PCR reaction preparation. DNA/RNA is removed within seconds after use. The solution contains a non-alkaline and non-carcinogenic agent. **PDS** is intended for use at PCR cabinets and laminars (e.g. UVT-S-AR), lab devices — BioMagPure 12, TS-100, pipettors — Assist series pipettes, etc.

PDS is effective against amplicon, plasmid, or genomic DNA and RNA from most surfaces except light or non-ferrous metals (e.g. aluminium, copper, lead, nickel, tin, titanium, zinc etc.).

The use of **PDS** both before and after PCR analysis is fast, easy and ideal to maintain a clean work area, thereby saving time and expenses.

PDS is heat resistant and stable for several years.

The decontamination solution is also available in 10 l containers — **PDS-10L**.



ORDERING INFORMATION:

Cat. number

PDS-250, DNA/RNA decontamination solution, spray 250 ml BS-040107-DK

PDS-10L, DNA/RNA decontamination solution, 10 l BS-040107-FK

Laboratory furniture

DESCRIPTION

A T-4, table



B T-4L, table



C LF-1, laboratory chest of drawers



Modular design of laboratory furniture provides flexibility and ease of use.

A T-4, table for — UVC/T-AR, UVC/T-M-AR, UVT-B-AR	
Maximum load	50 kg
Drawers	1
Mobility	Wheels with brakes
Material	Laminated particle board
Overall dimensions (WxDxH)	800 × 600 × 745 mm
Weight	23 kg

B T-4L, table for — UVT-S-AR	
Maximum load	75 kg
Drawers	1
Mobility	Wheels with brakes
Material	Laminated particle board
Overall dimensions (WxDxH)	1290 × 600 × 770 mm
Weight	36 kg

C LF-1, laboratory chest of drawers	
Drawers	5
Mobility	Wheels with brakes
Material	Laminated particle board
Overall dimensions (WxDxH)	300 × 450 × 705 mm
Weight	28 kg

UVT-S-AR on T-4L table with two LF-1 laboratory chests



ORDERING INFORMATION:

	Cat. number
T-4, table	BS-040101-BK
T-4L, table	BS-040107-BK
LF-1, laboratory chest of drawers	BS-050101-BK

Ultrapure water systems: Labaqua series

DESCRIPTION

Labaqua ultrapure systems are multi-purpose water purification systems. The Labaqua systems produce ultrapure and pure water directly from tap water.

Ultrapure (Grade 1) water is dispensed through the point-of-use filter on the front panel. Pure (Grade 2) water is dispensed directly from the storage tank.

Labaqua ultrapure water can be used for the most demanding applications, including, but not limited to: Inorganic trace analysis, Liquid chromatography, Cell culture, Molecular biology.

With resistivity of 18.2 Mega — Ohm \times cm (0.055 μ S/cm), ultrapure water produced by a Labaqua system exceeds requirements of all relevant standards (ISO 3696 Grade 1, ASTM Type I, CLSI Type I). Purified water is collected in a storage tank. An integrated recirculation system ensures consistent quality of water and reduces total organic carbon (TOC) to very low levels: <5ppb.

Pure water produced by the Labaqua systems complies with ISO 3696 Grade 2 water requirements and can be used for labware washing, wet chemistry methods, flame spectrophotometers, etc.

All cartridges and filters are easily accessible, and no tools are required to replace them. The Labaqua system can be installed on a laboratory bench or mounted on a wall.

FEATURES:

- **Volumetric dispense** — enables the user to set accurate dispensing volume for each dispense cycle. The dispense volume can be set either from the keyboard or by using the "teaching" mode.
- **Water quality** — embedded recirculation loop ensures stable premium water quality and enables practical elimination of Total Organic Carbon (TOC).
- **Low running costs** — performance of the deionization and polishing modules is constantly monitored. Monitoring algorithm enables cutting running costs, as replacement of the modules is requested only when service life is close to the end.
- **Total organic carbon (TOC) monitor** — organic contaminants may not affect the conductivity of water, so conductivity sensors cannot be used for TOC monitoring. Therefore, a special TOC monitoring module is needed to measure TOC level.
- **Color graphic LCD display** — system component status is reflected on the display in an intuitive colour pattern (Green/Yellow/Red).
- **System flowchart** — shows all component status and water quality parameters at a glance.



The Labaqua systems include:

- Boost pump
- Pre-filter set
- Reverse osmosis module
- Deionization module
- Final stage polishing module
- 30 L storage tank with an integrated Grade 2 dispensing valve
- Recirculation system

Model specific modules:

- **Labaqua Trace** — Point-of-use microfilter
- **Labaqua HPLC** — Point-of-use microfilter, TOC monitor
- **Labaqua Bio** — Point-of-use ultrafilter, UV sterilization module, TOC monitor



Ultrapure water systems: Labaqua series

SPECIFICATIONS

Purified water specifications	Labaqua Trace		Labaqua HPLC	Labaqua Bio
Ultrapure (Grade 1) water resistivity			18.2 MΩ × cm	
Ultrapure (Grade 1) water conductivity			0.055 μS/cm	
Pure (Grade 2) water resistivity			>10 MΩ × cm	
Pure (Grade 2) water conductivity			<0.1 μS/cm	
TOC	<30 ppb		<5 ppb	
RNase	—		—	<0.01 ng/ml
DNase	—		—	<4 pg/ml
Bacteria	<1 CFU/ml			<0.01 CFU/ml
Endotoxins	<0.15 EU/ml			<0.001 EU/ml
Particles >0.22 μm			<1/ml	
Deionization module life (standard module)			1 m ³	
Dimensions (W×D×H)			320 × 560 × 620 mm	
Storage tank			30 l	
Feed water pressure			0.8–4 bar	
Feed water conductivity			<1,300 μS/cm	
Weight	24 kg		25 kg	26 kg
Nominal operating voltage			230 V, 50/60 Hz	
Power consumption			130 W	

	Application	Labaqua Trace	Labaqua HPLC	Labaqua Bio
General laboratory applications	Glassware rinsing	+	+	+
	Laboratory washers	+	+	+
	Autoclaves	+	+	+
	Electrochemistry	+	+	+
	Wet chemistry	+	+	+
	Spectrophotometry	+	+	+
	Buffer and media preparation	+	+	+
	Reagent preparation	+	+	+
Inorganic analysis methods	Flame atomic absorption spectrophotometry	+	+	+
	Graphite atomizer atomic absorption spectrophotometry	+	+	+
	Plasma mass-spectrometry (ICPMS)	+	+	+
	Plasma spectrophotometry (ICPOES)	+	+	+
	Ion chromatography	+	+	+
Organic analysis methods	Liquid chromatography (HPLC/ UHPLC)		+	+
	Gas chromatography		+	+
	Total organic carbon measurements		+	+
Molecular Biology	Flow cytometry			+
	Cell and tissue culture			+
	Molecular biology			+

Ultrapure water systems: Labaqua series

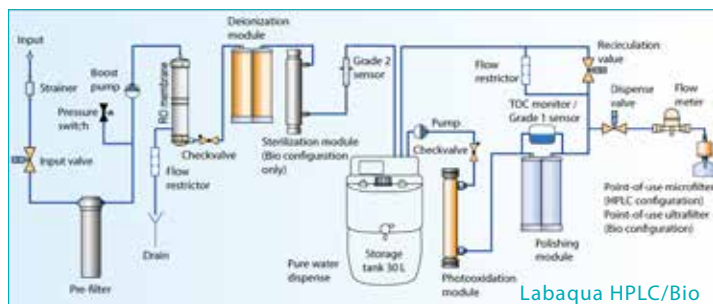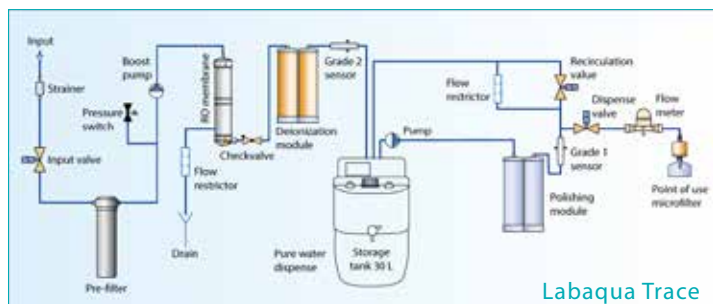

Optional universal remote dispenser



ORDERING INFORMATION

Labaqua Trace include 30 l tank, power cord

Labaqua HPLC include 30 l tank, power cord

Labaqua Bio include 30 l tank, power cord

Water purification system sets including different storage tanks are available

Optional accessories:

External pre-filter set (polyphosphate/carbon/1 μm) with manometer

External pre-filter set (carbon/1 μm) with manometer

Storage tank "Comfort" with base, tap and multipoint level switch, 60 l

Storage tank "Comfort" with base, tap and multipoint level switch, 100 l

Universal remote dispenser set with 3 m supply hose and water distribution module

External input pressure reducer

Replacement parts

Internal prefilter set

Deionization module

Polishing module

RO membrane

Point-of-use microfilter – 0.22 μm non sterile

Point-of-use microfilter – 0.22 μm sterile

Point-of-use ultrafilter

UV bulb 254 nm

UV bulb 185 nm

0.22 μm air vent filter for the storage tank

Air filter for storage tank

Filter set for BS-070104-LK (Polyphosphate, Carbon/PP, PP 1 μm)

Filter set for BS-070104-KK (Carbon/PP, PP 1 μm)

Cat. number

BS-070105-A02

BS-070104-A02

BS-070106-A02

please inquire

BS-070104-LK

BS-070104-KK

BS-070102-SK

BS-070102-FK

BS-070104-JK

10175

BS-070104-AK

BS-070104-IK

BS-070104-BK

BS-070104-NK

BS-070104-EK

BS-070104-FK

BS-070104-GK

BS-070104-RK

BS-070104-DK

BS-070102-AK

BS-070104-PK

410223

410222

DENSITOMETERS, PHOTOMETER



DEN-1
Densitometer



DEN-1B
Densitometer



DEN-600
Photometer

DEN-1 and DEN-1B, McFarland Densitometers

Densitometers **DEN-1** and **DEN-1B** are designed for measurement of cell suspension's turbidity in the range:

0.0–6.0 McFarland units
(0–180 × 10⁷ cells/ml);

Densitometers provide the opportunity to measure solution turbidity in a wider range (up to 15.0 McFarland units), however, it is necessary to remember that, in this case, the standard deviation values increase.

A densitometer is used for measurement of cell concentration (bacterial, yeast cells) during the fermentation process, determination of microorganism sensitivity to antibiotics, microorganism identification using various test-systems, for measurement of absorption at the definite wavelength, as well as for quantitative estimation of the colour solution concentration, absorbing green light.

The operation principle is based on the measurement of optical density with digital presentation of results in McFarland units. The unit is calibrated at the factory (for operation with 16 mm diameter glass tubes) and keeps calibration without power supply. However, if necessary, it is possible to calibrate the unit by 2–6 points in 0.0–6.0 McFarland unit range. Both commercial standards offered by Biosan and the cell suspensions prepared in a laboratory can be used for calibration.

Following polymer microparticles calibration kits and glass tubes are available on request:

- **CKG16** for glass tubes with diameter 16 mm, set of 0.5; 1.0; 2.0; 3.0; 4.0 McFarland Turbidity Standards (latex particles)
- **Calibration kit** for glass tubes with diameter 18 mm, set of 0.5; 1.0; 2.0; 3.0; 4.0; 5.0 McFarland Turbidity Standards (BaSO₄)
- **Calibration kit** for glass tubes with diameter 12 mm, set of 0.0 (blank); 0.5; 2.0; 3.0 McFarland Turbidity Standards (latex particles)
- **Glass sample tubes** without lid (diameter 16 mm, height 100 mm), which are suitable for working with **DEN-1**, **DEN-1B** factory calibration

Up to date information on calibration kits can be found on the website: www.biosan.lv

Two versions of the product are available:

1. **DEN-1** powered from external energy supply;
2. **DEN-1B** powered both from external energy supply and batteries (AA).

Basic Plus
Product Class



Premium
Product Class



Adapter A-16



DEN-1B rear side with calibration controls



DEN-1 and DEN-1B, McFarland Densitometers

	DEN-1	DEN-1B
Light source	LED	
Wavelength	$\lambda = 565 \pm 15 \text{ nm}$	
Measurement range	0.00–15.00 McF	
Display resolution	0.01 McF	
Accuracy	(0.0–6.0 McF) $\pm 3\%$	
Measurement time	1 s	
Sample volume	not less than 2 ml	
Tube external diameter	18 mm (without adapter) 16 mm (using included A-16 adapter) 12 mm (using optional A-12 adapter)	
Possibility to restore factory calibration settings		
Display	LCD	
Overall dimensions (W×D×H)	165×115×75 mm	
Weight	0.7 kg	
Independent power supply	—	3×AA batteries
Input current/power consumption	12 V, 7 mA/0.1 W	
External power supply	Input AC 100–240 V, 50/60 Hz, Output DC 12 V	
Standard set	Adapter A-16 External power supply	Adapter A-16 External power supply and 3×AA batteries



Application of **DEN-1** for determining microbial cells concentration of supernatant in tubes during centrifugation. Turbidity is determined in McFarland units.



DEN-1B

ORDERING INFORMATION:

Cat. number



DEN-1 with **A-16** adapter

BS-050102-AAF

DEN-1B with **A-16** adapter

BS-050104-AAF

Optional accessories:

A-12 adapter for 12 mm tubes

BS-050102-IK

CKG16 for glass tubes with diameter 16 mm, set of 0.5; 1.0; 2.0; 3.0; 4.0 McF

BS-050102-BK

Calibration kit for glass tubes with diameter 18 mm, set of 0.5; 1.0; 2.0; 3.0; 4.0; 5.0 McF

70900

Calibration kit for glass tubes with diameter 12 mm, set of 0.0; 0.5; 2.0; 3.0 McF

21255

Glass sample tubes 16 mm diameter without lid (16 × 100 × 0.8 mm), 78 pcs.

BS-050102-LK

Fits DEN-1, DEN-1B factory calibrated

DEN-600, DEN-600m Photometers NEW

DESCRIPTION

DEN-600/DEN-600m is a compact, portable, rechargeable battery-powered photometer. It comprises of 600 nm wavelength optical system, which enables to apply —

1. OD₆₀₀ method estimates the total number of cells.
2. McFarland (McF) turbidity measurement method.
3. Bradford protein assay method for protein concentration measurement.
4. Other methods that can be adjusted or optimized using 600 nm wavelength.

The device serves as an inexpensive alternative to a spectrophotometer commonly used for these applications. Because **DEN-600/DEN-600m** is battery powered and compact, it can be comfortably located in a biosafety cabinet, anaerobic chamber or quickly moved to another lab room. Additionally, the vessel holding mechanism allows accommodating round bottom, conical vials or falcon tubes, therefore enabling to measure the absorbance (Abs) and turbidity in Abs, OD and McFarland units.

USB connectivity and DEN software allow for data transfer, data processing and calculation, software calibration for Bradford protein assay method or a custom calibration for a specifically

SPECIFICATIONS

Light source	LED
Photodetector	Silicone photodiode
Measurement wavelength (λ)	600 nm ±10 nm
Vessel type	Cuvettes, round bottom tubes, falcon tubes
Battery type	LiPo
PC system requirements:	Intel/AMD Processor, 1 GB RAM, Windows Vista/7/8/10/11, USB
Dimensions (W×D×H)	120 × 145 × 65 mm
Weight	0.5 kg
External power supply	Input AC 100–240 V 50/60 Hz, Output DC 12 V

Measurement modes	Absorbance	McFarland
Measurement range	0–3.0 Abs	0–16.00 McF
Resolution	0.001 Abs	0.01 McF
Accuracy	±0.006 @ 1 Abs	±0.1 @ 0–8 McF
Repeatability	±0.003 @ 1 Abs	±0.05 @ 0–8 McF




IVD DEN-600m
IVD version



COMMON APPLICATIONS:

- Cell concentration measurement
- Cell growth data estimation
- Log phase estimation for microbial cells induction
- Competent cell preparation
- Bradford protein assay method
- Antibiotic susceptibility testing
- Inhibitory tests

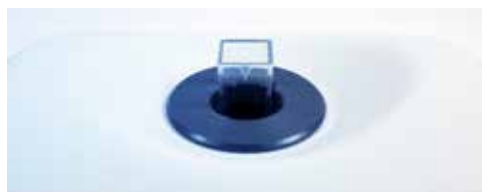

ORDERING INFORMATION

	Cat. number
DEN-600	BS-050109-AAA
DEN-600m IVD version	BS-050112-A01
Optional accessories:	
Verification set for Abs. Certified reference material, neutral density glass filter set of 4 Abs calibration points - 0.3532, 1.0512, 2.0425, 2.927 (the values may vary slightly from batch to batch)	BS-050109-AK
CKG16 for glass tubes with diameter 16 mm, set of 0.5; 1.0; 2.0; 3.0; 4.0 McF	BS-050102-BK
Glass sample tubes 16 mm diameter without lid (16 × 100 × 0.8 mm), 78 pcs.	BS-050102-LK
Calibration kit for glass tubes with diameter 12 mm, set of 0.0; 0.5; 2.0; 3.0 McF	21255

ASPIRATORS, PIPETTES



Assist
Pipette series



Assistboy
Pipette controller



FTA-2i
Aspirator with Trap Flask

FTA-2i, Aspirator with Trap Flask

DESCRIPTION

Aspirator with trap flask **FTA-2i** is designed for aspiration or removal of alcohol, buffer and liquid from reaction vessels (e.g. during DNA/RNA purification or other macromolecule reprecipitation techniques).

The device can be applied for routine operations of cells washing from culture medium and resuspension in a buffer. Aspirator operation principle is based on creating negative pressure in trapping flask using built-in microcompressor. The collecting tip is connected with polyethylene tube to the trapping flask. Liquid is removed from the reaction vessel when the collecting tip is in contact with the solution. A tube holder-organizer is conveniently located at **FTA-2i** right-hand side; it accommodates two 1.5–2 ml tube slots (e.g. for hydrochloric acid solution and distillate) necessary for collecting tip washing and storing, so that a tip can be re-used.

FTA-2i is equipped with a level sensor that detects excess liquid with consequential prevention of the overflow by automatically switching off the pump with a sounding alarm indication.

The devices come, as standard, with a vacuum regulation control knob that allows to select a preferable aspiration speed smoothly.

Additionally, a hand operator can be purchased for a more comfortable usability of the new accessories. The autoclavable hand operator features a pressure-sensitive button that can control the aspiration speed.

Common applications:

Removal and disposal of liquid from various reaction vessels

SPECIFICATIONS

Aspiration speed	up to 10 l/min (air)
Vacuum regulation	-200 to -800 mbar (adjustable)
Trap flask	2 l, polypropylene (autoclavable)
Liquid level sensor type	Invasive
Timer sound signal	yes
Overflow protection	Motor stops, light and sound signal
Filtration: Hydrophobic microbiologic filter 2200/02	eliminates risk of contamination from the trap flask by bacteria, viruses and infected particles
Filter pore diameter	0.027 micron
Input current/power consumption	12 V, 1 A/10.8 W
External power supply	Input AC 100–240V 50/60 Hz; Output DC 12 V
Dimensions (WxDxH)	185 × 290 × 390 mm
Weight*	1.85 kg

* — Accurate within ±10%.



Product video



LIQUID LEVEL
SENSOR

OVERFLOW
PROTECTION

2 LITRES

HYDROPHOBIC
MICROBIOLOGIC
FILTER


TUBE / TIP
HOLDER-ORGANIZER
(1.5 & 2 mL)

FTA-2i, Aspirator with Trap Flask

Optional accessories:

HAS-1, hand operator set

- ① Handheld vacuum controller;
- ② 1-channel adapter (with ejector) for 200 μ L tips;
- ③ 1-channel adapter with 125 mm stainless steel pin;
- ④ 1-channel adapter with 40 mm stainless steel pin;
- ⑤ 8-channel adapter (with ejector) for 200 μ L tips;
- ⑥ 8-channel adapter with 35 mm stainless steel pin;
- ⑦ Stand SH-6.

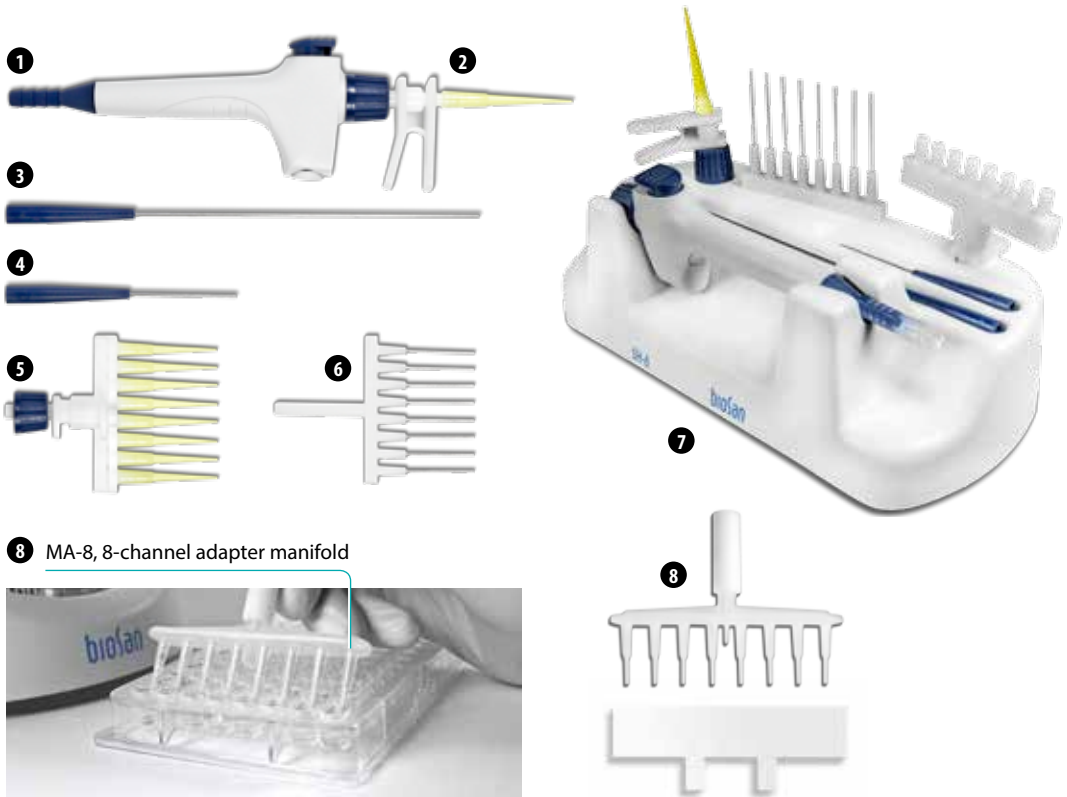

⑧ MA-8, 8-channel adapter manifold



ORDERING INFORMATION:

FTA-2i, with 2l trap flask, universal adapter **MA-U** (for 200/1000 μ L single use tips)

Cat. number

BS-040120-A02

Optional accessories:

HAS-1, hand operator set

BS-040118-PK

MA-8, 8-channel adapter manifold

BS-040108-BK

Extended tubing 2 m long, with fittings and MA-U adapter

BS-040120-DK

Replacement parts:

Suction microbiologic hydrophobic filter

BS-040120-S10

MA-U, universal adapter for 200/1000 μ L single use tips

BS-040118-AK

FTA-1, Aspirator with Trap Flask

DESCRIPTION

Aspirator with trap flask **FTA-1** is designed for aspiration/removal of alcohol/buffer remaining quantities from micro-test tube walls during DNA, RNA purification and other macromolecule reprecipitation techniques.

The device can also be used for routine operations of cells washing from culture medium and resuspension in a buffer. Aspirator operation principle is based on creating negative pressure in trapping flask using built-in microcompressor. The collecting tip is connected with polyethylene tube to the trapping flask. Liquid is removed from the microtest tube when the collecting tip touches the solution surface. A tube holder-organizer is conveniently located at **FTA-1** right-hand side; it accommodates two tubes (e.g. for hydrochloric acid solution and distillate) necessary for collecting tip washing and storing, so that a tip can be reused.

① Suction microbiological hydrophobic filter type 2200/02: Suction microbiologic filter eliminates the risk of contamination with bacteria, viruses and infected particle from patient to suction pump or central vacuum distribution. Suction microbiological filter is hydrophobic with very high bacterial blocking efficiencies, up to 99.99999% particles bigger than 0.027 µm (which is smaller than Hepatitis A, B and C).

SPECIFICATIONS

Vacuum	-500 mbar
Trap flask volume	1 l
Dimensions with trap flask (WxDxH)	160 × 210 × 340 mm
Weight with trap flask	1.7 kg
Input current/power consumption	12 V, 300 mA / 3.6 W
External power supply	Input AC 100–240 V; 50/60 Hz; Output DC 12 V



ORDERING INFORMATION:

Cat. number

FTA-1 with 1l trap flask BS-040108-AAG

Optional accessories:

MA-8 BS-040108-BK

Replacement parts:

Suction microbiologic hydrophobic filter BS-040108-S25

Basic Plus
Product Class

① Suction hydrophobic
microbiological filter



Product
video

Optional 8-channel adapter manifold MA-8



I LOVE PIPETTE



Assist, pipette series

The Assist series pipettes are single, 8 or 12 channel variable volume pipettes designed to measure and transfer volumes.

Single-channel pipettes are produced in ten ranges of volumes from 0.1 μl to 10,000 μl depending on the model.

Multichannel pipettes are produced in four ranges of volumes: 0.5–10 μl , 5–50 μl , 20–200 μl , 50–300 μl .

The pipettes are equipped with an analogue counter which shows the pipetting volume. The volume setting is done by turning the pipetting pushbutton knob or the black adjustment knob in the right direction. The volume range is shown on the pipetting pushbutton.

DESCRIPTION

Common pipettes usage depending on the volume

AP2, AP10, AP8-10, AP12-10	Measurement and transfer of micro-volumes, DNA sequencing and enzyme-assay applications.
AP20, AP50, AP100, AP200, AP250, AP1000, AP8-50, AP12-50, AP8-200, AP12-10, AP8-300, AP12-300	Measurement and transfer of general aqueous solution, acids and bases.
AP5000, AP10000	Measurement and transfer of large volumes.


Pipette:	Volume (μl)	Colour code	Fit to tips	Cat. number
Single channel:				
AP2	0.1–2.0	●	10 μl	BS-010501
AP10	0.5–10.0	●		BS-010502
AP20	2–20	●	200 μl	BS-010503
AP50	5–50	●		BS-010504
AP100	10–100	●		BS-010505
AP200	20–200	●		BS-010506
AP250	50–250	●	300 μl	BS-010507
AP1000	100–1,000	●	1,000 μl	BS-010508
AP5000	500–5,000	○	5,000 μl	BS-010509
AP10000	1,000–10,000	○	10,000 μl	BS-010510
Multichannel:				
AP8-10 AP12-10	0.5–10	—	10 μl	BS-010511 BS-010512
AP8-50 AP12-50	5–50	—	200 μl	BS-010513 BS-010514
AP8-200 AP12-200	20–200	—		BS-010515 BS-010516
AP8-300 AP12-300	50–300	—	300 μl	BS-010517 BS-010518
Sets:				
AP10, AP20, AP200, AP1000, 4 position stand, demo tips				BS-010519
AP10, AP100, AP1000, AP5000, 4 position stand, demo tips				BS-010520

Assist, pipette series

Features:

- Contoured shape of the handle and light weight;
- Proven accuracy and precision;
- UV resistant & fully autoclavable;
- 5 & 10 ml shaft protected by filter;
- Available in 8 and 12-channel version;
- Colour coded for easy volume identification;
- The adjustable ejector height system – to accommodate virtually all brands of tips;
- Dual volume setting using the pushbutton or the thumbwheel;
- Soft spring system for smooth, effortless pipetting.

Colour coded for easy volume identification

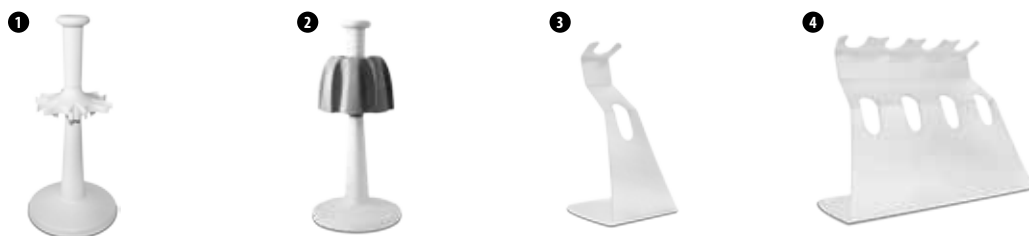


ORDERING INFORMATION:

Pipette stands:	Cat. number
1 Carousel stand (rotating) for 6 pipettes	BS-010522
2 Multiple stand (fixed) for 8 pipettes	BS-010523
3 1-position stand	BS-010524
4 4-position stand	BS-010525

Pipette tips features:

- The tip is made of imported high-quality medical grade polypropylene, and the filter element is made of ultra-high molecular weight polyethylene, which has good hydrophobicity.
- The filter element effectively protects the pipette from the sample to ensure the safety of sample aspiration.
- The filter element can reduce the hazard of the gas residue in the tip body to the operator during the pipetting process and improve the repeatability of the experiment.
- No inner surface coating, will not pollute the sample
- Have precise tick marks
- DNase, RNase and Pyrogen – free
- Can be autoclaved
- Suitable for conventional pipettes such as Biosan (Assist series) Eppendorf and Gilson.
- Provided in sterile racks.

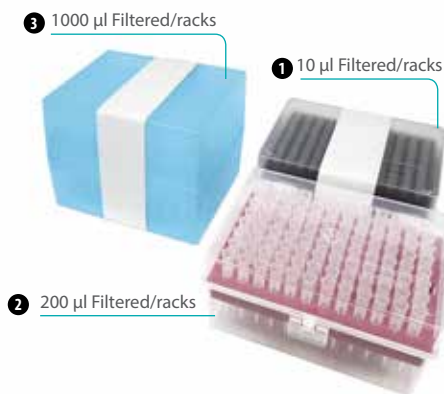

ORDERING INFORMATION:

Pipette tips:	Size (HxØ)	Packaging	Color	Compatibility	Cat. number
1 10 µl Filtered/racks	32 × 0.4 mm	96 tips/rack	Natural	Universal	BHZ01R1W-FS
2 200 µl Filtered/racks	59 × 0.5 mm	96 tips/rack	Natural	Universal	BHZ03R1W-FS
3 1000 µl Filtered/racks	86 × 0.8 mm	96 tips/rack	Natural	Universal	BHZ05R1W-FS

Assistboy, pipette controller



Assistboy pipette controller is a device intended for pipetting liquids with the use of measuring pipettes. It can work with all types of glass or plastic serological pipettes in the volume range from 0.5 ml to 100 ml.

Controller is equipped with an exchangeable filter membrane which protects shaft mechanism from aggressive liquid fumes.

Two dispense modes permit selection of dispensing intensity depending on the user's needs. The selected setting of the pipette controller mode is shown on display.

Safe and efficient work

- Protected by a PTFE filter blocking any liquid from entering the unit
- Autoclavable filter, the pipette holder and the nosepiece
- UV resistant body for safe sterilization
- Powerful, environmentally friendly 3 Ni-MH batteries enable many hours of continuous work
- LCD display showing battery charge level

Speed and working mode adjustment

- Function buttons for SPEED and working MODE control in a reach of a thumb
- Additional speed adjustment by the pressure applied to the trigger buttons

Working comfort

- Suitable for glass & plastic volumetric pipettes 0.5–100 ml
- Ergonomically shaped handle
- Well located function buttons
- Convenient charging stand

Charging stand



ORDERING INFORMATION:

Cat. number 

Assistboy with charging stand

BS-010521



BIOPROCESSING:

CO₂ INCUBATOR, SHAKER-INCUBATORS, PERSONAL BIOREACTORS



S-Bt Smart Biotherm
Compact CO₂ Incubator



ES-20/80
Shaker-Incubator



RTS-1 and RTS-1C
Personal Bioreactors

RTS-1 and RTS-1C, Personal bioreactors



Reverse-Spin®
Innovative Mixing Technology



USB
connection



Product
video

Users articles: biosan.lv/report

SPECIFICATIONS

	RTS-1	RTS-1C
<i>E.coli</i> BL21 Factory calibration measurement range, in OD ₈₅₀ : at 10–20 ml volume at 20–30 ml volume		0–10 OD (0–19 OD ₆₀₀ equivalent) 0–8 OD (0–15.2 OD ₆₀₀ equivalent)
Factory calibration measurement precision	±0.3 OD ₈₅₀	
Mass transfer coefficient k _a (h ⁻¹)	Up to 350 ±26 h ⁻¹ at 5 ml	
Measurement Wavelength (λ)	850 ±15 nm	
Light source	LED	
Real time measurement	1–60 min	
Temperature setting range	+25°C ... +70°C (increment 0.1°C)	+4°C ... +70°C (increment 0.1°C)
Bottom control range point	5°C above ambient	15°C below ambient
Top control range point	70°C	
Stability	±0.°C	
Sample temperature accuracy: 20–45°C <20°C >45°C	±1 ±2 ±3	
Sample temperature heating/cooling rate	0.7°C/min	
Sample volume	10–30 ml	
Speed control range	50–2,000 rpm (increment 10 rpm)	
Speed control precision	±15 rpm	
Reverse-Spin Time	1–60 s (increment 1 s)	
Display	LCD	
Minimum PC requirements	Intel/AMD Processor, 1 GB RAM, Windows Vista/7/8/8.1/10, 2.0 USB port	
Optimal PC requirements	Intel/AMD Processor, 3 GB RAM, Windows 7/8/8.1/10, 2.0 USB port	
Overall dimensions (WxDxH)	130 × 212 × 200 mm	
Weight	1.7 kg	2.2 kg
Input current/power consumption	12 V DC, 3.3 A/40 W	12 V DC, 5 A/60 W
External power supply	Input AC 100–240 V 50/60 Hz; Output DC 12 V	

* — Highest k_a (h⁻¹) is achieved at 5 ml working volume which is optimal for aerobic cultivation

** — Conversion coefficients from OD₈₅₀ to OD₆₀₀ vary between strains and phases of growth



Reverse-Spin® Technology — Innovative Principle of Microbial Cultivation
on web page biosan.lv/rts-tech

RTS-1 and RTS-1C, Personal bioreactors

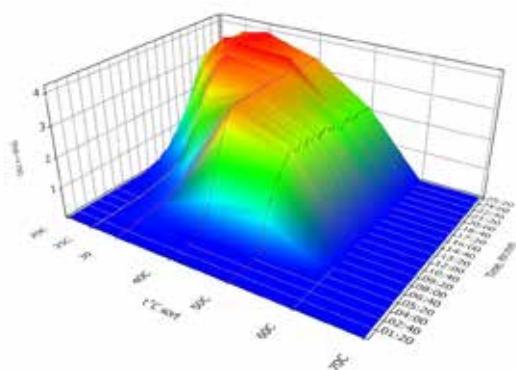
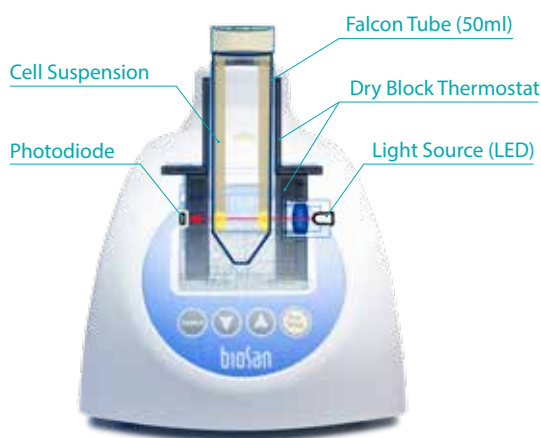


Figure 1. 3D graph of *E.coli* BL21 growth kinetics showing the effect of different temperatures in 7 parallel RTS bioreactors.

TYPICAL APPLICATIONS

- Fermentation real-time growth kinetics
- Clone candidate screening
- Protein expression
- Temperature stress and fluctuation experiments
- Media screening and optimization
- Growth characterization
- Inhibition and toxicity tests
- Strain quality control

ORDERING INFORMATION

RTS-1C including TubeSpin® Bioreactor 50, TPP®, 20 pcs.

RTS-1 including TubeSpin® Bioreactor 50, TPP®, 20 pcs.

Optional accessories:

TubeSpin® Bioreactor 50, TPP®, 20 pcs.

TubeSpin® Bioreactor 50, TPP®, 180 pcs.

USB 2.0 Hub 10 × ports

RTS-1 and RTS-1C are personal bioreactors that utilize patented Reverse-Spin® technology that applies non-invasive, mechanically driven, low energy consumption, innovative type of agitation where cell suspension is mixed by the single-use falcon bioreactor tube rotation around its axis with a change of direction of rotation motion resulting in highly efficient mixing and oxygenation for aerobic cultivation. Combined with a near-infrared optical system, it is possible to register cell growth kinetics non-invasively in real-time.

FEATURES

- Reverse-Spin® mixing principle in 50 ml falcon tubes allows to achieve high $k_L a$ (h^{-1}) up to 450, which is essential for efficient aerobic cultivation;
- Individually controlled bioreactor accelerates optimization process;
- Possibility to cultivate microaerophilic and obligate anaerobic microorganisms (not strict anaerobic conditions);
- Reverse-Spin® mixing principle enables non-invasive biomass measurement in real-time ;
- Near-infrared optical system makes it possible to register cell growth kinetics;
- Free of charge software for storage, demonstration and analysis of data in real-time;
- Compact design with a low profile and small footprint for personal application;
- Temperature control for bioprocess applications;
- Active cooling for rapid temperature control, e.g. for temperature fluctuation experiments;
- Task profiling for process automatization;
- Cloud data storage possibility to remotely monitor the process of cultivation while at home or using a mobile phone.

SOFTWARE FEATURES

- Real-Time cell growth logging;
- 3D graphical representation of OD or growth rate over time over unit;
- Pause option;
- Save/Load option;
- Report option: PDF and Excel;
- Connect up to 10 units (recommended) simultaneously to 1 computer;
- Remote monitoring option (requires internet connection);
- Cycling/Profiling options;
- User manual calibration possibility for most cells.

Cat. number

BS-010160-A04

BS-010158-A04

BS-010158-AK

BS-010158-CK

BS-010158-BK



RTS-1 and RTS-1C, Personal bioreactors

Recommendations for creating personal settings for cultivation of microorganisms.

Points that should be considered:

1. The growth rate directly depends on the tube's rotation speed, since it is directly proportional (in the range from 1,500 to 2,500 rpm) with the rate of saturation of the medium with oxygen.
2. Naturally, with aerobic metabolism, the change in OD over time will also proportionally increase depending on same as above.
3. This will also affect the specific growth rate $\Delta OD/\Delta t$.
4. As well as the time for the growth curve to reach the stationary growth phase during aerobic fermentation (the higher the tube rotation speed, the faster the culture's exit to the stationary phase)
5. The saturation of the medium with oxygen will depend on the frequency of switching the tube rotation to the opposite (RST) (the more often the direction of rotation of the tube is reversed, the higher the oxygen mass transfer)
6. OD $\lambda=850$ — this wavelength is used to measure microorganism cell concentration because nutrient media and microorganism cells have colour. This must be taken into account when monitoring the specific dynamics of microorganism growth. In order to go into the "shadow" region (independent of the colour of the medium and microorganism), we offer the near infrared (not visible to the human eye) light scattering measurement range of — 850 nm. Since we are still in the sensitive range and, at the same time, are independent of the natural colouration of microorganism colonies. The conversion factor OD_{850}/OD_{600} is about 2.

It is known that the aerobic bacterial growth is influenced by efficient gas exchange. Figure 2 a-c, serves as an example of growth optimization and illustrates the relationship between RST and gas exchange. As RST decreased the specific growth rate, and biomass yield increased, the highest aeration and optimal growth conditions for *E. coli BL21* optimized at 2000 RPM 1 s RST.

Cell growth depending on rotation intensity

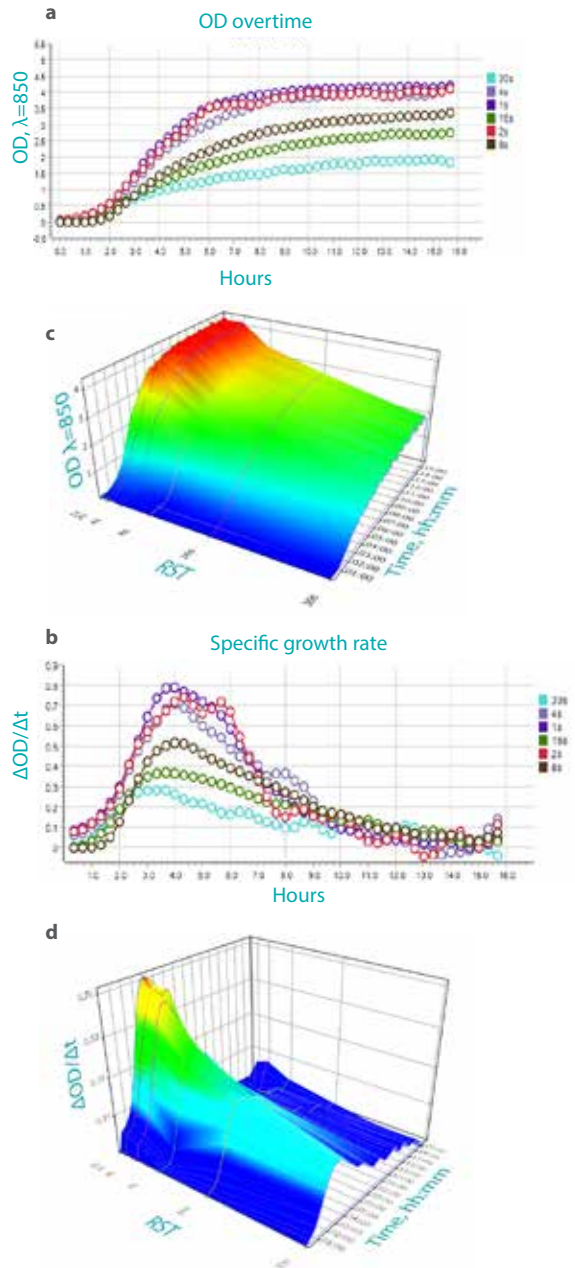


Figure 2, a-c. Influence of Reverse Spin Time (RST) on growth kinetics of *E. coli BL21* in OD_{600} . (a-c) Biomass growth; (b-d) Specific growth rate; throughout cultures were grown in 50 ml TPP Bioreactor tubes, 30% filling volume, 2,000 RPM, RST 1, 2, 4, 8, 16, 30 seconds, LB medium and 37°C temperature, to convert OD_{850} to OD_{600} simply multiply OD_{850} by 1.9.

RTS-1 and RTS-1C, Personal bioreactors

$k_{i,a}$ (h⁻¹) results in RTS-1/C

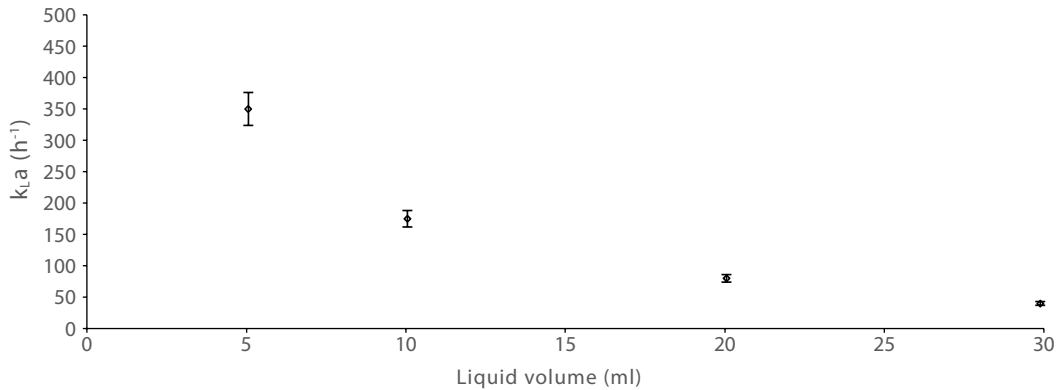


Figure 3. Determination of $k_{i,a}$ in 50 ml TPP Bioreactor tubes. The bioreactor vessels were filled with 5, 10, 20, 30 ml deionized water, and measurements were made by non-invasive O₂ sensors and optics (PreSens, Regensburg, Germany) at 37°C using the gassing-out method. Mean and standard deviation of at least five independent experiments are shown.

The $k_{i,a}$ was measured in 5, 10, 20, 30 ml of deionized water in 50 ml TPP Bioreactor tubes at agitation rate of 2,000 rpm and 1 s RST, this agitation rate was found optimal for Reverse-Spin® mixing principle during initial optimization studies. Over the working volume range, the $k_{i,a}$ increased with the decrease of liquid volume (Figure 3). At the smallest working volume of 5 ml, the highest $k_{i,a}$ of 350 ± 26 h⁻¹ was reached.

Cells successfully cultivated

Saccharomyces cerevisiae, *Pichia pastoris*, *Yarrowia lipolytica*, *Bacillus subtilis*, *Escherichia coli*, *Lactobacillus acidophilus*, *Bifidobacterium bifidum*, *Pseudomonas aeruginosa*, *Hybridoma*, *Jurkat* and CHO cells.

Types of recommended tubes

For aerobic microorganisms, it is recommended to use tubes that are supplied by TPP — TubeSpin® Bioreactor 50ml. For obtaining optimal results growing aerotolerant anaerobes, it is required to seal the screw cap of TPP TubeSpin® Bioreactor 50ml by tape or purchase TPP TubeSpin® 50 ml falcon tubes without the membrane filter. It is also possible to use other manufacturer tubes of the same type, e.g. Corning® 50 ml Mini Bioreactor, but the device rotor must be modified. It is possible to request this specific modification.

Factory calibration particle size and calibration coefficients 600nm/850nm

Factory calibration of the instrument is designed for rod-shaped bacteria size of *E.coli BL21*. In case of exceeding this size, the measurement system will not work correctly. Optical density OD₈₅₀ to OD₆₀₀ conversion coefficient of the factory calibration is equal to 1.9.

Factory calibration growth phase influence on measurement accuracy

During the growth transition of *Escherichia coli* culture from exponential growth to the stationary phase, many morphological and physiological changes occur, including cell volume decrease and cell shape change. Therefore, if cells were taken for referent measurement using a spectrophotometer at different stages from the stationary phase, then the correctness of measurement will be worse than specified.

Conversion rate coefficient of user calibration

Optical density OD₈₅₀ to OD₆₀₀ nm conversion rate coefficient depends on the cell size and volume. Therefore, the coefficient will be different for other cell sizes. The device can be calibrated at desired reference wavelength to meet user's needs, e.g. OD₆₀₀.

Do you want to test this system?

We can provide demo units for 50% of the price for testing or creating an application note. For such, inquiries please contact our R&D department directly at igor@biosan.lv.

RTS-8 and RTS-8 Plus, Multi-channel Bioreactors

DESCRIPTION

RTS-8 and **RTS-8 plus** are multi-channel bioreactors that utilise patented Reverse-Spin® technology that applies non-invasive, mechanically driven, low energy consumption, innovative type of agitation where cell suspension is mixed by the single-use falcon bioreactor tube rotation around its axis with a change of direction of rotation motion resulting in highly efficient mixing and oxygenation for aerobic cultivation.

Combined with a near-infrared measurement system, it is possible to register cell growth kinetics and additionally on **RTS-8 plus** fluorescence and luminescence measurement systems used to register pH and O₂ non-invasively in real-time. For pH and O₂, innovative single-use sensor spots are used inside the tubes. Although O₂ supply is one of the major issues in the cultivation of aerobic organisms, especially in oxygen-limited conditions, adequate methods for real monitoring of dissolved oxygen were missing, and sufficient O₂ supply was usually assumed. Innovative non-invasive oxygen sensors integrated into falcon tubes now enable online oxygen monitoring and give new insights into metabolic activities.

The pH is one of the major issues in the cultivation of cells, yeast or bacteria. Cultivation vessels, that are sensor-limited, are widely applied in academic and industrial bioprocess development. As adequate methods for real monitoring of pH were not available, cumbersome at-line sampling was used, lacking high data density and interfering with growth. Non-invasive real-time pH measurement provides new insights into metabolic activity and changes in metabolic pathways.

SOFTWARE FEATURES

- Real-Time cell growth logging;
- 3D graphical representation of OD or growth rate over time over unit;
- Pause option;
- Save/Load option;
- Report option: PDF and Excel;
- Remote monitoring option (requires internet connection);
- Cycling/Profiling options;
- User manual calibration possibility for most cells.



only RTS-8 Plus

- Real-Time pH and O₂ measurement and logging

Tube for RTS-8 Plus with sensor



Smart Plus
Product Class



USB
connection



Reverse-Spin®

Innovative Mixing Technology

FEATURES

- Parallel cultivation of 8 tube bioreactors enables to save time and resources for bioprocess optimization;
- Individually controlled bioreactor accelerates optimization process;
- Possibility to cultivate microaerophilic and obligate anaerobic microorganisms (not strict anaerobic conditions);
- Reverse-Spin® mixing principle enables non-invasive biomass measurement in real time;
- Near-infrared optical system makes it possible to register cell growth kinetics;
- Free of charge software for storage, demonstration and analysis of data in real time;
- Compact design with low profile and small footprint for personal application;
- Individual temperature control for bioprocess applications;
- Active cooling for rapid temperature control, e.g. for temperature fluctuation experiments;
- Task profiling for process automatization;
- Cloud data storage to remotely monitor the process of cultivation while at home or using a mobile phone.



only RTS-8 Plus

- Non-invasive O₂ and pH measurement allows for accurate monitoring of metabolic activities cultivation while at home or using a mobile phone

Advantages of the sensor spots:

- They are small;
- Their signal does not depend on the flow rate of the sample;
- They can be physically divided from the measuring system which allows a non-invasive measurement;
- They can be used in disposables;
- Therefore, they are ideally suited for the examination of small sample volumes, highly parallelized measurements in disposables, and biotechnological applications.

RTS-8 and RTS-8 Plus, Multi-channel Bioreactors

TYPICAL APPLICATIONS

- Fermentation real time growth kinetics;
- Clone candidate screening;
- Protein expression;
- Temperature stress and fluctuation experiments;
- Media screening and optimization;
- Growth characterization;
- Inhibition and toxicity tests;
- Strain quality control;
- Initial bioprocess optimization studies.



Product
video

SPECIFICATIONS

	RTS-8	RTS-8 Plus
Light source		Laser
Measurement wavelength (λ)		850 ±15 nm
Measurement range		0–100 OD ₆₀₀
<i>E. coli</i> factory calibration measurement range		0–50 OD ₆₀₀
<i>S. cerevisiae</i> factory calibration measurement range		0–75 OD ₆₀₀
Achievable user calibration measurement error (range 0.1–6 OD ₆₀₀)		±0.3
Achievable user calibration measurement error (range 6–50 OD ₆₀₀)		≤5%
Achievable user calibration measurement error (range 50–75 OD ₆₀₀)		≤10%
Measurement periodicity per hour		1–60 (increment 1 min)
Temperature setting range		+15°C ... +60°C
Temperature control range		+15°C below ambient ... +60°C (increment 0.1°C)
Temperature stability		±0.3°C
Sample temperature accuracy (20–37°C)		±1°C
Tube sockets		8
Sample working volume range		3–50 ml
Speed control range		150–2,700 rpm (increment 1 rpm)
Display		LCD
Dimensions (W×D×H)		350 × 690 × 300 mm
Weight		20 kg
Nominal operating voltage		AC 230 V, 50 Hz
Power consumption		3.15 A / 500 W
O ₂ sensor* / pH sensor**	— / —	+ / +

*O₂ sensor

Range	0–100%
Accuracy	±0.05% O ₂ at 0.2%, ±0.4% O ₂ at 20.9%
Drift	<0.03% O ₂ within 30 days
Temperature range	up to 40°C
Response time (t ₉₀)	<6 s
Storage stability	18 months

**pH sensor

Range	4.0–8.5 pH
Accuracy	±0.10 pH at pH 7
Drift	<0.005 pH per day
Temperature range	up to 40°C
Response time (t ₉₀)	<120 s
Storage stability	18 months

ORDERING INFORMATION

	Cat. number
RTS-8	
Calibration <i>E.coli</i>	BS-010168-A01
Optional calibration	
Calibration <i>S.Cerevisiae</i>	BS-010168-A09
Calibration <i>E.coli</i> and <i>S.Cerevisiae</i>	BS-010168-A10
Including TPP TubeSpin® Bioreactor vessels 50 ml, 20 pcs	

Optional accessories:

TubeSpin® Bioreactor 50, TPP®, 20 pcs.	
TubeSpin® Bioreactor 50, TPP®, 180 pcs.	
USB 2.0 Hub 10 × ports	
Sterile TPP TubeSpin® Bioreactor vessels, 50 ml, with pH and O ₂ sensors	

	Cat. number
RTS-8 Plus	
Calibration <i>E.coli</i>	BS-010170-A01
Optional calibration	
Calibration <i>S.Cerevisiae</i>	BS-010170-A08
Calibration <i>E.coli</i> and <i>S.Cerevisiae</i>	BS-010170-A11
Supports pH and pO ₂ measurement, including TPP TubeSpin® Bioreactor vessels 50 ml, 20pcs and sterile TPP TubeSpin® Bioreactor vessels 50 ml with pH and O ₂ sensors, 10 pcs	

BS-010158-AK
BS-010158-CK
BS-010158-BK
200001368



ES-20, Orbital Shaker–Incubator

DESCRIPTION

The **ES-20** is a compact bench-top Shaker-Incubator used for mixing of biological liquids and incubation and cultivation of biological liquids according to the operator set program.

Built-in microprocessor thermocontroller provides constant temperature control in the incubator chamber. Forced heated air circulation inside the transparent plexiglass chamber guarantees even temperature distribution. Dismountable construction makes transportation easy.

Orbital shaking is controlled by the digital tachometer (rpm) and Digital time setting regardless of the temperature. The unit is equipped with a direct-drive system, ensuring the most reliable, stable long-time operation (up to 30 day nights).

The **ES-20** is extremely easy to operate, with a very straightforward setup of temperature, speed and time, using the two-line set-up and status display, which clearly indicates both set and actual values for each of the three parameters.

DIFFERENT INTERCHANGEABLE PLATFORMS

ALLOW USING **ES-20** FOR:


- Growing cell cultures in flasks and other laboratory glassware
- Extracting tissue samples at physiological temperatures
- Other sample preparation processes

DESCRIPTION

Temperature setting range	+25°C ... +42°C
Speed control range	50–250 rpm
Temperature control range	5°C above ambient ... +42°C
Setting resolution	0.1°C; 10 rpm
Temperature stability	±0.5°C
Temperature accuracy	±0.5°C
Temperature uniformity	±0.5°C
Orbit	10 mm
Display	LCD, 2 × 16 signs
Digital time setting	1 min–96 h/non-stop (1 min increment)
Timer sound signal	yes
Plexiglas walls thickness	7 mm
Maximum load	2.5 kg
Overall dimensions (W×D×H)	340 × 340 × 435 mm
Dimensions of the inner chamber	305 × 260 × 250 mm
Weight	13.2 kg
Nominal operating voltage	230 V, 50/60 Hz or 120 V, 50/60 Hz
Power consumption	160 W (0.7 A)/ (230/120 V) 170 W (1.6 A)

Basic Plus
Product Class

Ø 10 mm
orbit

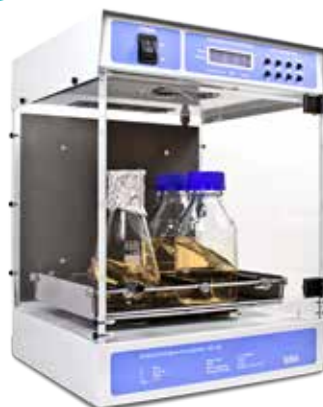

Product
video



Heat up time for **ES-20**:

from 25°C to 42°C

ORDERING INFORMATION:

Cat. number

ES-20 without platform BS-010111-AAA

Optional accessories:

Platforms:	UP-12	BS-010108-AK
	PP-4	BS-010108-BK
	P-12/100	BS-010108-EK
	P-6/250	BS-010108-DK
	P-16/88	BS-010116-BK

Description of all platforms
for **ES-20** on page 22

ES-20/60, Orbital Shaker–Incubator

DESCRIPTION

SPECIFICATIONS

Premium
Product Class

Ø 20 mm
orbit



Orbital Shaker-Incubator **ES-20/60** for biotechnological and pharmaceutical laboratories is a professional category equipment designed to cultivate microorganisms and eukaryotic cells, including animal, plant and insect cells. It is also possible to cultivate thermophilic bacteria in **ES-20/60** shaker-incubator.

Shaker is equipped with a direct-drive mechanism for platform motion. It provides a reliable and stable operation for the long term experiments needed for cell growth.

Shaker–Incubator **ES-20/60** provides smooth or intensive mixing in flasks installed on the platform.

Built-in noiseless thermoresistant brushless fan provides precise temperature distribution inside the chamber (adjustable for up to +80°C). The inner chamber is made of stainless steel. State-of-the-art motor, newest thermal insulation materials, soft-start of the platform motion and temperature **PID–control** decrease the energy consumption and make the Shaker-Incubator highly energy efficient despite its relatively large size.



Heat up time for **ES-20/60**

from 25°C + 90 min to 80°C



Temperature setting range	+25°C... +80°C
Speed control range	50–250 rpm
Temperature control range	10°C above ambient ... +80°C
Setting resolution	0.1°C; 10 rpm
Temperature stability	±0.5°C
Temperature accuracy	±0.5°C
Temperature uniformity	±0.5°C
Orbit	20 mm
Display	LCD, 2 × 16 signs
Digital time setting	1 min–96 h/non-stop (1 min increment)
Timer sound signal	yes
Maximum load	8 kg
Overall dimensions (W×D×H)	590 × 525 × 510 mm
Dimensions of the inner chamber	460 × 400 × 310 mm
Weight	41.1 kg
Nominal operating voltage	230 V, 50/60 Hz or 120 V, 50/60 Hz
Power consumption	450 W (2 A)/450 W (4.5 A)

ORDERING INFORMATION:

Cat. number



ES-20/60 without platform

BS-010135-AAA

Platforms cat. numbers for ES-20/60 can be found on page 23

ES-20/80 and ES-20/80C Orbital Shaker–Incubator

DESCRIPTION

ES-20/80, ES-20/80C shakers-incubators for biotechnological and pharmaceutical laboratories is a professional category equipment. The typical applications include microbial and cell culture cultivation, protein expression, solubility studies, general mixing, as well as other various applications in the fields of biology and chemistry. The unit is equipped with a newly developed triple eccentric mechanism for platform motion that provides supreme balancing characteristics, superior reliability and quiet operation. The achieved stability of the unit during vigorous mixing allows for stacking installation of up to 3 units which enables to save space.

The new display and easy to use user interface provide a clear and intuitive control of parameters and also allow data logging, storage and display over time. Additional features like out-of-balance sensor and automatic thermostat failure detection make this shaker-incubator an advanced and safe product. Bluetooth connectivity to PC allows for data management, data logging, parameter control and profiling in the included software.

A combination of an electric heating element and a built-in heat-resistant brushless fan provides precise temperature distribution inside the chamber up to +80 °C. In **ES-20/80C** shaker-incubator cooling of the chamber is provided by additional Peltier modules allowing to cool until 12.5 °C below the ambient temperature. Additionally, excellent sample temperature uniformity is achieved. The inner chamber is made of stainless steel. State-of-the-art motor, thermal insulation materials and parameter PID-control decrease the energy consumption and make the shaker-incubator highly energy efficient despite its relatively large size.

Stacking kit for 3× **ES-20/80, ES-20/80C**



Smart Plus
Product Class

Ø 20 mm
orbit



Bluetooth®
connection



50-400 rpm

Heat up time for **ES-20/80**

from 25°C + 75 min to 80°C

Smart Plus
Product Class

Ø 20 mm
orbit

NEW



Bluetooth®
connection



50-400 rpm

Heat up time for **ES-20/80C**

from 4°C + 70 min to 80°C

ES-20/80 and ES-20/80C, Orbital Shakers-Incubators

	ES-20/80	ES-20/80C
Temperature setting range	+25°C... +80°C	4°C... +80°C
Speed control range	50–400 rpm (increment 10 rpm)	
Temperature control range	5°C above ambient ... +80°C	12.5°C below ambient... +80°C
Setting resolution	0.1°C; 10 rpm	
Temperature stability	—	±0.1°C at 37°C
Temperature accuracy	±0.1°C at 37°C	
Temperature uniformity	±0.3°C at 37°C	±0.2°C at 37°C
Orbit	20 mm	
Display	TFT, 5 inches	
Digital time setting	1 min–96 h/non-stop (1 min increment)	
Timer sound signal	yes	
Maximum load	10.6 kg	
Data transfer	Bluetooth®	
Stacking	up to 3*	
PC system requirements	Intel/AMD Processor, 1 GB RAM, Windows Vista/7/8/8.1/10/11, USB, Bluetooth	
Overall dimensions (WxDxH)	620 × 530 × 510 mm	620 × 620 × 510 mm
Inner chamber dimensions (WxDxH)	460 × 400 × 310** mm	460 × 400 × 325** mm
Weight	48 kg	50 kg
Nominal operating voltage	230 V, 50/60 Hz	
Power consumption	500 W (2.2 A)	

* Additional stacking kit required

** Height is measured from platform surface

ORDERING INFORMATION

Cat. number 

ES-20/80 with software, without platform

BS-010135-AAA

ES-20/80C with software, without platform

BS-010173-A01

Optional accessories:

USB Bluetooth® adapter

BS-010425-FK

Stacking kit for 2× **ES-20/80**, **ES-20/80C**

BS-010167-OK

Stacking kit for 3× **ES-20/80**, **ES-20/80C**

BS-010167-PK

Platforms cat. numbers for **ES-20/80**, **ES-20/80C** can be found on next page



Platforms cat. numbers for **ES-20/80**, **ES-20/80C** can be found on next page

Platforms for ES-20/80 and ES-20/80C

Platform	Description	Dimensions (Working Area)	Cat. number
HSP-30/100 	Platform with 30 tight fit clamps for 100–150 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010167-KK
HSP-16/250 	Platform with 16 tight fit clamps for 250–300 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010167-MK
HSP-9/500 	Platform with 9 tight fit clamps for 500 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010167-NK
HSP-6/1000 	Platform with 6 tight fit clamps for 1,000 ml flasks	360 × 400 mm (360 × 400 mm)	BS-010167-LK
PP-400 	Flat platform with non-slip silicone mat	360 × 400 mm (360 × 400 mm)	BS-010135-FK
UP-168 	Universal platform for different flasks (Clamps ordered separately)	360 × 400 mm (360 × 400 mm)	BS-010135-JK
+			
HSC-50 HSC-100 HSC-250 HSC-500 HSC-1000 	Tight fit clamp for 50, 100, 250, 500, 1,000 ml flask (for UP-168)	Ø 50 mm Ø 65 mm Ø 85 mm Ø 105 mm Ø 130 mm	BS-010167-DK BS-010167-EK BS-010167-FK BS-010167-JK BS-010167-IK
+			
NEW SPML 	Set of 3 double-sided adhesive strips as an alternative for regular flask clamps (for UP-168)	390 × 80 × 3 mm	BS-010135-MK
+			
TR-21/50 	Test tube rack for 50 ml with 21 drillings (for UP-168)	340 × 124 mm (2 per platform)	BS-010135-KK
+			
TR-44/15 	Test tube rack for 15 ml with 44 drillings (for UP-168)	340 × 124 mm (2 per platform)	BS-010135-LK

Premium
Product Class



NEW CTR-6, CO₂ Tube Roller

CO₂ Tube Roller **CTR-6** provides regulated rocking and rolling of maximum up to 6 rollers and is designed for use specifically in CO₂ incubators. **CTR-6** is specifically designed for use in harsh environments such as CO₂ and humidity and provides reproducible results for cell culture growth. Possibility to remove rollers makes the unit flexible and allows for performing various procedures and techniques in various cultivation vessels. The specially designed remote controller allows for protection of electronics from CO₂ incubator environment, as well as, does not interfere with the experiment.

DESCRIPTION

SPECIFICATIONS

Speed control range*	5–80 RPM
Increment	1 RPM
Digital time setting	1 min – 96 h or non-stop
Increment	1 min
Maximum load	3 kg
Tilt angle	4°
Dimensions (WxDxH)	310 × 262 × 80 mm
Input current/power consumption	12 V, 415 mA/ 5 W
Weight, accurate within ±10%	3 kg
External power supply	Input AC 100–240 V; 50/60 Hz Output DC 12 V

* Maximum speed depends on the load on the platform and the shape of the vessels



ORDERING INFORMATION

CTR-6, CO₂ Tube Roller

Cat. number 

BS-010174-A01

S-Bt Smart Biotherm, Compact CO₂ Incubator

DESCRIPTION

S-Bt Smart Biotherm is designed for work in the areas of cell biology (operations with animal cell cultures and tissues), molecular biology (DNA/RNA reaction analysis, hybridization reactions), biotechnology (synthesis of target proteins and other molecules), immunology (synthesis of antibodies and other proteins of the immune system). Unit provides six-sided heating: the heating elements are located on the walls and on the door, thus providing excellent uniform temperature distribution, regardless of external factors, such as ambient temperature and positioning of the device.

Built-in infrared CO₂-sensor allows accurate control of the CO₂ level. The sensor makes measurement non-sensitive to changes in temperature and humidity inside the incubator.

The chamber is made of stainless steel with smoothed seams to minimize contamination and to facilitate cleaning.

S-Bt is equipped with a UV air recirculation system — 1 UV lamp and a fan are mounted behind the rear wall, providing decontamination of the working volume.

A convenient access port is built in the wall of the incubator for easy output of wire sensors or devices' installed inside. The access port is heated independently to prevent the formation of condensate.

Unit is equipped with error tracing and alarm systems, which significantly lower potential risks during operation.

Unit is equipped with a "black box" system that records temperature, humidity and CO₂ levels, as well as statuses for door opening, UV lamp, fan and errors, to the internal memory.

Bluetooth® connection to PC is available.

SPECIFICATIONS

Chamber Material	Stainless steel (1 mm)
Temperature setting range	+25°C ... +60°C
Temperature stability	±0.1°C
Temperature uniformity @37°C	±0.3°C
Timer sound signal	yes
Working volume	46 l
Number of shelves	3 (max 6)
Inner door	Glass
Relative humidity	>90% @ 37°C
Humidity delivery	Water bath
CO ₂ control range*	0–20%
CO ₂ sensor	Infrared sensor
Temperature and CO ₂ level input	Digital
UV lamp	1 × 6 W, TUV G6T5
Data transfer	Wireless
Access port	1 (Ø 26 mm)
Working voltage	230V, 50/60 Hz; 115 V, 50/60 Hz
Power consumption	600 W
Weight	37.7 kg
Dimensions (W×D×H)	500 × 560 × 550 mm
Inner chamber dimensions (W×D×H)	350 × 330 × 390 mm

* — At set temperature from ambient to 50 °C



Bluetooth®
connection



Product
video

APPLICATION AREAS

- **Cell biology:** operations with animal cell cultures and tissues;
- **Molecular biology:** DNA/RNA reaction analysis, hybridization reactions;
- **Biotechnology:** synthesis of target proteins and other molecules;
- **Immunology:** synthesis of antibodies and other proteins of immune system.

FEATURES

- Six-sided heating provides uniform distribution of the temperature inside the chamber;
- Infrared CO₂ sensor, non-sensitive to temperature and humidity changes;
- UV recirculation system for decontamination cycles;
- Bluetooth data transfer to PC;
- «Black box» parameter logging system;
- Error tracing and alarm system;
- Separately heated lockable port for chamber access for cables.

S-Bt Smart Biotherm, Compact CO₂ Incubator

Simple CO₂ tank connection




Air UV recirculation system in the chamber




Gas purification filter



PC software

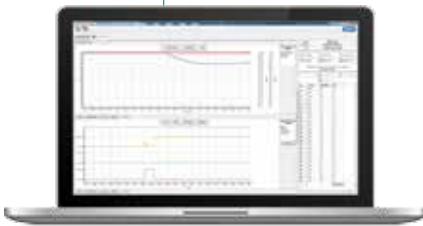

Incubator stacking device

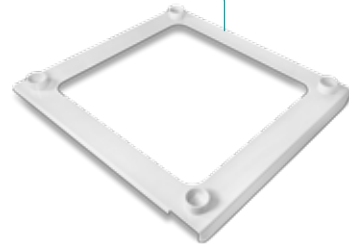

ORDERING INFORMATION

Cat. number

S-Bt Smart Biotherm, PC software included + RS6, rack with 3 shelves

BS-010425-A01

S-Bt Smart Biotherm, PC software included + RS2, rack for CPS-20/CTR-6 installation

BS-010425-A10

Optional accessories:

CPS-20, CO₂ Shaker

BS-010172-A01

CTR-6, CO₂ Tube Roller

BS-010174-A01

Shelf

BS-010425-AK

USB Bluetooth® adapter

BS-010425-FK

Incubator stacking device

BS-010425-CK

CPS-20, CO₂ Shaker NEW

DESCRIPTION

CO₂ Shaker **CPS-20** provides regulated orbital motion of the platform and is designed for use specifically in CO₂ incubators. **CPS-20** is specifically designed for unnecessary harsh environments such as CO₂ and humidity and provides reproducible results for cell culture growth. A choice of five interchangeable platforms provides the possibility of performing various procedures and techniques in various cultivation vessels.

CPS-20 incorporates a brushless motor with a guaranteed service life of up to 35,000 hours and a direct-drive system, ensuring the most reliable, stable long-time operation (up to 30 day nights). The unit is equipped with a triple eccentric mechanism for platform motion that provides supreme balancing characteristics, superior reliability and quiet operation. The specially designed remote controller allows for the protection of electronics from a CO₂ incubator environment, as well as, the remote control minimizes interference with the incubator environment and the ongoing experiment.

SPECIFICATIONS

Speed control range	50–250 rpm (increment 10 rpm) max. speed depends on the load and vessels' shape
Digital time setting	1 min–96 h/non-stop (increment 1 min)
Digital speed control	+
Maximum continuous operation time	168 h
Orbit	20 mm
Maximum load	3 kg
Overall dimensions (WxDxH)	255 × 255 × 100 mm
Weight	3.4 kg
Input current/power consumption	470 mA / 5.7 W
External power supply	Input AC 100–240 V; 50/60 Hz Output DC 12 V



RS2, rack for CPS-20 installation



Platforms for CPS-20

Platform	Description	Dimensions (mm)	Working Area (mm)	Cat. number
UP-12 	Universal platform with adjustable bars for different types of flasks, bottles and beakers with silicone mat	285 × 220 × 40	270 × 185 × 40	BS-010108-AK
Bio PP-4 	Flat platform with silicone mat for Petri dishes, culture flasks, agglutination cards	255 × 255	230 × 230	BS-010116-AK
P-12/100 	Platform with clamps for flasks, 100–150 ml (12 places)	250 × 190	250 × 190	BS-010108-EK
P-6/250 	Platform with clamps for flasks, 250–300 ml (6 places)	250 × 190	250 × 190	BS-010108-DK
P-16/88 	Platform with spring holders for up to 88 tubes up to 30 mm diameter (e. g. 10 ml, 15 ml, 50 ml tubes)	275 × 205 × 75	275 × 205 × 75	BS-010116-BK

S-Bt Smart Biotherm, Compact CO₂ Incubator

ORDERING INFORMATION

CPS-20, CO₂ Shaker

RS2, rack for CPS-20/CTR-6 installation

Cat. number

BS-010172-A01

BS-010425-HK

LAB DIAGNOSTICS:

DNA/RNA PURIFICATION, IMMUNODIAGNOSTICS



3D-IW8
Inteliwasher



PST-60HL
Plate Shaker-Thermostat



HiPo MPP-96
Microplate Photometer

Automatic DNA/RNA extraction

Liquid Handling

Assist Series
Pipettes



Incubation

TS-100



TS-100C



Sample resuspension
and droplet spin down

FVL-2400N



Automatic
extraction

BioMagPure 12 Plus
and reagents



Manual DNA/RNA extraction using magnetic beads technology

Sample resuspension and droplet spin down

FVL-2400N



Sample preparation in UV-Cabinet for PCR

UVC/T-M-AR, or similar, see UV-Cabinets for PCR



Mixing and resuspension

MPS-1



V-1 plus



Multi Bio RS-24



Capture of magnetic beads

MagSorb-16



Centrifugation

Microspin 12



Vacuum aspiration

FTA-1



FTA-2i



Incubation

TDB-120



TS-100C



Reagents

Reagents for extraction*



* — Information about current offers on the products of other manufacturers are available in the corresponding sections of our site www.biosan.lv/en/products

BioMagPure 12 Plus, Compact Bench-Top Robotic Workstation For Automated Nucleic Acid Purification

DESCRIPTION

The **BioMagPure 12 Plus** consists of a compact bench-top robotic workstation for automated nucleic acid purification. Usage of pre-filled reagent cartridges and disposable consumables enable a true walk-away automation and high-quality nucleic acid extraction solution. Proven magnetic separation technology makes purification efficient, easy to use, reliable, safe and cost-effective.

BioMagPure 12 Plus has an ingeniously designed polygonal reaction chamber with patented parts that ensure high efficiencies of lysis and elution through a large contact area of magnet and heating element allowing maximisation of magnetic bead recovery, minimisation of the residues of magnetic beads and alcohols in the final elute product. Specific formation of reaction chamber ensures unrivalled mixing ability and exclude conventional mixing by tip or pipetting thus eliminates cross-contamination possibility.

Reagent kit contains everything for extraction procedure performance, including all necessary plastics, pre-filled reagent cartridges, incubation buffers and solutions for sample pre-treatment (if needed),

With the flexibility of processing 1–12 samples per run, the **BioMagPure 12 Plus** is tailor-made to fit small clinics and early-stage laboratories. By occupying minimal counter space and greatly reducing technician man-hours, this series allows organizations to operate facilities in a much more cost-effective fashion.



Product video


3 easy steps

- 1

LOAD
- 2

RUN
- 3

OBTAIN

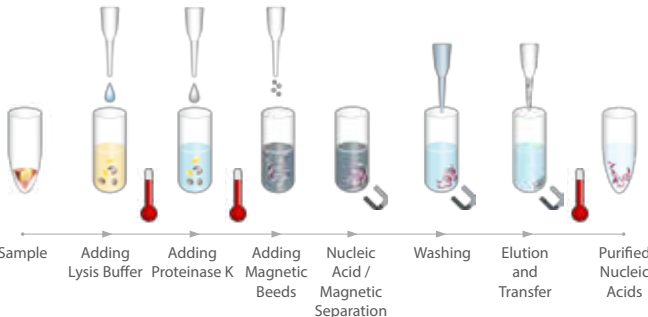

SPECIFICATIONS

Processing time	45–60 min
Processing capability	1–12 samples per run
Extraction technology	magnetic particle separation technology
Protocol	programmed by scanning a barcode
Protocol input	barcode scanner
Sample volume	10–2,000 µl (depending on the kit)
Elution volume	50–400 µl
Connection to PC	not required
Display	LCD (20 × 4)
Certification	CE IVD
Nominal operating voltage	110–240 V, 50/60 Hz
Dimensions (W×D×H)	560 × 590 × 510 mm
Weight	55 kg

Features:

- Advanced magnetic bead technology;
- Reaction chamber with patented parts;
- Piercing-pin system for elimination of cross-contamination;
- Walk-away automation;
- Reliable quality;
- No PC required;
- Ready-to-use reagent cartridges;
- 3 easy steps: LOAD-RUN-OBTAIN.

Reagents for BioMagPure 12 Plus



ORDERING INFORMATION:

Name	Description	Cat. number
BioMagPure 12 Plus	Compact Bench-Top Robotic Workstation For Automated Nucleic Acid Purification.	BS-060202-AAA
Blood DNA Extraction Kit 200	Blood DNA Extraction Kit is used with the BioMagPure 12 Plus instrument to extract DNA from 10-400µl mammalian whole blood, suspension of mammalian blood cells.	BS-060201-AK
Blood DNA Extraction Kit 1200	Blood DNA Extraction Kit is used with the BioMagPure 12 Plus instruments to extract gDNA from 400-1000µl mammalian blood, suspension of mammalian blood cells.	BS-060201-BK
Viral Nucleic Acid Extraction Kit	Viral Nucleic Acid Extraction Kit is used with the BioMagPure 12 Plus instrument to extract Viral DNA or RNA from human biological specimens such as serum, plasma, and other cell-free fluids.	BS-060201-CK
Tissue DNA Extraction Kit	BioMagPure 12 Plus Tissue DNA Extraction Kit is used with the BioMagPure 12 Plus instrument to extract genomic DNA from a variety of animal tissues, swab samples and bloodstain.	BS-060201-DK
Cultured Cell DNA Extraction Kit	Cultured Cell DNA Extraction Kit is used with the BioMagPure 12 Plus instrument to extract genomic DNA from culture cells and buffy coat.	BS-060201-EK
Bacterial DNA Extraction Kit	Bacterial DNA Extraction Kit is used with the BioMagPure 12 Plus instrument to extract genomic DNA from both Gram-positive and Gram-negative bacteria.	BS-060201-FK
HPV DNA Extraction Kit for Swab	HPV DNA Extraction Kit is used with the BioMagPure 12 Plus instrument for DNA extraction of the Human Papillomavirus (HPV) from cervical cell samples which collected by cervical brush or genital swab in liquid-based Medium (e.g. Hologic Thinprep PreservCyt®, BD Surepath™, etc.) or other STM (sample transport media) preservation solutions(e.g. QIAGEN DNA PAP Cervical sampler, Roche Cobas® PCR Cell Collection Media, HybriBio cell preservation solution, etc.).	BS-060201-GK
TB DNA Extraction Kit	TB DNA Extraction Kit is used with the BioMagPure 12 Plus instrument for extraction of genomic DNA of <i>Mycobacteria spp.</i> (e.g. <i>Mycobacterium tuberculosis</i>) from different specimen.	BS-060201-IK
FFPE DNA Extraction Kit	FFPE DNA Extraction Kit is used with the BioMagPure 12 Plus instrument for extraction of genomic DNA from FFPE (Formalin-Fixed, Paraffin-Embedded) tissue samples. Providing good quality, high integrity DNA for Molecular diagnosis and research works	BS-060201-JK
Forensic DNA Extraction Kit	Forensic DNA extraction kit is used to extract and isolate genomic DNA from forensic samples.	BS-060201-KK
Viral/Pathogen Nucleic Acids Extraction Kit A	Viral/Pathogen Nucleic Acids Extraction Kit A is used with the BioMagPure 12 Plus instrument to extract Viral and bacterial DNA/RNA from cell-free samples, such as serum, plasma, and other cell-free body fluids.	BS-060201-LK
Viral/Pathogen Nucleic Acids Extraction Kit B	Viral/Pathogen Nucleic Acids Extraction Kit B is used with the BioMagPure 12 Plus instrument to extract viral and bacterial DNA/RNA from swab samples (cell-rich samples).	BS-060201-MK
Viral RNA Extraction Kit	Viral Nucleic RNA Extraction Kit is used with the BioMagPure 12 Plus instrument to extract Viral RNA from human biological specimens such as serum, plasma, and other cell-free fluids.	BS-060201-NK
Plant DNA Extraction Kit	Plant DNA Extraction Kit is used with the BioMagPure 12 Plus instrument to extract genomic DNA from plants (leaf, seeds and spores) and fungal tissues. Up to 100 mg of tissue can be used for purification.	BS-060201-OK
Total RNA Extraction Kit	Total RNA Extraction Kit is used with the BioMagPure 12 Plus instrument to extract total RNA from whole blood, blood cells, animal tissue, plant tissue, yeast or cultured cells.	BS-060201-PK
Viral Nucleic Acid Large Volume Extraction Kit	Viral Nucleic Acid Large Volume Extraction Kit is used with the BioMagPure 12 Plus instrument to extract Viral DNA or RNA from human biological specimens such as serum, plasma, and other cell-free fluids.	BS-060201-QK
CFC DNA Extraction Kit Large Volume	CFC DNA Extraction Kit Large Volume - is used with the BioMagPure 12 Plus instrument to extract circulating DNA from plasma serum or cell-free body fluids sample volume ranging up to 5 ml.	BS-060201-RK


MagSorb-16, Magnetic rack for manual nucleic acid extraction

DESCRIPTION

MagSorb-16 is a magnetic rack designed to isolate nucleic acids utilizing magnetic particles. High-performance permanent neodymium magnets are built into the rack.

MagSorb-16, consists of two convenient quick-detachable tube holders that hold up to 16 tubes (1.5–2 ml) and a magnetic stand on which the tube holders are placed.

Magnetic beads based NA extraction methods offer quick and efficient separation of genetic material from cellular leftovers. Wide variety of commercially available extraction kits manufactured by companies like ThermoFisher Scientific, Promega, Qiagen etc. are available in the market. Our **MagSorb-16** magnetic rack can be used together with any manufacturers magnetic beads based NA extraction reagents that allow working in 1.5–2.0 ml tube format.

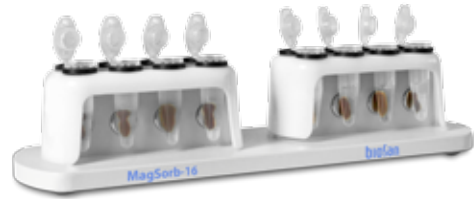

SPECIFICATIONS

Number of places in stand	16
Tube's volume	1.5–2 ml
Type of tubes	Eppendorf or equivalent

ORDERING INFORMATION

MagSorb-16, magnetic rack

Cat. number

BS-010601



PCR Analysis

Sample
resuspension and
droplet spin down



FVL-2400N
(tubes)



CVP-2
(PCR-plates)

Real time
PCR detection



BioQuant-96

Reagents



abTes COVID-19 qPCR kit
(Reagents are CE-IVD certified — among WHO listed kits.)

BioQuant-96, Real-time PCR detection system



Product
video



DESCRIPTION

BioQuant-96 is the newest product of Biosan Molecular diagnostic product family.

It has adopted innovative thermoelectric refrigeration technology, brand-new light source and light path design. Detection from the top allows using different consumables – 0,2 tubes, 8-tube strips and semi-skirted and non-skirted 96-well PCR plates. The unique constant current power and 6-zone independent temperature control method ensure more rapid, correct and stable fluorescence quantitative analysis, while maintaining its excellent performance in lowest possible energy consumption.

Device is available in 5-channel and 6-channel configuration. Meanwhile, it has been added with functions including independent temperature control, low temperature storage of sample at 4°C and FAST mode for more faster cycling (confirm reagent compatibility with fast mode).

BioQuant-96 is comprehensively realizing automatic gain setting and improving user experience. It will fully meet the demand of scientific research laboratories and registered as I (A) class Medical device in EU it also will fully meet the requirements of any diagnostic laboratory.



SPECIFICATIONS

Sample capacity	0.2 ml single tube (transparent cap), 96 × 0.2 (0.1) ml Plate (transparent cap), 8-strip tubes (transparent cap)
Dynamic range	1~10 ¹⁰ copies
Excitation wavelength	300–800 nm
Emission wavelength	500–800 nm
Channels	5
Temperature setting range	4~105°C (increment: 0.1°C)
Heating/cooling rate	6.0°C/s / 5.5°C/s
Temperature accuracy	≤±0.1°C
Temperature fluctuation	≤±0.3°C
Temperature uniformity	≤±0.3°C
Temperature control mode	BLOCK/Tube simulation mode (automatic control based on sample volume)
Sample volume	5~100 µl
Gradient temperature range	1~36°C
Hot-lid temperature range	30~110°C (adjustable default: 105°C, Automatic Hot-lid)
Repeatability	5%
Scan mode	Entire plate or designated line
Program	Max. 20 segments for each program, max. 99 cycles
Continuous operation	+
Scan period	5.5 s
Software	BioQuant-96 Software included
Operating system	Windows 7/8/10/11
Minimum PC requirements	Memory: 512 MB, HD: 10 GB, CPU: Pentium® 4, Virtual Memory: ≥1,000 MB
Nominal operating voltage	100–240 V; ~50/60 Hz
Power consumption	600 W
Dimensions (WxDxH)	490 × 290 × 391 mm
Number of sockets	USB adapter, Bluetooth adapter
Certification	CE-IVD
Interface	Built-in Touchscreen

ORDERING INFORMATION

BioQuant-96

Cat. number

BS-050110



abTes COVID-19 qPCR kit, Reagents

COVID-19 is a single-stranded, positive-sense RNA virus that caused a nationwide outbreak of serious respiratory problems across the world. This RNA virus is capable of person-to-person transmission, spreading very rapidly and affecting almost every country. It is genetically highly related to SARS and MERS coronaviruses.

Together with our partners, AITbiotech, we are participating in the global fight against COVID-19, providing reliable Real-Time PCR kits to detect nCoV-2019 virus in clinical samples.

Kit is CE-IVD certified
and is among WHO listed kits.

The **abTES COVID-19 qPCR I Kit** is a qualitative real-time polymerase chain reaction (qPCR) kit that enables simultaneous detection of two COVID-19 — specific signature regions from its non-structure polypeptide (orf1a) in a single reaction. It also includes detection of the human housekeeping gene, GAPDH, as an Internal Control (IC) to identify possible PCR inhibitions from sample processing.


The kit contains all the necessary PCR reagents for rapid, sensitive, and specific detection using target-specific primers and double-labelled hydrolysis probes. This kit has been validated on samples extracted from sputum, nasopharyngeal and throat swabs.

Amount per kit — **100 tests**

Sensitivity: 2.2 copies/ μ l or 11 copies/rxn

Validated sample types (as recommended by CDC):

- Nasopharyngeal wash/aspirate or nasal aspirate;
- Nasopharyngeal & oropharyngeal swabs;
- Bronchoalveolar lavage;
- Tracheal aspirate;
- Sputum;
- Serum.

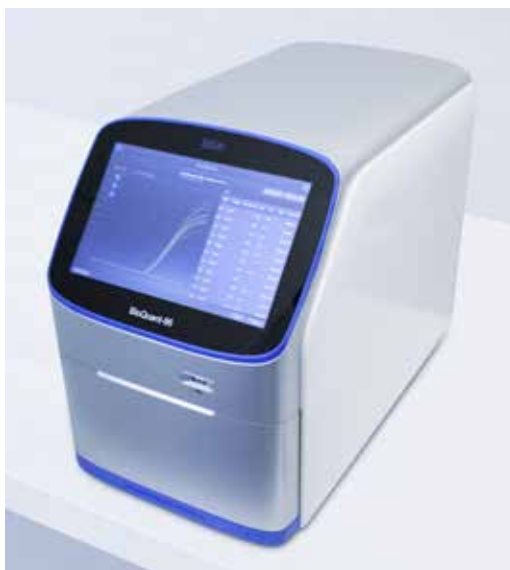

Enzyme-linked immunosorbent assay (ELISA)

Sample preparation
in laminar flow cabinet



Biological Safety Cabinet class II



PST-60HL



PST-60HL-4



PSU-2T

Incubation



IW-8



3D-IW8

Washing (Automated)



FTA-1 with MA-8



FTA-2i

Washing (Manual)



HiPo MPP-96
Microplate Photometer
with QuantAssay software



Reagents of ELISA*

Reading and Analysis

* — Information about current offers on the products of other manufacturers are available in the corresponding sections of our site www.biosan.lv/en/products

IW-8, IW-8m Intelispeed Washers



DESCRIPTION


Intelispeed Washer IW-8/IW-8m is designed to wash standard flat-bottom (two point aspiration) and U-shape (only in single point aspiration) 96 well plates and microstrips. The unit is fully programmable, ensuring multi-step solution ripening, aspiration (aspiration, combination of aspiration/liquid dispensing and soaking, as well as soaking cycle during a particular time).

The unit has 100 user-defined programs. Standard version is supplied with an 8-channel washing head for dispensing/aspiration, three bottles for washing and rinsing solutions, a waste bottle and bottle with filter. Optional 4-channel washing solution weight logger, **4 CHW Logger** is available.

The unit is designed for washing standard 96-well plates during analyses.

The unit provides:

- Washing mode;
- Rinsing mode;
- Mixing mode;
- Single point, two point aspiration;
- Possibility of additional solution mixing during time gap between two work cycles;
- Possibility to use microtest plates by different manufacturers, ensured by automated plate set up (adjusting to different depths of plate wells);
- Plate and strip washing mode;
- User-defined programs with adjustable parameters;
- Saving work programs.

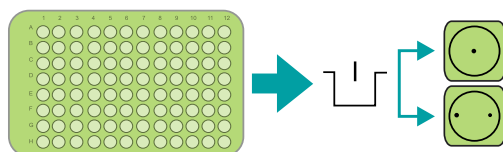

ORDERING INFORMATION:

IW-8

IW-8m IVD version

4 CHW Logger

Cat. number

BS-060106-AA1

BS-060108-A01

BS-060102-AK

IW-8, Intelispeed Washer



4-channel washing solution weight logger, **4 CHW Logger** provides automatic control of rinsing solutions and waste volume. The washer shows the remaining volume for each bottle as a percentage and gives a warning message in case of low solution volume or full waste bottle when **4 CHW Logger** is connected.

**4 CHW Logger Specifications:**

Max. loading per scale cup	2 kg
Dimensions (WxDxH)	267 × 252 × 97 mm
Weight	3 kg

Choice of 3 washing liquid bottles	
Minimum dispense volume	25 µl
Maximum dispense volume	1,600 µl
Dispense increment	25 µl
Dispensing accuracy	±2.5%
Allowed residual liquid volume not more than 2 µl in plate well	
Number of wells washed simultaneously	8
Number of washing cycles for each channel	1–15
Timer sound signal	yes
Aspiration time	0.2–3 s
Aspiration/dispensing speed	3 levels
Max. number of channels in a program	2
Soaking time	0–300 s (increment 10 s)
Shaking time	0–150 s (increment 5 s)
Number of washed rows	1–12
Time of plate single wash (350 µl), not more	45 s
Number of programs	101
Plate platform and washing head movement	automated
Indication of operation modes	8-line LCD
Dimensions (WxDxH)	375 × 345 × 180 mm
Weight with accessories	9.6 kg
External power supply	DC 12 V, 5 A
Consumed power	22 W

The unit is designed for use in closed laboratory rooms at temperatures from +4 to +40°C and relative humidity up to 80% at +31°C decreasing linearly to 50% relative humidity at 40°C

3D-IW8, Inteliwasher



4 CHW Logger

3D-IW8



Product
video

DESCRIPTION

Inteliwasher **3D-IW8** series microplate washer is designed to wash various types of standard 96-well microtitre plates, microstrips as well as microarrays on FastFRAME (rectangular well shape). It is suitable for washing wells with different bottom shapes: flat, U-shape and V-shape. The unit is fully programmable, ensuring multi-step solution ripening, aspiration (aspiration, combination of aspiration/liquid dispensing and soaking, as well as soaking cycle during a particular time). Dispense system of liquid dosage for each channel separately.

The unit provides:

- Washing mode;
- Rinsing mode;
- Mixing mode;
- Single point, two point, circular (circle or rectangular path) aspiration;
- Possibility of additional solution mixing during time gap between two work cycles;
- Possibility to use microtest plates by different manufacturers, ensured by automated plate set up (adjusting to different depths of plate wells);
- Round-bottom plate and strip washing mode;
- Possibility of user-defined programs with adjustable parameters.



ORDERING INFORMATION:

Cat. number

3D-IW8


BS-060102-AAI

4 CHW Logger

BS-060102-AK

3D-IW8, Inteliwasher


SPECIFICATIONS

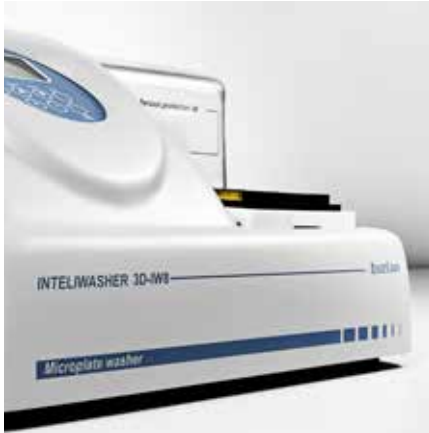

The unit has 50 programs divided into 5 following aspiration categories (see figure below):

- 1 **Type 1** (1.0–1.9) **IPF96 U/V** is intended for round and V-shape immunoplates, 1 point aspiration.
- 2 **Type 2** (2.0–2.9) **IPF96 FLAT-2** is intended for flat-bottom shape immunoplates, 2 point aspiration.
- 3 **Type 3** (3.0–3.9) **IPF96 FLAT-C** is intended for rectangular shape immunoplates, full-circle aspiration direction.
- 4 **Type 4** (4.0–4.9) **FastFRAME-2** is intended for multi-slide plate* with rectangular wells, 2 point aspiration.
- 5 **Type 5** (5.0–5.9) **FastFRAME-C** is intended for multi-slide* plate with rectangular wells, full-square aspiration direction.

* — The **Fas tFRAME** multi-slide plate or analog plate of another manufacturer, that is compatible with standard 25 × 76 mm (1 × 3 inch) glass slides.

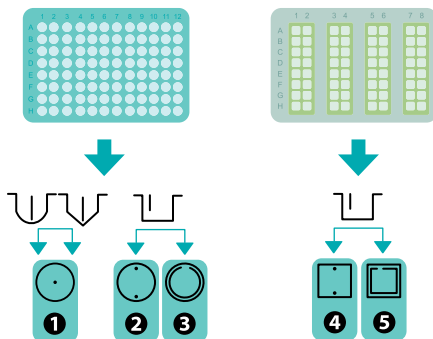

Minimum dispense volume	25 µl
Maximum dispense volume	1,600 µl
Dispense increment	25 µl
Dispensing accuracy	±2.5%
Allowed residual liquid volume in plate well, not more	2 µl
Number of wells washed simultaneously	8
Number of washing cycles	1–15
Timer sound signal	yes
Aspiration time	1–3 s
Final aspiration time	1–3 s
Aspiration/dispensing speed	3 levels
Max. number of channels in a program	2
Choice of 3 washing liquid bottles	
Soaking time	0–300 s (increment 10 s)
Shaking time	0–150 s (increment 5 s)
Number of washed rows	1–12
Time of one plate wash (300 µl), not more	45 s
Number of programs	50
Plate platform and washing head movement	automated
Indication of operation modes	LCD, 8-line
Dimensions (W×D×H)	375 × 345 × 180 mm
Weight with accessories	9.9 kg
External power supply	Input AC 100–240 V 50/60 Hz, Output DC 12 V
Input current/ power consumption	12 V, 1.8 A / 22 W

The unit is designed for use in closed laboratory rooms at temperatures from +4°C to +40°C and relative humidity up to 80% at +31°C decreasing linearly to 50% relative humidity at 40°C.

4-channel washing solution weight logger, **4 CHW Logger**, provides automatic control of rinsing solution and waste volumes. The washer shows the volume for each bottle as a percentage and gives a warning message in case of low solution volume or full waste bottle when **4 CHW Logger** is connected.

4 CHW LOGGER SPECIFICATIONS:

Max. loading per scale cup	2 kg
Dimensions	267 × 252 × 97 mm
Weight	3 kg

ImagerBio C, Combined Colorimetric MicroArray and ELISPOT Imager NEW



DESCRIPTION

ImagerBio C is a compact and portable instrument designed for imaging colorimetric Microarrays and Elispots. The instrument by standard operates in bright-field and darkfield illumination modes with bottom and top illuminations respectively. The instrument is able to work with 96-well plate, a 12x8 well strip, as well as microscopic glasses (placed in a special adapter). The device is equipped with a highly sensitive CMOS camera. The instrument has an automatic mechanism of 96 well plate ejection outside the body.

Unit uses various autofocus methods.

The instrument can be operated by a 20V portable lithium-ion battery, which makes the device usable in the field studies or point-of-care.

The compact design allows the device to be carried in the hand luggage of the aircraft.


SOFTWARE FEATURES:

- Working in Microarray or Elispot interface
- End-User interface, with simplified functional, where user just chooses assay, loads samples and gets results. Visual validation step is optional
- Automated image acquisition, analysis and report generation (up to 12x12 array in 96 well plate in under 3 minutes)
- Automatic array finder via machine learning and image recognition
- Grid lay-outting (Manual and Automatic)
- Password protected Assay Developers interface with full access to the vast software parameters
- Analysis of images by the average/median intensity of the spots
- Qualitative/quantitative analysis of the arrays
- Creating Qualitative/quantitative analysis assays
- Quantitative assays with 4/5 parameter logistics functions, etc.
- Setting multi level interpretation thresholds for different type of samples in the same well (e.g. tolerance to egg and lettuce)
- Reports available in PDF, CSV, EXCEL files
- Control of camera exposure, gain, XYZ kinematics.


Analysis Report

Experiment information:

Sample: AJ
 Well name: A20 0
 Assay: Covid19 beta2.1 avg DF in cell
 Template: OVID.1.2

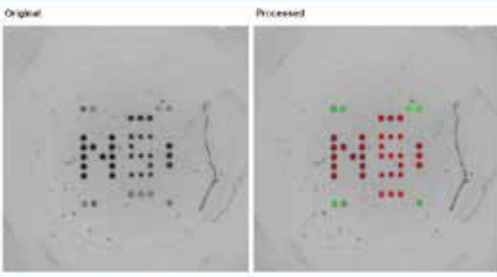

Test name	Result	ResultC
13p4xRBC	Positive	47.4
25p4x2	Positive	77.3
30p4xNTD	Positive	74.3
43p4x1	Positive	70.4
5h4xocapico	Positive	83.3
Test Diagnosis:		Result
Covid19 beta2.1 avg DF in cell		Positive

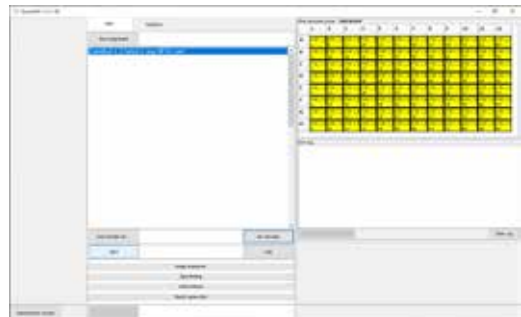

ImagerBio C, Combined Colorimetric MicroArray and ELISPOT Imager

SPECIFICATIONS

Illumination	Channels Brightfield, Darkfield (bottom, top illumination), Chemiluminescence
Arbitrary units measurement range	0 to 65,535
NIST Certified OD measurement range (for brightfield)	0.1 to 2.0 OD
NIST Certified Diffused reflectance meas. range (for darkfield)	2 to 99
Plates and Vessels	96 well plate / 12 × 8 well strip / 4 microscope slides
Light sources	LED
Lifetime of the light source	>10,000 hours
Data interface for unit controls / camera	USB 2.0 / USB 3.0
Camera	CMOS
Standard image resolution	1280×1280 pixels
Standard Sample Resolution	6 μm per pixel
Per Request Sample Resolution	5–7 μm per pixel
Image formats	png or tiff, 16 bit
Focus	Manual, Automatic, adjustable via PC
Software	Included
PC requirements (recommended to use NVIDIA GPU for microarrays, only for Elispot not necessary)	CPU: Intel i7, RAM: 8 GB Video card: Nvidia GTX 1050 Ti 4GB, or better (Capability only with Nvidia cards) SSD: 256 GB, OS: Windows 10/11 (64 bit)
Size (W × D × H)	330×345×150 mm
Weight, w/o power supply	not more than 6 kg
Power supply	Input AC 100–240 V 50/60 Hz, Output DC 20–24 V, 2.5A

Software developers interface, Assay Editor interface and Analytes template in well



ORDERING INFORMATION

ImagerBio C, Combined Colorimetric MicroArray and ELISPOT Imager

Cat. number

BSM000101-A01



ImagerBio F, Fluorescent MicroArray Imager NEW




DESCRIPTION

ImagerBio F is a compact and portable fluorescent imager designed for imaging fluorescent microarrays. The instrument by standard operates on 2 channels: CY3 and CY5. The instrument is able to work with 96-well plate, a 12x8 well strip, as well as microscopic glasses (placed in a special adapter). The device takes images from the below, so the wells must have a transparent bottom. The device is equipped with a highly sensitive CMOS camera capable of taking more sensitive images in a shorter exposure times, which saves operating time. Image filtering goes through fluorescent filters. The instrument has an automatic mechanism of tray ejection outside the body.

The instrument can operate with a 20V portable lithium-ion battery, which makes the device usable in the field studies or point-of-care.

The device is available in 2, 3 or 4 channel configuration: DAPI, FITC, CY3, CY5. Qdot options are also available.


The compact design of 2/3 channel version allows the device to be carried in the hand luggage of the aircraft.

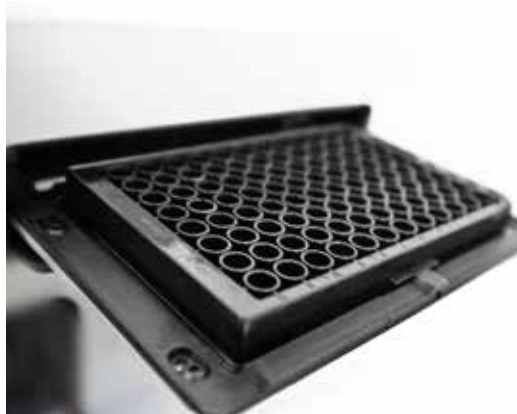

SOFTWARE FEATURES:

- Automatic array finder via machine learning and image recognition
- Grid lay-outing (Manual and Automatic)
- Analysis of images by the average/median intensity of the spots
- Qualitative/quantitative analysis of the arrays
- Creating Qualitative/quantitative analysis assays
- Quantitative assays with 4/5 parameter logistics functions, etc.
- Setting multi level interpretation thresholds for different type of samples in the same well (e.g. tolerance to egg and lettuce)
- Report exports to PDF, CSV, EXCEL
- Control of camera exposure, gain, XYZ kinematics.

Color overlap of CY3 and CY5 channels

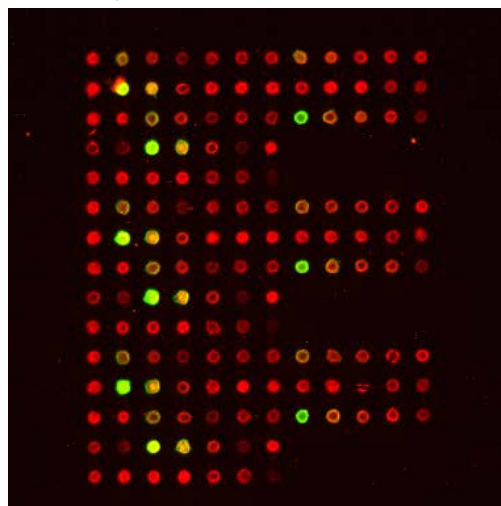

ImagerBio F, Fluorescent MicroArray Imager

SPECIFICATIONS

Fluorescence channels	CY3, CY5
Optional Fluorescence channels	DAPI, FITC, CY3, CY5, Qdots
Arbitrary units measurement range	0 to 65,535
Vessels	96 well plate / 12x8 well strip / 4 microscope slides
Light source	Laser
Lifetime of the light source	>10,000 hours
Data interface for unit controls/camera	USB 2.0 / USB 3.0
Camera	3 MP, CMOS
Standard image resolution	1280 x 1280 pixels
Resolution	6 µm per pixel, 5–7 per user request
Data interface for camera control	USB 3.0
Image formats	png/tiff 16 Bit, or other on request
Focus	Manual, adjustable via PC
Exposure	Controllable, up to ~10 s
Software	Include
PC requirements	CPU: Intel i7, RAM: 8 GB Video card: Nvidia GTX 1050 Ti 4GB, or better (Capability only with Nvidia cards) SSD: 256 GB, OS: Windows 10/11 (64 bit)
Size (2/3 channel) (W x D x H)	330x345x150 mm
Weight (2 channel), w/o power supply	not more than 8 kg
Power supply	Input AC 100–240 V 50/60 Hz, Output DC 20–24 V, 2.5A

Software developers interface, Assay Editor interface and Analytes template in well



ORDERING INFORMATION

- Imager Bio F2**, Fluorescent MicroArray Imager, 2 channels
- Imager Bio F3**, Fluorescent MicroArray Imager, 3 channels
- Imager Bio F4**, Fluorescent MicroArray Imager, 4 channels

Cat. number

- BSM000102-A01
- BSM000103-A01
- BSM000104-A01



HiPo MPP-96, Microplate Photometer

DESCRIPTION

Microplate Photometer HiPo is a compact tabletop device for measuring optical density — results of ELISA and micro-biological studies in 96-well microplates. Photometer is controlled and outputs data via computer. An extensive range of additional interference filters is available (with average increment of 10 nm).

The device is supplied with specialized software **QuantAssay**. Features of **QuantAssay** software:

- ELISA assays of any complexity can be carried out via robust assay editor with help of Assay Wizard
- Quantitative assay, includes up to 20 standards
- Avidity/Affinity assays
- Multiplex assays with up to 7 assays on one plate
- Qualitative assay includes up to 11 controls
- BestFit function for selecting the best calibration curve
- User-friendly interface: get your results in 3 clicks
- Save, load and export results
- LIMS export integration
- Creates visual reports



SPECIFICATIONS

Detection mode	Absorbance
Light source	LED, self-calibrating
Photodetector	8 silicon photodiodes
Plate type	96-well microplates (including strip-well microplates)
Reading Speed	5–8 s per wavelength
Measurement modes	Endpoint, Kinetic
Measurement channels	8
Reference channel	1
Measurement range	0–4.3 OD
Resolution	0.0001 OD
Wavelength range	400–700 nm
Wavelength selection	up to 8* filters on wheel standard filters 405, 450, 492 and 620 nm
Shaking	4 amplitudes, 4 speeds
Software	QuantAssay / compatible with IDEXX xChekPlus™
PC system requirements	Intel/AMD Processor, 1 GB RAM, Windows Vista/7/8/10/11, USB
Overall dimensions (WxDxH)	140 × 300 × 130 mm
Weight	4.6 kg
External power supply	Input AC 100–240 V 50/60 Hz, Output DC 12 V



Accuracy (405, 450, 492, 620 nm)

0.000–2.000 OD	≤ (0.5 % ±0.010 OD) typical
2.000–3.000 OD	≤ (1 % ±0.010 OD) typical

Precision / Reproducibility (405, 450, 492, 620 nm)

0.000–2.000 OD	≤ (0.5 % ±0.005 OD)
2.000–3.000 OD	≤ (1.0 % ±0.005 OD)

ORDERING INFORMATION:

Cat. number

HiPo MPP-96 BS-050108-A02

Optional accessories:

OD Plate, Verification tool BS-050108-AK

Additional filters* On request

* — It is possible to install up to 4 additional filters on request. Additional filters are available in two specifications: optical absorption not less than 3.5 OD or 4.3 OD

Quant Assay, Software for HiPo MPP-96



Software video is available on the website

ELISA assays of any complexity can be carried out via robust assay editor with help of **Assay Wizard**:

Qualitative assay includes up to 11 controls; Results can be outputted as Positive/Negative or Positive/Gray Zone/Negative; Gray zone can be set as symmetric and non-symmetric; Positivity ratio can be outputted

Avidity/Affinity results be outputted as Positive/Negative or Positive/Gray Zone/Negative; Avidity index margins can be easily set; Avidity Index can be outputted

	Margin	Result
If AI <	0.30	+
If AI >=	0.30 and 0.50	++
If AI >=	0.50	+++

User-friendly interface: get your results in 3 clicks: Choose an assay, a template and press Play

Choose an assay
 Quantitative ▶

Choose a Template
 Usual template ▶


Save, load and export results
 Creates reports: Excel, PDF, CSV



Quantitative assay includes up to 20 standards; User can choose Standard/Reverse type of curves

BestFit function for selecting the best calibration curve from following models:

4/5 Parameters logistics, Piece-wise linear, Linear, Index/Logarithm/Exponent/Cubic regression models

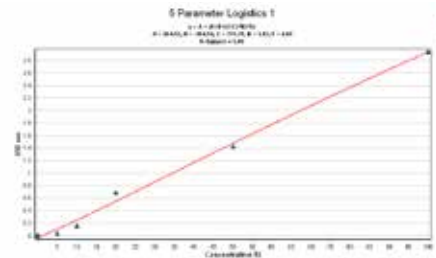

Install up to 7 assays on one plate by using **multiplex**

	1	2	3	4	5	6	7
A	Samp1	Samp1	Samp1	Samp1	Samp1	Samp1	Samp1
B	Samp2	Samp2	Samp2	Samp2	Samp2	Samp2	Samp2
C	Samp3	Samp3	Samp3	Samp3	Samp3	Samp3	Samp3
	0	1	2	3	4	5	6

Easy fill of the samples

Name: Samp 2 Bkg P₁ N₁ Std X

Group: 2 Reset

	1	2	3	4	5	6	7	8	9
A	Samp1	Samp1	1.296	1.368	1.915	1.814	1.581	1.633	2.592

PDF report contains: Experiment information, Results table, List of variables and it's calculations, Interpretation parameters


OD Plate, Verification Instrument for MPP-96 HiPo



DESCRIPTION **OD Plate** is the measurement verification instrument for microplate photometer MPP-96 HiPo. The instrument is designed to verify the accuracy and precision of measurements of the photometer at six levels of nominal optical density: 0.3; 0.6; 1.0; 2.0; 3.0; 4.0 OD. The instrument is supplied with the following verification wavelength range: 405–700 nm.

Instrument is provided in a shockproof container with an USB flash drive containing:

- Copy of measurement results
- User manual

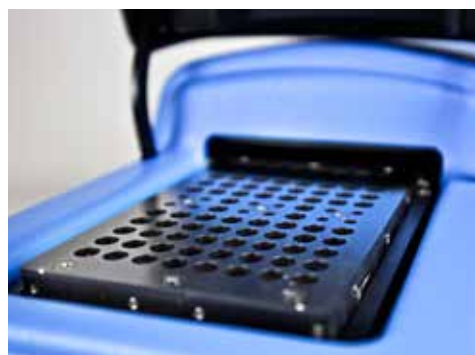

SPECIFICATIONS


Nominal optical density levels	0.3; 0.6; 1.0; 2.0; 3.0; 4.0 OD (±0.1 OD)
Verification wavelength range	405, 414, 450, 480, 492, 515, 540, 550, 560, 568, 580, 594, 620, 630, 650, 690, 700 nm
Instrument dimensions	128 × 86 × 12 mm
Net weight	0.3 kg

ORDERING INFORMATION:

OD Plate, Verification tool

Cat. number

BS-050108-AK

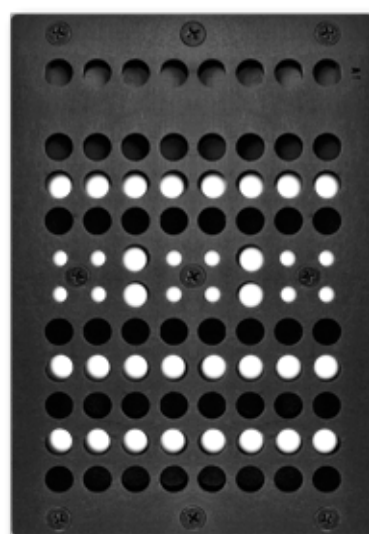

General Information

Safety

All Biosan laboratory equipment meets the requirements of International Standard IEC 61010-1: Safety requirements for electrical equipment for measurement, control and laboratory use, and applicable specific parts e.g. IEC 61010-2-010: particular requirements for laboratory equipment for the heating of materials, IEC 61010-2-020: particular requirements for laboratory centrifuges, IEC 61010-2-051: particular requirements for laboratory equipment for mixing and stirring.

CE Mark

All Biosan laboratory equipment bears a CE mark to indicate that it meets the requirements of all applicable European Directives.

Compliance with the Low Voltage Directive is demonstrated by meeting EN 61010 (as indicating in paragraph on safety) and the EMC Directive by meeting EN61326-1: EMC requirements for electrical equipment for measurement, control and laboratory use. Some products also fall within the scope of IVD Directive.

Electrical Supplies

All standard Biosan laboratory equipment is available for voltages within the range 220–240 V, 50 or 60 Hz. Most of the equipment is also available for voltages 100–120 V, 50 or 60 Hz.

Quality

The Biosan Quality Management System complies with the requirements of LVS EN ISO 9001:2015, the scope of certification is development, production, sales and service of laboratory equipment.

High quality customer service and readiness to meet ever growing customer requirements to modern equipment are the main goals of ISO 9001 compliance (certified since 2004). Since 2022 Biosan Quality Management System complies with the requirements of LVS EN ISO 13485:2016, the scope of certification is development, design, production, service and distribution of medical devices: device for measuring optical density and automatic microplate washer, sales, storage and distribution of active and non active nonimplantable medical devices.

Environmental Conditions

Biosan laboratory equipment is designed for operation in cold rooms, incubators (excluding CO₂ incubators) and closed laboratory rooms at ambient temperature from +4°C to +40°C in a non-condensing atmosphere and maximum relative humidity 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C.

Guarantee and After Sales Service

Biosan equipment is reliable, designed and built to provide years of trouble-free service. Most Biosan equipment is guaranteed for two years against faulty materials and workmanship and premium product class up to 3 years upon registration in our website support section. Warranty terms and conditions are indicated in the product manual. All Grant standard laboratory equipment is guaranteed for three years against faulty materials and workmanship. Local distributors and service centres provide necessary technical assistance within and outside the warranty period. Biosan technical support team provides direct support offering the best solution for assistance upon receipt of request via e-mail service@biosan.lv or forms available at Technical Support section of Biosan web-site.

World Wide Availability and Support for Biosan Laboratory Equipment

Biosan laboratory equipment and specialist technical support is available world-wide. Please, visit multilingual (English, French, German, Italian, Latvian, Russian, Spanish) web-site <http://www.biosan.lv> for further product information (videos, brochures, manuals, articles), placing enquiries and locating your locally appointed distributor or contact customer service direct at service@biosan.lv.

As Biosan is committed to a continuous program of improvement, specifications may be changed without notice.



Basic Plus



Premium



Smart Plus

PRODUCT CLASS FEATURES

	Basic Plus	Premium	Smart Plus
Designed to complete basic sample preparation tasks	●	●	●
Designed to complete sophisticated sample preparation tasks		●	●
Advanced specifications and special features		●	●
PC interface for logging, control, programming, alarms, online monitoring functions			●
Modern Bioform design	●	●	●
Small footprint	●	●	●
Low power consumption	●	●	●
Safe 12V DC	●	●	●
High quality	●	●	●
2 year warranty + 3rd year purchased via distributors	●		●
2 year warranty + 3rd year for free upon product registration		●	

Applications and Articles



Reverse-Spin® Technology — Innovative Principle of Microbial Cultivation.

biosan.lv/rts-tech



Investigation of the effect of aeration on growth dynamics, respiratory rate and pH changes of the aerobic bacterium *E. coli* BL21 cultivated in RTS-8 PLUS single-use bioreactor.

biosan.lv/rts-growth



Development and evaluation of DNA amplicon quantification.

biosan.lv/uv-box



UVR-M and UVR-Mi, UV Air Recirculators Test Report.

biosan.lv/uvr-test



Germicidal and Antiviral decontamination of air by UV irradiation and UV recirculator method.

biosan.lv/uv-effect



how to choose

A PROPER SHAKER, ROCKER, VORTEX



Medical-Biological
Research & Technologies

Sample volume
 $10^3 \dots 10^2$ ml

Erlenmeyer flask
and Cultivation flask



Sample volume
 10^1 ml

Petri dishes, vacutainers
and tubes up to 50 ml



Sample volume
 $10^0 \dots 10^{-3}$ ml

PCR plates, microtest plates
and Eppendorf type tubes



PSU-20i,
Orbital Shaker

ES-20/80,
Orbital Shaker-Incubator



Applications:

- Microbiology
- Extraction
- Cell cultivation



PSU-10i,
Orbital Shaker



ES-20,
Orbital
Shaker-Incubator

Applications:

- Agglutination
- Gel staining/destaining



MR-12,
Rocker-Shaker



Multi RS-60,
Programmable rotator

Bio RS-24,
Mini-Rotator



RTS-1 and RTS-1C,
Personal bioreactor



MR-1,
Mini Rocker-Shaker



Multi Bio 3D,
Mini Shaker

Applications:

- Agglutination
- Extraction
- Blot hybridisation
- Gel staining/destaining



Multi Bio RS-24,
Programmable rotator

Applications:

- Microbiology
- Extraction
- Cell cultivation
- Hematology



V-1 plus,
Vortex



MSV-3500,
Multi Speed Vortex

Applications:

- Nucleic acid Analysis
- Molecular Analysis
- Protein Analysis
- Genomic Analysis



PST-60HL-4,
Thermo-Shaker



PST-60HL,
Thermo-Shaker



MPS-1,
Multi Plate Shaker



CVP-2,
Centrifuge vortex for PCR plates

TS-100, TS-100C, TS-100C Smart
Thermo-Shakers

V-32,
Multi-Vortex



PST-100HL,
Thermo-Shaker

TS-DW,
Thermo-Shaker
for deep well
plates



Applications:

- ELISA Analysis
- Genomic Analysis
- Hybridization
- Immunology



PSU-2T,
Mini-Shaker




COVID-19

The global coronavirus (COVID-19) outbreak marks the necessity of fast and reliable sample preparation as well as safe working environment. To facilitate products selection here are shown ready product lines to rapidly start COVID-19 sample analysis.

Preventive

Decontamination

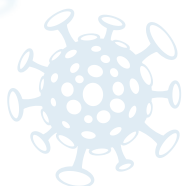

UVR-M/UVR-Mi
UV Cleaner-Recirculators

Air



PDS-250, PDS-10L
DNA/RNA decontamination solution

Surface

ELISA based detection methods

Steps

Incubation



PST-60HL-4
Plate Shaker-Thermostat



PST-60HL
Plate Shaker-Thermostat



PSU-2T
Mini-shaker
for immunology

Washing



Assist
Pipette series



Assistboy
Pipette
controller



FTA-2i
Aspirator
with Trap Flask



3D-IW8
Inteliwasher

Detection and analysis



HiPo MPP-96
Microplate
Photometer

+



Quant Assay
Software for
HiPo MPP-96

NA based detection methods

Manual preparation



UVT-S-AR
DNA/RNA UV-cleaner box

Working area

Resuspension of probes and reagents



V-1 Plus
Personal Vortex



MSC-3000
Centrifuge/Vortex Multispin



TDB-120
Dry block thermostat

Sample lysis



Assist
Pipette series



FTA-2i
Aspirator with Trap Flask

Sample wash



MagSorb 16
Magnetic Rack for Manual Nucleic Acid Extraction



Microspin 12
High-speed Mini-centrifuge

Elution of NA

Automated preparation



TS-100C Smart
Programmable Thermo-Shaker with cooling for microtubes and PCR plates



Biomagpure 12 Plus
Compact Bench-Top Robotic Workstation for Automated Nucleic Acid Purification (+ reagents)

Sample detection and analysis

Preparation



CVP-2
Plate centrifugation and vortexing



MSC-6000
Tube centrifugation and vortexing

Real time PCR detection



BioQuant-96
Real time PCR instrument



abTes COVID-19 qPCR kit
(Reagents are CE-IVD certified – among WHO listed kits.)

Manual and automatic DNA/RNA extraction

Liquid Handling

Assist Series
Pipettes



Incubation

TS-100, TS-100C or TS-100C Smart



Sample resuspension
and droplet spin down

FVL-2400N



or

Mixing and
resuspension

MPS-1



Automatic
extraction

BioMagPure 12 Plus
and reagents

Centrifugation

Microspin 12



Sample preparation
in UV-Cabinet for PCR

UVC/T-M-AR



Centrifugation

FTA-1



FTA-2i



Incubation

TDB-120



Reagents

Reagents for extraction*



* — Information about current offers on the products of other manufacturers are available in the corresponding sections of our site www.biosan.lv/en/products

Manual DNA/RNA extraction using magnetic beads technology

Sample resuspension and droplet spin down

FVL-2400N



Sample preparation in UV-Cabinet for PCR

UVC/T-M-AR, or similar,
see UV-Cabinets for PCR



Mixing and resuspension

MPS-1



V-1 plus



Multi Bio RS-24



Capture of magnetic beads

MagSorb-16



Centrifugation

Microspin 12



Vacuum aspiration

FTA-1



FTA-2i



Incubation

TDB-120



TS-100C



Reagents

Reagents for extraction*



* — Information about current offers on the products of other manufacturers are available in the corresponding sections of our site www.biosan.lv/en/products

PCR Analysis

Sample resuspension
and droplet spin down



FVL-2400N
(tubes)



CVP-2
(PCR-plates)



Real time
PCR detection

BioQuant-96



Reagents

abTes COVID-19 qPCR kit
(Reagents are CE-IVD certified
— among WHO listed kits.)



Reagents

* — Information about current offers on the products of other manufacturers are available in the corresponding sections of our site www.biosan.lv/en/products

Enzyme-linked immunosorbent assay (ELISA)

Sample preparation
in laminar flow cabinet



Biological Safety Cabinets class II



PST-60HL



PST-60HL-4



PSU-2T

Incubation



IW-8



3D-IW8

Washing (Automated)



FTA-1 with MA-8



FTA-2i

Washing (Manual)



HiPo MPP-96
Microplate Photometer
with QuantAssay software

+



Reagents of ELISA*

Reading and Analysis

* — Information about current offers on the products of other manufacturers are available in the corresponding sections of our site www.biosan.lv/en/products

Microbial Cell Cultivation

Sample preparation
in laminar flow cabinet

Biological Safety Cabinet class II



ES-20



ES-20/60

Cultivation



ES-20/80



ES-20/80C

CO₂ cultivation



S-Bt Smart Biotherm

Measurement



DEN-1



DEN-1B



DEN-600

or

Cultivation and
Real-time OD
Measurements and
Logging



RTS-1



RTS-1C



RTS-8 and RTS-8 plus



Biosan, SIA

Ratsupites iela 7 k-2,

Rīga, LV-1067, Latvia

Phone: **+371 674 261 37**

E-mail: [**marketing@biosan.lv**](mailto:marketing@biosan.lv)

[**www.biosan.lv**](http://www.biosan.lv)

